

ApneaLink Plus featuring EasySense technology

Sleep-disordered breathing (SDB) is recognized as a serious health problem that impacts approximately 43 million US adults. More than 80% remain undiagnosed, and many barriers prevent patients from getting access to therapy.¹

Now, the detection of this chronically debilitating condition has been made easier with the ApneaLink Plus, a Type III portable monitoring device, the latest addition to the ApneaLink family of diagnostic products and accessories.

The ApneaLink devices provide you with a cost-effective, easy-to-use method of diagnosing or screening patients for obstructive sleep apnea (OSA) in the home. The device reports apneas, hypopneas, flow limitation, snoring, blood oxygen saturation and the probability of Cheyne–Stokes respiration (CSR) breathing patterns within the recording.

The ApneaLink and ApneaLink Plus have been validated in several studies worldwide

Validation of MicroMESAM² as screening device for sleep-disordered breathing

(Wang Y, Teschler T, Weinreich G, Hess S, Wessendorf TE, Teschler H)

MicroMESAM-generated flow-time curves correspond well with pneumotachograph-generated curves, producing automated AHI³ that are highly sensitive in detecting SDB.

Validation of ApneaLink as screening device for Cheyne–Stokes respiration.

(Weinreich G, Armitstead J, Töpfer V, Wang YM, Teschler, H)

The study demonstrated that the screening classifier was able to detect CSR with high diagnostic accuracy. Hence, the ApneaLink equipped with the CSR classifier is an appropriate screening tool that may help to prioritize patients with CSR for a polysomnography (PSG).

Validation of the ApneaLink for the screening of sleep apnea: A novel and simple single-channel recording device

(Erman MK, Stewart D, Einhorn D, Gordon N, Casal E)

The ApneaLink device provides reliable information, is simple to use, and is highly sensitive and specific in calculating AHI when compared with the AHI obtained from a full PSG.

¹ Young et al. *Am J Respir Crit Care Med* 2002

² Distributed by ResMed as the ApneaLink in the US

³ Values reported as AHI for MicroMESAM (ApneaLink) are actually RDI values (AHI plus flow limitation). See full translated article for details.



RESMED

Ordering Information and Product Codes

ApneaLink Plus Complete Set (includes respiratory effort and oximetry)

- ApneaLink Plus recorder device
- Program CD software
- Quick software setup guide
- USB download cable
- ResMed EasySense respiratory effort sensor
- 2 reusable belts
- XPOD oximeter sensor clip
- ResMed XPOD oximeter
- 3 single-use oximeter sensors
- 3 nasal cannulas
- Carrying case
- 2 AA batteries

US, Latin America and Canada **22328**

ApneaLink Plus Basic Set (includes respiratory effort)

- ApneaLink Plus recorder device
- Program CD software
- Quick software setup guide
- USB download cable
- ResMed EasySense respiratory effort sensor
- 2 reusable belts
- 1 nasal cannula
- Carrying case
- 2 AA batteries

US, Latin America and Canada **22319**

ApneaLink Basic Set

- ApneaLink recorder device
- Program CD software
- Quick software setup guide
- USB download cable
- 1 reusable belt
- 3 nasal cannulas
- Carrying case
- 2 AA batteries

US, Latin America and Canada **22302**

ApneaLink Oximetry Accessories Kit

- ResMed XPOD oximeter
- XPOD oximeter sensor clip
- 3 single-use oximeter sensors

US, Latin America and Canada **22304**

ApneaLink Plus Accessories (for use with ApneaLink Plus only)

- ResMed EasySense Respiratory effort sensor **22333**
- Belt, reusable stretch (required) **629052**

ApneaLink Plus and ApneaLink Software

(compatible with all devices) **22326**

Other Optional Accessories and Disposables

- Belt, single use (24 pack) (for use with ApneaLink only) **70406**
- Nasal cannulas (25 pack) **70388**
- Nasal/oxygen cannulas (25 pack) **22976**

Oximetry Components

- Oximeter soft sensor, reusable (recommended) **70413**
- Oximeter sensor – single-use **22337**
- XPOD oximeter sensor fixation clip **22306**

RESMED

ApneaLink™ / ApneaLink Plus

Portable Monitoring Devices

TECHNICAL SPECIFICATIONS

Signal Recording

- Respiratory effort
- Respiratory flow
- Breathing sounds
- Blood oxygen saturation
- Pulse
- Battery voltage

Sampling Rates for the Channels

- Respiratory flow / breathing sounds: 100 Hz
- Blood oxygen saturation: 1 Hz
- Pulse: 1 Hz
- Battery: 1 Hz
- Respiratory effort: 10 Hz

Signal Processing

- Signal recording: 20 Bit
- Signal storage: 16 Bit

Internal Memory

- Storage capacity: 15 MB
- Recording period: 8 hours minimum

Power Supply to Recorder

- 2 batteries: LR 6 / Mignon / AA / 1.5 V / at least 2100 mAh
- 2 NiMH rechargeable batteries: Mignon / AA / 1.2 V / at least 2100 mAh

Dimensions (length x width x height)

- Recorder: 4.6" x 2.4" x 1.2" (125 x 60 x 30 mm)
- Pulse oximeter: 2.1" x 0.8" x 0.6" (53 x 20 x 15 mm)

Weight (recorder without batteries)

- Recorder (without batteries): Approx. 50 g (1.8 oz)
- Pulse oximeter: approx. 30 g (1.1 oz)

Operating Conditions

- Temperature: 68°F to 104°F (20°C to 40°C)
- Humidity: 10% to 90% RH (non-condensing)

Shipment/Storage Conditions

- Temperature: -4°F to +122°F (-20°C to +50°C)
- Humidity: 10% to 90% RH

Operating/Storage Air Pressure

- 800 hPa to 1060 hPa

Effective Range

- Flow sensor: -10 hPa to +10 hPa
- SpO₂: 70 to 100%
- Pulse: 18 to 300 bpm

Accuracy (No Movement)

- SpO₂: +/- 3 digits
- Pulse: +/- 3 digits

Interfaces

- Nasal pressure cannula: Luer connection
- Pulse oximeter: 3-pin binder plug
- Computer: Full speed USB 1.1



Simple, fast and easy to use.

New features for clearer diagnoses, saving you time and money.

ResMed.com

ResMed's ApneaLink devices are the easy choice in OSA diagnosis

The ApneaLink™ improves patient care by providing easy access to treatment while helping you grow your sleep apnea business.

ApneaLink features:

Automatic analysis derives apnea-hypopnea index (AHI), hypopnea index (HI), flow limitation, snoring and oxygen desaturation index (ODI)

Validated results meet AASM and CMS definitions for hypopnea scoring guidelines

Cheyne-Stokes probability detection determines when to refer patients for further in-lab diagnosis

Results can be scored manually for more detailed patient data

Your business logo can be added to increase brand awareness

Email summary and signal reports can be sent to referral physicians or other relevant parties

Patient instructions can be printed when programming the device

Extended report contains additional overview of respiratory data



Feature Comparison	ApneaLink	ApneaLink Plus
Apnea-hypopnea index	•	•
Risk indicator	•	•
Apnea index	•	•
UAI (Unclassified apnea index)	•	•
OAI (Obstructive apnea index)	•	•
CAI (Central apnea index)	•	•
MAI (Mixed apnea index)	•	•
Hypopnea index	•	•
Flow lim br without sn (FL)	•	•
Flow lim br with sn (FS)	•	•
Snoring events	•	•
ODI (Oxygen desaturation index)	Optional with oximetry	•

ResMed's ApneaLink Plus is simple, fast and easy to use

The ApneaLink Plus with EasySense technology, a unique respiratory effort sensor, is a simple, low-cost portable home sleep test diagnostic device that records up to four channels of information: respiratory effort, pulse, oxygen saturation and nasal flow.

- Effort belt with EasySense respiratory effort sensor
- Enhanced recorder light status; improved start/stop button
- Simple, easy-to-use component connectors
- Same robust design as ApneaLink

ApneaLink Plus additional features:

Differentiation of apneas leads to clearer diagnosis and more accurate, effective reports

New prescription page streamlines process for health care professionals

AHI graphic and risk indicator can highlight either AHI or RI

Configurable analysis parameters allow for the adjustment of obstructive and central apnea thresholds

Five measurements of oxygen saturation including ≤ 89 and ≤ 88 allow for accurate billing documentation



ApneaLink Report

Simple, cost-effective and reliable results.

YOUR LOGO HERE

ApneaLink - Report of 10/6/2011 12:14 PM

Treating physician		Referral to	
<input type="text"/>		<input type="text"/>	
Patient data			
First name:	OSA	Patient ID:	
Name:	ApneaLink Plus	DOB:	8/23/2008
Street:		Size:	0 ft 0 in
City, ST, Zip:		Weight:	0.00 lbs
Phone:		BMI:	kg/m ²
Recording		Evaluation	
Date:	8/23/2008	Start:	11:43 PM
Start:	11:33 PM	End:	3:42 AM
End:	3:44 AM	Duration:	3 h 59 min
Duration:	4 h 11 min		

AHI*

Points evaluation from AHI (see Clinical Guide for more details)

Indices		Normal	Result
AHI*	62	< 5 / h	Average breaths per minute [bpm]: 9.95
RI*	65	< 5	Breaths: 2376
Apnea index:	45	< 5 / h	Apneas: 179
UAI:	0		Unclassified apneas: 0 (0%)
OAI:	36		Obstructive apneas: 143 (80%)
CAI:	6		Central apneas: 25 (14%)
MAI:	3		Mixed apneas: 11 (6%)
Hypopnea index:	17	< 5 / h	Hypopneas: 66
% Flow lim. Br. without Sn (FL):	19	< Approx. 60	Flow lim. Br. without Sn (FL): 451
% Flow lim. Br. with Sn (FS):	13	< Approx. 40	Flow lim. Br. with Sn (FS): 304
			Snoring events: 1689
ODI Oxygen Desaturation Index*:	54	< 5 / h	No. of desaturations: 215
Average saturation:	89	94% - 98%	Saturation [] 90% : 147 min (61%)
Lowest desaturation:	66	-	Saturation [] 85% : 15 min (6%)
Lowest saturation:	66	90% - 98%	Saturation [] 80% : 9 min (4%)
Baseline Saturation:	92	%	Saturation [] 89% : 108 min (45%)
			Saturation [] 88% : 68 min (28%)
Minimum pulse frequency:	58	50 - 70 bpm	
Maximum pulse frequency:	83	60 - 90 bpm	
Average pulse frequency:	68	bpm	
Proportion of probable CS epochs:	0	0%	

Analysis status: Analyzed automatically

Analysis parameters used (Default)
 Apnea [20%: 10s; 80%: 20s; 20%: 60%: 8%]; Hypopnea [70%: 10s; 100%: 1.0s]; Snoring [6.0%: 0.3s; 3.5s; 0.5s]; Desaturation [4.0%]; CSR [0.50]

Comments

OSA ApneaLink Plus - 10/6/2011 12:14 PM
 * See Clinical Guide for abbreviations and ResMed standard parameters Firmware version: 05.0000R01 Software version: 9.00

Detailed Signal View

Full night view

Signals:

Nasal flow

Snore

Adjustable signal amplitude

Respiratory effort

Oxygen saturation

Pulse

Selectable time

Expanded view


Signals (Selectable):

Nasal flow

Respiratory effort

Oxygen saturation

Pulse

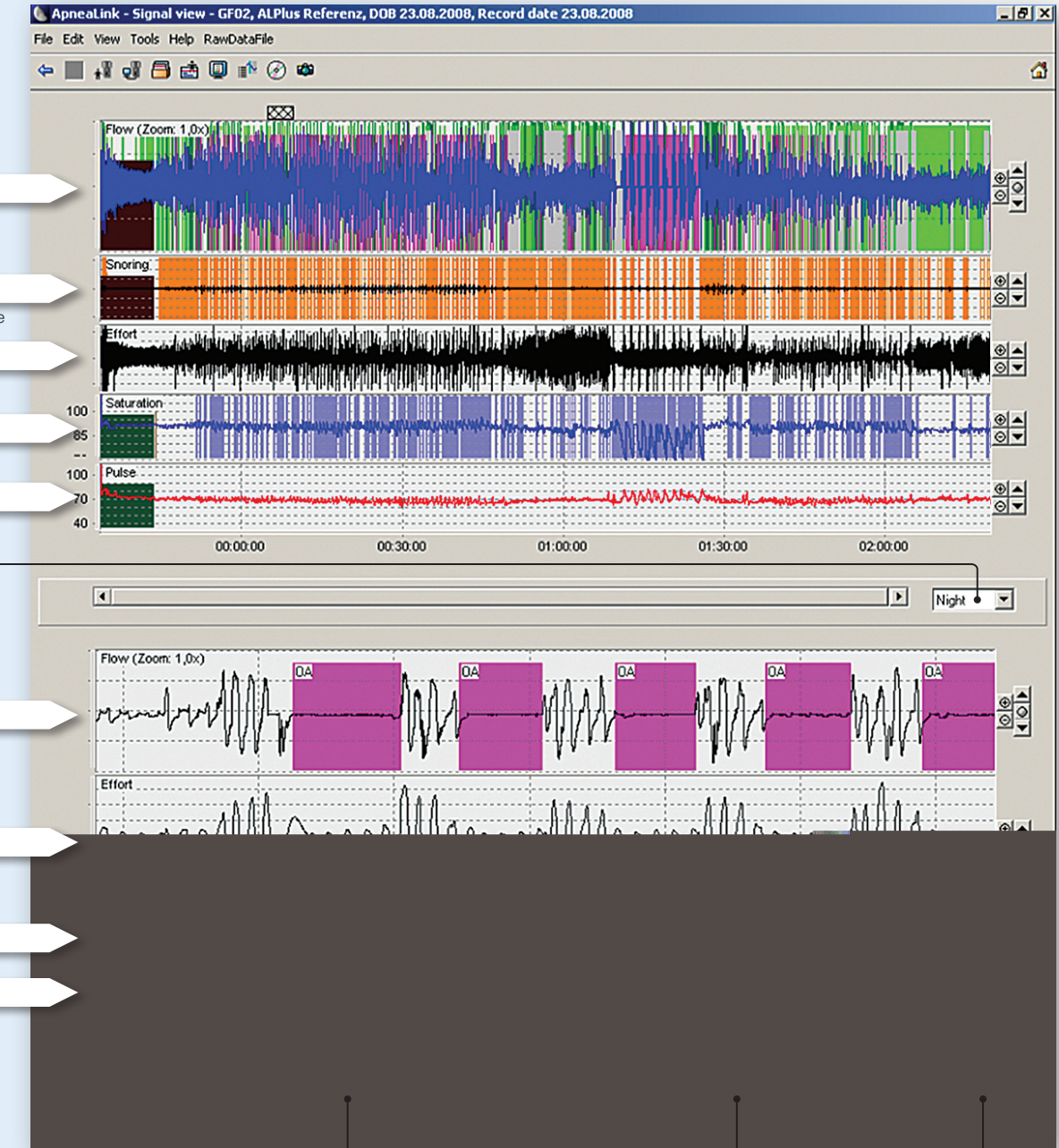

Time select bar

Event type selection

Selectable time