

## Mirage Quattro™

### FULL FACE MASK

#### FOURTH GENERATION TECHNOLOGY

**Wide fit range:** fourth generation full face mask from ResMed fits more than 95% of users.\*

**Mirage™ dual-wall cushion with spring air technology:** unique cushion accommodates movement during sleep and provides an effective solution for jaw drop.

**MicroFit dial:** new forehead support adjustment offers personalised fitting and comfort.

**Light, quiet and easy to use:** these benefits are optimised with streamlined mask design, improved vent design and a range of user-friendly features.

**MicroFit dial**  
provides better fit, seal and comfort for over 95% of users\*

**Set-and-forget headgear slots**  
save optimal fitting settings

**Improved built-in vents**  
provide minimal noise output

**Moving cushion membrane at nasal bridge**  
accommodates different types of nasal structures and prevents leak into the eyes


**Cushion clip with tabs**  
allows easy removal of cushion from frame

**Mirage dual-wall cushion with spring air technology**  
enhances seal where it's needed most

**Moving cushion membrane at chin**  
maintains seal and accommodates jaw drop and movement

**Set-and-forget headgear clips**  
allow convenient removal of mask

**Simplified elbow with anti-asphyxia valve (AAV)**  
few elbow parts for easy assembly and disassembly; AAV enhances safety if the airflow to the mask is impeded

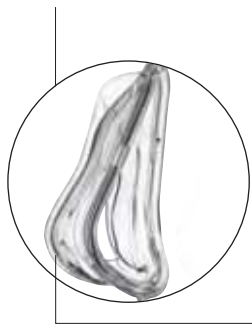

The membrane as seen on a side view of the mask cushion





## TECHNOLOGY HIGHLIGHT

### MicroFit dial

A choice of 24 positions at your fingertips—providing a comfortable, effective fit for over 95% of users\*

- **Adjustable:** can be adjusted while on the user's face
- **Precise:** each turn of the dial adjusts the forehead support in small increments
- **Progressive:** moves in towards the face, allowing the user to sense when the fit is most comfortable for them
- **Works in combination with the cushion:** as the dial adjusts the angle of the mask to the face, the cushion adapts to enhance the seal where it's needed most—without the need to overtighten the headgear

## USER BENEFITS

### Personalised comfort

- MicroFit dial adjusts the mask to suit user's profile
- Dual-wall cushion with spring air technology adapts to facial contours and accommodates movement during sleep (including jaw drop)

### Light and quiet

- Streamlined design feels light on the face
- Improved vent design reduces noise

### Making it easier

- MicroFit dial and spring air cushion technology help achieve good mask fit—first time, every time
- Fewer parts for easy handling, cleaning, assembly and disassembly

## PROFESSIONAL BENEFITS

### Fit range increased to over 95% of users\*

- MicroFit dial and spring air cushion technology adjust the mask to suit more facial structures
- New XS mask size

### Customer satisfaction and compliance

- Comprehensive user support package with travel bag, mask care DVD, quick-fitting guide and user guide
- Increased therapy acceptance through first-time fit
- Personalised fit, comfort and seal encourage long-term compliance
- Simplified parts for quick and easy assembly and disassembly

### Cost and time efficiencies

- Reduced inventory: 33% fewer product codes than Ultra Mirage™ Full Face Mask and fewer mask parts
- Reduced fitting and training time: improved cushion technology seals quickly
- Reduced callbacks: ease-of-use and fitting

## PRODUCT CODES

### Mirage Quattro Full Face Mask

#### Europe 1<sup>1</sup>

X-Small	61204
Small	61205
Medium	61206
Large	61207

#### Europe 2<sup>2</sup>

X-Small	61208
Small	61209
Medium	61210
Large	61211

#### Europe 3<sup>3</sup>

X-Small	61228
Small	61229
Medium	61230
Large	61231

#### France

X-Small	61220
Small	61221
Medium	61222
Large	61223

#### UK

X-Small	61224
Small	61225
Medium	61226
Large	61227

#### Asia-Pacific

X-Small	61216
Small	61217
Medium	61218
Large	61219

#### ANZ

X-Small	61212
Small	61213
Medium	61214
Large	61215

<sup>1</sup> English, German, French, Spanish, Italian, Portuguese and Dutch

<sup>2</sup> English, Swedish, Norwegian, Danish, Finnish, Hungarian, Polish and Greek

<sup>3</sup> English, Russian, Turkish, Czech, Estonian, Slovenian, Slovakian, Bulgarian, Romanian, Latvian and Lithuanian

\* Source: internal ResMed anthropometric data