

Swift™ LT

NASAL PILLOWS SYSTEM

LESS IS MORE

Minimal contact, optimal stability

Freedom of movement and side-sleeping

Super-quiet performance

Flexible fit, personalised comfort

**First nasal pillows
with dual-wall technology**


Dual-wall nasal pillows
seal on contact for a gentle, comfortable seal; thin outer layer enhances comfort, thick inner layer provides stability

Flexible cushion base
enables the nasal pillows to move in multiple directions without compromising seal

360° rotating elbow
provides flexibility in tube position

Quiet vents
diffuse and direct airflow quietly away from the user and partner

Flexible lightweight tube
minimises 'pull' on the mask, increasing stability

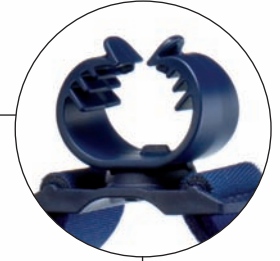

Unique removable tube retainer on top buckle offers three lock-in positions to set the desired tube direction with a simple twist action

Back strap with soft buckle offers a low-profile easy-to-adjust, flexible fastener

Stabilising elements are flexible and adjust to individual facial widths while keeping the Swift LT gently in place

Click-fit adjustor helps select the most comfortable angle for the nasal pillows on the face and keeps them in place

Narrow soft silicone frame makes sleeping on the side more comfortable



TECHNOLOGY HIGHLIGHT

Minimal contact, optimal stability



How it works

- Swift LT offers a unique combination of dual-wall nasal pillows on a flexible cushion base applying minimal contact on the user's face.
- The click-fit adjuster and stabilising elements hold the cushion and pillows in a comfortable and stable position.
- The nasal pillows can move vertically and sideways, without compromising seal.
- Dual-wall pillows have a soft thin outer layer for comfort and a thicker inner layer for stability.

The benefits

- Swift LT maintains seal with minimum facial contact.
- It adapts to movement quickly.
- It provides stability, whether the user is still or moving.
- The mask does not need to feel tight to seal effectively—it works best when the nasal pillows and cushion base are in light contact with the face.

USER BENEFITS

Freedom of movement with stability

- Dual-wall pillows and flexible cushion base maintain seal while adapting to movement.
- The flexible, lightweight tubing minimises pull on the mask and enhances stability.
- Minimalism reinforces the sense of freedom—it's as much to do with how the user feels as with how the mask looks on their face.

Personalised fit

- The click-fit adjuster enables the user to position the nasal pillows at the most suitable angle for their facial profile.

Choice of sleeping positions

- Swift LT users can sleep on their side or their back, with the option to use the tube retainer, positioning the tube as desired.

Quiet performance for quiet sleep

- Innovative vent design delivers super-quiet mask performance—only 25 decibels.*

PROFESSIONAL BENEFITS

Quick to fit

- Swift LT is quick and easy to fit, saving time for clinicians and sleep professionals.

Reduced inventory, 95% fit range**

- The system comes packed with three pillow sizes to suit most users, reducing inventory.

Easy product handling

- With very few parts, Swift LT is easy to clean and assemble.

Convenient to recommend or prescribe

- User guide and online support facilitate training for users and staff.

Patient compliance

- Designed for patient satisfaction, comfort and ease-of-use to improve compliance.

* Testing conducted per ISO 3744:1994 at 10 cm H₂O

** Based on ResMed's anthropometric studies

PRODUCT CODES

Swift LT Nasal Pillows System

Europe ¹	60561
Europe ²	60562
Europe ³	60563

¹ English (UK), Dutch, French, German and Italian (Switzerland)

² English, Danish, Finnish, Norwegian and Swedish

³ English, Czech, Dutch (Belgium), French (Belgium), German (Belgium), Greek, Hungarian, Italian, Polish, Portuguese, Russian, Slovak, Spanish and Turkish