

Swift™ FX

NASAL PILLOWS SYSTEM

There's Nothing Hard About It

The Swift FX is so soft, light and simple, it feels like it's barely there! Minimal contact on the face and unobtrusive design make therapy less intimidating and easy to accept.

- Simple, soft and stable
- Stays in place when patient moves
- Designed for optimal compliance
- Instant, intuitive fitting
- Stable at high pressures
- Helps patients sleep better

FOR THE PROFESSIONAL

Intuitive, almost instant fit

- Once the right pillows size is selected, it requires minimal or often, no adjustment.

Stable seal at high pressures

- The system stays in place and maintains seal effectively even at higher treatment pressures*

Easier therapy acceptance by patients

- Designed for instant appeal, Swift FX is non-intimidating to facilitate therapy acceptance.

Reduced inventory, wide (>95%) fit range**

- Packed with three pillow sizes to suit most users. An Extra-Small size is also available.

Easy product handling

- With very few parts, Swift FX is easy to clean and assemble.

Convenient to recommend or prescribe

- User guide and online support enable easy training for users and staff.

FOR THE PATIENT

Soft comfort all around

- From the tip of the nasal pillows to the soft, secure sides with Soft Wraps ...the entire system enhances comfort

Stable during movement

- Dual-wall pillows and flexible cushion base maintain seal while adapting to movement.
- Spring-flex lightweight tubing minimises pull on the mask, enhancing stability.

'Barely there' experience

- Minimalism reinforces the sense of freedom—in terms of how the user feels and looks.

Choice of sleeping positions

- Swift FX users can sleep on their side or their back, positioning the tube as desired.

Quiet performance for quiet sleep

- Innovative vent design diffuses airflow for quiet mask performance—only 25 dBA.***

*Tested up to 20 cm H₂O in clinical trials

**Based on ResMed's market research

***Declared dual numbered noise emission values (with uncertainty 3 dBA) in accordance with ISO 4871.





Swift FX
making hair management easy

PRODUCT CODES


Swift FX Nasal Pillows System

Europe ¹	61501
Europe ²	61502
Europe ³	61503

¹ English (UK), Dutch, French, German and Italian (Switzerland)

² English, Danish, Finnish, Norwegian and Swedish

³ English, Czech, Dutch (Belgium), French (Belgium), German (Belgium), Greek, Hungarian, Italian, Polish, Portuguese, Russian, Slovak, Spanish and Turkish

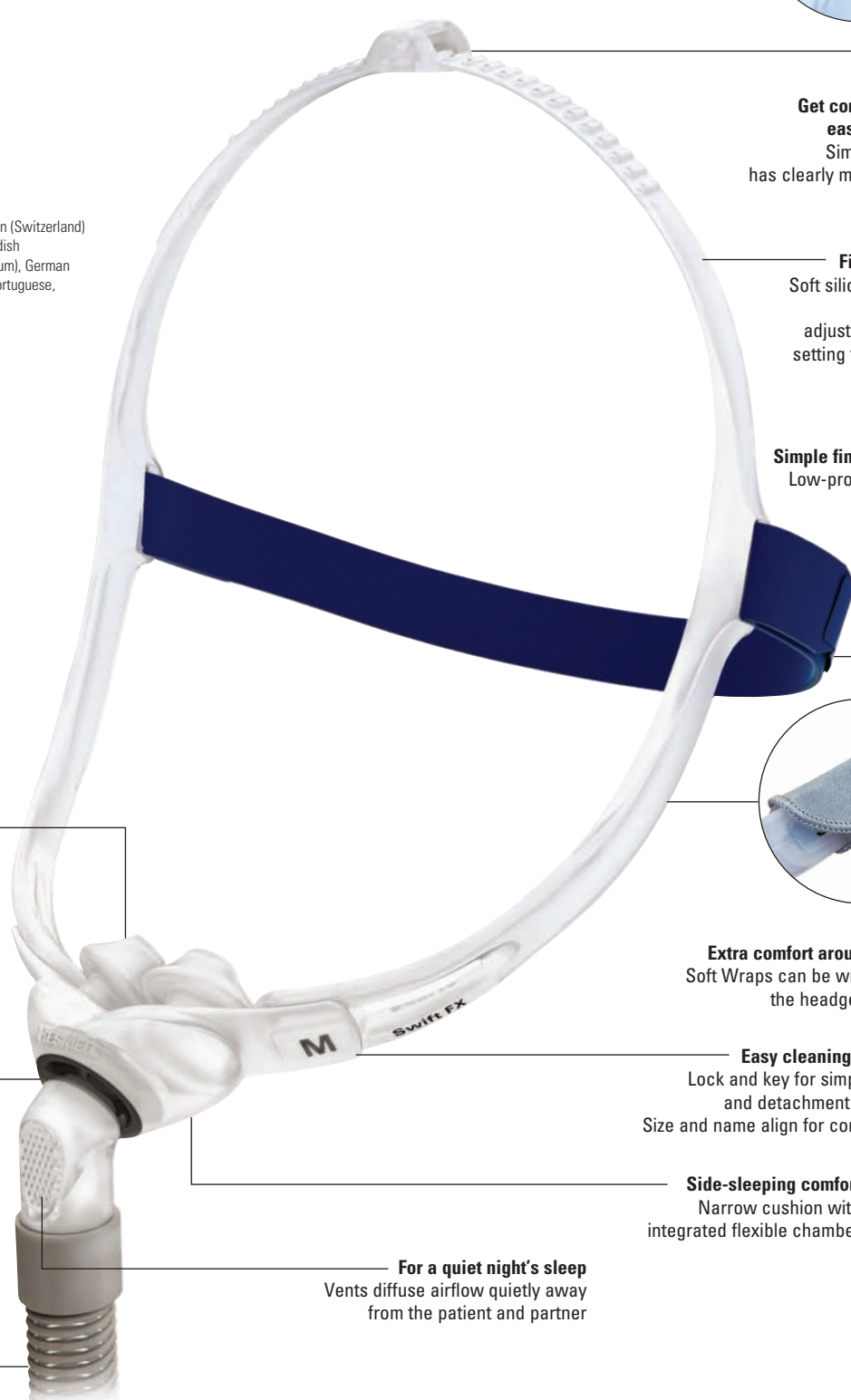

Get comfortable with easy adjustments
Simple top buckle has clearly marked notches

Fits like a glove
Soft silicone headgear needs minimal adjustment - default setting fits most users

Simple fine adjustments
Low-profile back strap with fasteners


Gentle, secure seal
Dual-wall nasal pillows seal on contact and gently disperse air to the nose



Freedom of movement, stability and comfort
Flexible pillows move in multiple directions

Ultra-smooth 360° rotating elbow swivels freely

Spring-flex lightweight tube stretches and bends easily

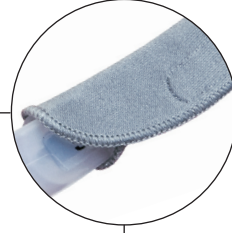

Extra comfort around the cheeks
Soft Wraps can be wrapped around the headgear side straps

Easy cleaning and assembly
Lock and key for simple attachment and detachment of side straps; Size and name align for correct assembly

Side-sleeping comfort
Narrow cushion with integrated flexible chamber

For a quiet night's sleep
Vents diffuse airflow quietly away from the patient and partner