

ResMed

LifeChoice **Activox™** 4L

Portable oxygen concentrator



The oxygen your patients need,
wherever they go

Portable therapy, quality of life.

The LifeChoice Activox 4L portable oxygen concentrator (POC) is your patient's passport to a life of greater freedom and peace of mind.

Reliable oxygen, outstanding battery life and a lightweight, low-fuss format means your patients won't feel weighed down by therapy. Instead, they can forget about their disease and focus on enjoying an active life.



Designed for mobility

Light and portable



- Weighs just 2.2 kg (4.8 lbs)
- User-friendly carry case with straps offers enhanced freedom
- Converts easily between four carrying options: briefcase, shoulder bag, hands-free backpack and hands-free waist pack

Powered to keep up

- Up to 10 hours of use with the powerful internal battery*
- An extra 4 hours of use with the optional external battery*
- AC/DC adapter delivers quick, efficient charging and operation; both at home and on the go

Simple, easy-to-use interface

- Transflective display screen makes it easy to operate at any time, even in direct sunlight and total darkness
- Control panel shows relevant day-to-day information (e.g. battery life, mode) as well as error messages and alarms
- Simplicity and ease of use inspires confidence and facilitates self-management

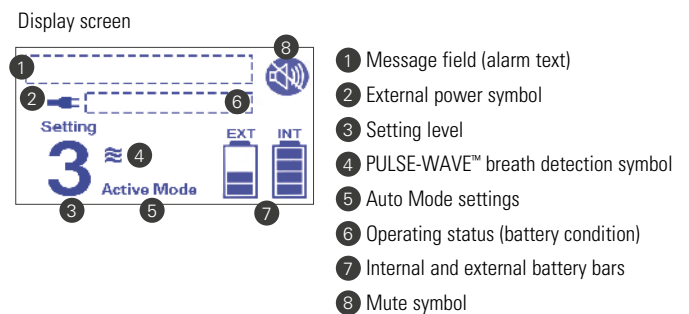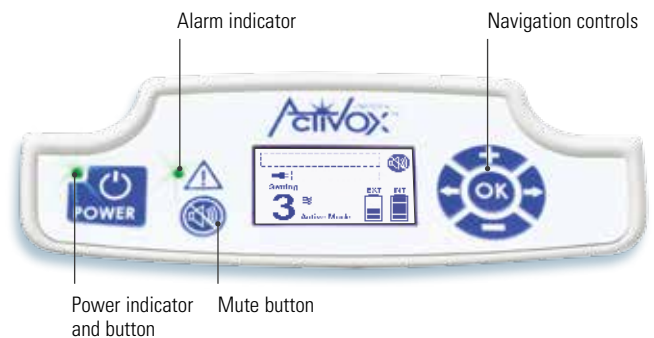




Simple therapy, maximum comfort

The LifeChoice Activox 4L is firmly focused on meeting your patients' needs and offering them the freedom, comfort and confidence that comes with portable therapeutic support.

Responsive algorithms deliver the right dose of oxygen, exactly when it's needed. The reliable, thoughtfully designed LifeChoice Activox 4L keeps pace as your patient's activity levels change during the day and their condition evolves over the long term.

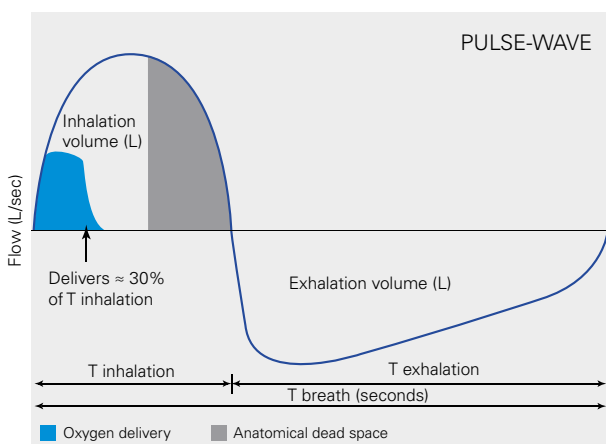

Efficient oxygen delivery

- PULSE-WAVE™ technology effectively delivers the required dose of oxygen while minimising waste and reducing the risk of nasal irritation and nose bleeds.
- Auto Mode senses changes in your patient's activity level and adapts to their breathing requirements.
- Auto Dose, a feature of Auto Mode, ensures treatment continuity. It provides a 'safety' dose of oxygen.



Right for everyone

The lightweight, attractive, user-friendly LifeChoice Activox 4L helps to increase patient satisfaction. With this reliable portable oxygen concentrator, call-outs and equipment failures are reduced, enabling you to focus on providing care.

Add in ResMed's service, warranty and maintenance support model, and you've got a package that generates growth potential for your business and delivers freedom to your patients.

A partnership to change lives with every breath

As a global leader in respiratory care devices and applications, ResMed is a partner you can rely on. We develop cutting-edge solutions for treating and managing respiratory conditions and are committed to improving quality of life and reducing healthcare costs. With ResMed, you have a partner who offers you the technology and expertise you need to make a positive change to your patients' lives.



Self confidence on the go

The LifeChoice Activox 4L package delivers on all fronts. The device is equipped with responsive algorithms, variable settings and relevant features to meet your patients' therapeutic requirements.

The carry bags, adjustable straps and charging options provide the freedom and flexibility to meet your patients' lifestyle needs. Comfort, ease of use, portability and effective therapy: it's all contained in the LifeChoice Activox 4L package.



LifeChoice Activox 4L Package

- A** LifeChoice Activox (4L model shown)
- B** 4-way carry case
Use as a backpack, shoulder bag, waist pack or briefcase
- C** Adjustable straps
For use with 4-way carry case
- D** DC charger
For charging on the go in vehicle
- E** AC charger
For charging at home or on the go through a standard wall outlet
- F** Accessory bag
To hold power supplies, extra straps and manual
- G** Standard 210 cm long single lumen nasal cannula (sold separately)
- H** External battery (sold separately)



Technical specifications

Oxygen delivery	PULSE-WAVE – Settings 1-4
Oxygen concentration	90% ± 3%
Weight	2.2 kg (4.8 lbs)
Size (W x H x D)	230 mm x 200 mm x 112 mm
FAA approved	•

Battery runtimes

Setting	1	2	3	4
Internal battery runtime**	10hr 15min	8hr 15min	5hr	4hr
External battery runtime**	4hr 45min	3hr 45min	2hr 45min	2hr 25min
Total runtime**	15hr	12hr	7hr 45min	6hr 25min
Charge time**	4 hours internal – 2 hours external			

Product codes

Activox P4L (package)	AOX-P4L-EU	AC charger	XYC089
Carry case	123-102	EU power cord	XYC494
Shoulder strap kit	XYC414	Accessory bag	XYC340
DC charger	XYC104	Adult nasal cannula (sold separately)	XYC262
		External battery (sold separately)	200122

* In setting 1

** Hours are approximate. Runtimes are based on breaths per minute.