




Product Information Guide Fall 2020

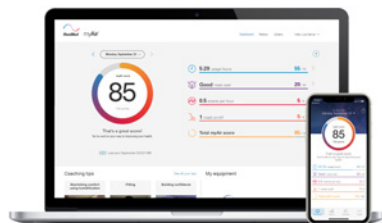

ResMed.com

1 Combined responses of those expressing a preference in an independent patient survey in the USA, UK, Germany, France & Australia, March 2017. Visit [ResMed.com/MaskBrand](https://www.resmed.com/maskbrand).

ResMed **Air** Solutions

The Science of Progress


Featuring connected devices and cloud-based software solutions, ResMed Air Solutions is a first-of-its-kind ecosystem that helps you get the most out of your sleep business. From early intervention to a robust resupply program, Air Solutions is helping patients, providers and physicians communicate and share data like never before.

Harnessing the power of digital health to transform the sleep industry – this is the Science of Progress.

ResMed ReSupply makes it easy to manage and monitor patients remotely with their preferred outreach method and helps support patient success in therapy with wellness checkups and supply orders.

GoScripts is one system in which payers and healthcare providers can send, receive, request and sign securely from anywhere completely electronically.

AirView™ is the easiest, most efficient way to manage your sleep business. This cloud-based data hub also features tools that are revolutionizing patient care.

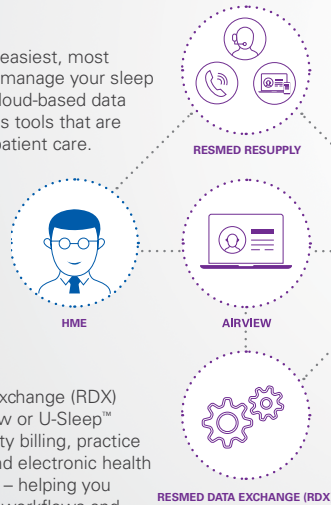

ResMed Data Exchange (RDX) connects AirView or U-Sleep™ data to third-party billing, practice management and electronic health record software – helping you streamline your workflows and improve patient care.



Air10™ devices are connected to AirView and myAir via built-in cellular connectivity – an industry first.

With myAir™, patients can track data and are coached to get more involved with sleep therapy.

"We love that the AirSense 10 has one menu for clinicians to set up the device – we're saving around 20 minutes per patient. It's also easy for patients to learn and we're getting less callbacks."

– The Medical Store
South Burlington, VT





21st century sleep therapy

From integrated cellular connectivity and humidification in Air10™ devices to the first therapy device designed specifically to treat sleep apnea in women, ResMed devices bring sleep therapy into the 21st century.

CPAP and APAP devices

DEVICE LINEUP



HCPCS codes for ResMed's CPAP/APAP devices

CODE	PRODUCT	RECOMMENDED REPLACEMENT SCHEDULE*
E0601	CPAP	1 per 5 years
E0562	Humidifier – heated	1 per 5 years
A4604	Tubing w/ integrated heating	1 per 3 months
A7037	Tubing	1 per 3 months
A7038	Filter, disposable	2 per month
A7046	Humidifier chamber – replacement	1 per 6 months
A9279	Monitoring feature/device	Contact payer

*Replacement cycles may vary depending on your patients' insurance coverage

CPAP

AirSense™ 10 series



AirSense™ 10 Elite

The AirSense 10 Elite is indicated for the treatment of obstructive sleep apnea (OSA) in patients weighing more than 66 lbs (30 kg). It is intended for home and hospital use.

- Premium CPAP therapy device with integrated humidifier, integrated cellular connectivity, patient access to myAir™ program, AutoRamp™ with sleep onset detection and EPR™ with Easy-Breathe technology
- Provides advanced event detection, including central sleep apnea and Cheyne–Stokes respiration detection, and RERA reporting as well as remote diagnostic information in AirView™ software
- Optional ClimateLineAir™ heated tubing with swivel connector and pulse oximetry monitoring

	US	Canada
AirSense 10 Elite with HumidAir™	37205	37400
AirSense 10 Elite with HumidAir and ClimateLineAir	37206	37401



AirSense™ 10 CPAP

The AirSense 10 CPAP is indicated for the treatment of obstructive sleep apnea (OSA) in patients weighing more than 66 lbs (30 kg). It is intended for home and hospital use.

- Basic CPAP therapy device with integrated humidifier, integrated cellular connectivity, patient access to myAir™ program, AutoRamp™ with sleep onset detection and EPR™ with Easy-Breathe technology
- Provides usage-compliance data, AHI and leak data reporting, and access to remote diagnostic information in AirView software
- Optional ClimateLineAir heated tubing with swivel connector

	US	Canada
AirSense 10 CPAP with HumidAir	37203	37398
AirSense 10 CPAP with HumidAir and ClimateLineAir	37204	37399

APAP

AirSense™ 10 series



AirSense™ 10 AutoSet™

The AirSense 10 AutoSet (self-adjusting) is indicated for the treatment of obstructive sleep apnea (OSA) in patients weighing more than 66 lbs (30 kg). It is intended for home and hospital use.

- Advanced AutoSet™ therapy device with integrated humidifier, integrated cellular connectivity, patient access to myAir™ program, AutoRamp™ with sleep onset detection, EPR™ with Easy-Breathe technology and AutoSet™ Response comfort setting
- Provides advanced event detection, including central sleep apnea and Cheyne–Stokes respiration detection, and RERA reporting, as well as remote diagnostic information in AirView™ software
- Optional ClimateLineAir™ heated tubing with swivel connector and pulse oximetry monitoring

	US	Canada
AirSense 10 AutoSet with HumidAir™	37207	37402
AirSense 10 AutoSet with HumidAir and ClimateLineAir	37208	37403

Remind your non-Medicare patients to check with their insurance companies for reimbursement eligibility. Replacement schedule is based on Centers for Medicare and Medicaid Services (CMS) Local Coverage Determination for positive airway pressure (PAP) devices for the treatment of obstructive sleep apnea.



AirSense™ 10 AutoSet™ for Her

The AirSense 10 AutoSet for Her (self-adjusting) is indicated for the treatment of obstructive sleep apnea (OSA) in patients (female patients with mild to moderate OSA when using AutoSet for Her treatment mode) weighing more than 66 lbs (30 kg). It is intended for home and hospital use.

- Advanced AutoSet therapy device with an FDA-cleared, dedicated AutoSet algorithm for female patients with mild to moderate OSA
- Includes integrated humidifier, integrated cellular connectivity, patient access to myAir™ program, AutoRamp™ with sleep onset detection, EPR™ with Easy-Breathe technology and AutoSet™ Response comfort setting
- Provides advanced event detection, including central sleep apnea (CSA) and Cheyne–Stokes respiration detection, RERA reporting and access to remote diagnostic information in AirView software
- Optional ClimateLineAir™ heated tubing with swivel connector and pulse oximetry monitoring

	US	Canada
AirSense 10 AutoSet for Her with HumidAir	37209	37404
AirSense 10 AutoSet for Her with HumidAir and ClimateLineAir	37210	37405

APAP

AirMini™

AirMini

AirMini is a self-adjusting system indicated for the treatment of obstructive sleep apnea (OSA) in patients (female patients with mild to moderate OSA when using AutoSet for Her treatment mode) weighing more than 66 lb (30 kg). Its small size makes it an ideal travel PAP.

- ActiveAir technology regulates flow while still increasing therapy pressure, allowing the system to be smaller, quieter and more efficient without compromising on therapy delivery or pressure swings.
- AirMini features HumidX™ and HumidX™ Plus, a waterless humidification* system providing comfortable therapy without the need for distilled water.




AirMini	38113
AirMini mask pack P10 (mask included)	38824
AirMini setup pack N20 (mask not included)	38823
AirMini setup pack F20/F30 (mask not included)	38825




HumidX HumidX Plus
38809 (3 pk) **38812** (3 pk)
38810 (6 pk) **38813** (6 pk)
38808 (50 pk) **38811** (50 pk)




38840 AirMini travel bag
38841 AirMini travel case

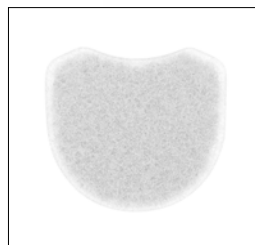

AirMini filter
38814 (1 pk)
38815 (2 pk)
38816 (12 pk)
38837 (50 pk)

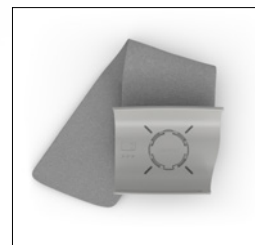

38842 AirMini mount system

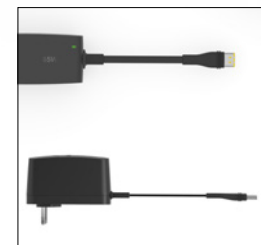

38839 AirMini DC/DC converter
38828 AirMini 20W AC adapter (US only)

* Not available with full face mask.

Accessories

Air10™ humidification and filters

HumidAir™ humidifier: Integrated humidifier for Air10 devices

Humidification is delivered in one easy-to-use system. The humidifier is intended for single patient use in the home environment and reuse in a hospital/institutional environment. Enhanced Auto option with default temperature and humidity settings and pre-set at levels designed to deliver optimal comfort.



37299 HumidAir standard tub
37479 HumidAir cleanable tub



37296 ClimateLineAir heated tubing
37357 ClimateLineAir Oxy heated tubing
36810 SlimLine™ tubing
14994 6' standard tubing
14987 6'6" standard tubing
14999 9' standard tubing (AirCurve only)
14986 6'6" autoclavable standard tubing



37335 HumidAir side cover – gray (AutoSet for Her, S, VAuto, ST, ASV)
37303 HumidAir side cover – charcoal (CPAP, Elite, AutoSet, ST-A)



Filters – standard
36850 (1 pk) **36851** (2 pk)
36852 (12 pk) **36853** (50 pk)
Filters – hypoallergenic
36855 (1 pk) **36856** (2 pk)
36857 (12 pk) **36858** (50 pk)
Filter cover
19768 Gray (AutoSet for Her, AirCurve)
19729 Charcoal (CPAP, Elite, AutoSet ST-A)

Accessories

S9™ humidification and filters

H5i™

- Integrates seamlessly with the S9
- Sleek and stylish design
- Easy-to-use operation and data management display
- Climate Control system automatically adjusts to environmental conditions for maximum comfort



36900 H5i + standard chamber
36803 H5i water chamber (standard)



36901 H5i + cleanable chamber
36800 H5i water chamber (cleanable)



36995 ClimateLine™ tubing
36996 ClimateLine MAX Oxy tubing



Filters – standard
36850 (1 pk) **36851** (2 pk)
36852 (12 pk) **36853** (50 pk)
Filters – hypoallergenic
36855 (1 pk) **36856** (2 pk)
36857 (12 pk) **36858** (50 pk)

Accessories

Power

ResMed Power Station (RPS) II kits

- All-in-one pack that provides power to ResMed ventilators and sleep therapy devices
- Provides up to 13 hours of run time on a single charge
- Weighs only 2 lbs and can be easily transported in the RPS II carry bag
- Ideal for travel or when additional mobility and freedom are required

RPS II kit for Stellar™	249002
RPS II kit for S9	249004
RPS II kit for Air10	249009
RPS II	24923
RPS II DC cable for Air10	37343
RPS II adapter for Air10	37342



37474 Air10 Power Supply Unit (90W) II (not compatible with RPS II or STA)
37344 Air10 power supply unit (90W) (RPS II and STA Only)



37297 DC/DC Converter (Air10)



36920 30W PSU (S9)
36821 90W PSU (S9)



36970 DC/DC Converter (S9)

Accessories

Faceplate and travel bag



Air10 Faceplate

- All three faceplates can be used on any AirSense 10 or AirCurve 10 device

Air10 Silver Faceplate	19736
Air10 White w/gray design Faceplate	19735
Air10 Black Faceplate	19734



Air10 travel Bag

Travel bag	37304
Travel bag for Her	37305

Accessories

Data management and oximetry



Air10 SD card
protective folder

37329 with card (10 pk)

37330 without card (10 pk)

S9 SD card protective folder

36840 with card (10 pk)

36845 without card (10 pk)

Side panel rubber SD card door

19766 Gray (AutoSet for Her,
AirCurve)

19728 Charcoal (CPAP, Elite,
AutoSet)



36931 SD card reader



37302 Air10 oximetry module

370004 Air10 oximetry

Complete Kit: Includes
oximetry module, XPOD,
XPOD clip, reusable soft
sensor (medium)



36940 S9 oximetry adapter

36950 S9 USB adapter

369100 S9 oximetry complete
kit: Includes oximetry module,
XPOD, XPOD clip, reusable
soft sensor (medium)



36944 Air10 and S9 oximetry
accessory kit

22374 XPOD LP oximeter

22371 XPOD LP clip

707560 Reusable soft sensor
(medium) – 3 m cable

1431002 Three-meter
disposable oximeter flex
sensor with single-use tapes
(25 pk)

Limited warranty information

For ResMed products

ResMed Ltd (hereafter 'ResMed') warrants that your ResMed product shall be free from defects in material and workmanship from the date of purchase for the period specified below.

Product	Warranty Period
AirTouch UltraSoft cushion	7 days
- Mask systems (including mask frame, cushion, headgear and tubing) – excluding single-use devices - Accessories – excluding single-use devices - Flex-type finger pulse sensors - Humidifier water tubs	- ClimateLine™ tubes - ClimateLineAir™ tubes - Standard tubes - SlimLine™ tubes 90 days
Batteries for use in ResMed internal and external battery systems	6 months
- Clip-type finger pulse sensors - CPAP and bilevel device data modules - Oximeters and CPAP and bilevel device oximeter adapters - Humidifier cleanable water tubs - Titration control devices	1 year
- CPAP, bilevel and ventilation devices (including external power supply units) - Humidifiers - Battery accessories - Portable diagnostic/screening devices	2 years

This warranty is only available to the initial consumer. It is not transferable.

If the product fails under conditions of normal use, ResMed will repair or replace, at its option, the defective product or any of its components.

This Limited Warranty does not cover: a) any damage caused as a result of improper use, abuse, modification or alteration of the product; b) repairs carried out by any service organization that has not been expressly authorized by ResMed to perform such repairs; c) any damage or contamination due to cigarette, pipe, cigar or other smoke; and d) any damage caused by exposure to ozone, activated oxygen or other gases.

Warranty is void on product sold, or resold, outside the region of original purchase.

Warranty claims on defective product must be made by the initial consumer at the point of purchase.

This warranty replaces all other expressed or implied warranties, including any implied warranty of merchantability or fitness for a particular purpose. Some regions or states do not allow limitations on how long an implied warranty lasts, so the above limitation may not apply to you.

ResMed shall not be responsible for any incidental or consequential damages claimed to have resulted from the sale, installation or use of any ResMed product. Some regions or states do not allow the exclusion or limitation of incidental or consequential damages, so the above limitation may not apply to you.

This warranty gives you specific legal rights, and you may also have other rights which vary from region to region. For further information on your warranty rights, contact your local ResMed dealer or ResMed office.

Bilevel devices

DEVICE LINEUP



ResMed E0470 and E0471 devices

E0470-Bilevel without a backup rate:

- AirCurve™ 10 VAuto
- AirCurve™ 10 S
- VPAP™ COPD

E0471-Bilevel with a backup rate:

- AirCurve 10 ST
- AirCurve 10 ASV
- AirCurve 10 ST-A
- Stellar 150

HCPCS codes for ResMed's bilevel devices

CODE	PRODUCT	REPLACEMENT SCHEDULE*
E0470	Bilevel	1 per 5 years
E0471	Bilevel w/ backup rate – non-invasive	1 per 5 years
E0472	Bilevel with backup rate used with invasive interface	1 per 5 years
E0562	Humidifier – heated	1 per 5 years
A4604	Tubing w/ integrated heating	1 per 3 months
A7037	Tubing	1 per 3 months
A7038	Filter, disposable	2 per month
A7046	Humidifier chamber – replacement	1 per 6 months
A9279	Monitoring feature/device	Contact payer

* Replacement schedule is based on Centers for Medicare and Medicaid Services (CMS) positive airway pressure and respiratory assist device Local Coverage Determinations. Commercial payer policies vary.

Bilevel

AirCurve™ 10 series



AirCurve™ 10 S

The AirCurve 10 S is indicated for the treatment of obstructive sleep apnea (OSA) in patients weighing more than 66 lbs (30 kg). It is intended for home and hospital use.

- Fixed-pressure bilevel device that is ideal for non-compliant CPAP patients
- Features integrated humidifier, integrated cellular connectivity, patient access to myAir™ program, Easy-Breathe waveform, Vsync leak compensation, adjustable TiControl™, trigger and cycle sensitivity
- Provides advanced event detection and reporting, as well as access to remote diagnostic information in AirView™ software
- Optional ClimateLineAir™ heated tubing with swivel connector and pulse oximetry monitoring

	US	Canada
AirCurve 10 S with HumidAir™	37213	37408
AirCurve 10 S with HumidAir and ClimateLineAir	37214	



AirCurve™ 10 VAuto

The AirCurve 10 VAuto is indicated for the treatment of obstructive sleep apnea (OSA) in patients weighing more than 66 lbs (30 kg). It is intended for home and hospital use.

- Auto-adjusting bilevel that uses the AutoSet™ algorithm and natural comfort of the Easy-Breathe waveform to treat non-compliant CPAP patients
- Features integrated humidifier, integrated cellular connectivity, patient access to myAir™ program and a variety of technologies that improve therapy comfort
- Provides advanced event detection and reporting, as well as access to remote diagnostic information in AirView software
- Optional ClimateLineAir heated tubing with swivel connector and pulse oximetry monitoring

	US	Canada
AirCurve 10 VAuto with HumidAir	37211	37406
AirCurve 10 VAuto with HumidAir and ClimateLineAir	37212	

Bilevel

AirCurve™ 10 series



AirCurve™ 10 ASV

The AirCurve 10 ASV is indicated for the treatment of patients weighing more than 66 lbs (30 kg) with obstructive sleep apnea (OSA), central and/or mixed apneas, or periodic breathing. It is intended for home and hospital use*

- Adaptive servo-ventilator targeting patients' recent minute ventilation
- Features integrated humidifier, integrated cellular connectivity (including access to remote diagnostic information in AirView™ software) and patient access to myAir™ program, Easy-Breathe waveform and Vsync leak compensation
- Choice of fixed EPAP in ASV mode or auto-adjusting EPAP in ASVAuto mode
- Optional ClimateLineAir™ tubing with swivel connector and pulse oximetry monitoring

	US	Canada
AirCurve 10 ASV with HumidAir™	37215	37410
AirCurve 10 ASV with HumidAir and ClimateLineAir	37216	



AirCurve™ 10 ST

The AirCurve 10 ST is indicated for the treatment of obstructive sleep apnea (OSA) in patients weighing more than 66 lbs (30 kg). It is intended for home and hospital use.

- Bilevel device with backup rate that provides patient-device synchrony so patients remain comfortable and well-ventilated, reducing the work of breathing
- Ideal choice for reliable, non-invasive ventilation, when alarms are not required
- Provides event detection and reporting, integrated cellular connectivity and patient access to myAir™ program, as well as access to remote diagnostic information in AirView™ software
- Optional ClimateLineAir™ tubing with swivel connector and pulse oximetry monitoring

	US	Canada
AirCurve 10 ST with HumidAir	37306	37412
AirCurve 10 ST with HumidAir and ClimateLineAir	37307	

*ASV therapy is contraindicated in patients with chronic, symptomatic heart failure (NYHA 2-4) with reduced left ventricular ejection fraction (LVEF ≤ 45%) and moderate to severe predominant central sleep apnea.

Bilevel

Non-invasive ventilation (NIV)



VPAP™ COPD

The VPAP COPD is indicated to provide non-invasive ventilation for patients weighing more than 30 lbs (13 kg) with respiratory insufficiency such as that associated with hypercapnic chronic obstructive pulmonary disease (COPD) or obstructive sleep apnea (OSA). It is intended for home and hospital use.

- Customized therapy settings make it easy to initiate therapy
- Integrated oxygen entrainment support provides oxygen and warmed, humidified air in one line
- Configurable audio and visual alarms provide added security
- Custom therapy settings optimized for COPD

S9 VPAP COPD

US Canada
36034 36034



AirCurve™ 10 ST-A

The AirCurve 10 ST-A is indicated to provide non-invasive ventilation for patients weighing more than 30 lbs (13 kg) with respiratory insufficiency or obstructive sleep apnea (OSA). It is intended for home and hospital use.

- Homecare NIV device optimized to treat patients with respiratory insufficiency
- Features iVAPS™ mode, which is designed to maintain a target alveolar ventilation by automatically adjusting pressure support. iVAPS mode is indicated for patients weighing more than 66 lbs (30 kg).
- Features integrated humidifier, integrated cellular connectivity (including access to remote diagnostic information in AirView™ software) and patient access to myAir™ program

	US	Canada
AirCurve 10 ST-A with HumidAir™	28210	28219
AirCurve 10 ST-A with HumidAir and ClimateLineAir	28211	
AirCurve 10 ST-A with optional Auto-EPAP		28224

Bilevel

Stellar 150

Stellar™ 150

ResMed Stellar™ 150 supports invasive and non-invasive ventilation for a range of adult and pediatric patients (66 lb / 30 kg and above). It offers reliable ventilation in one easy, accessible device. Stellar features an internal battery with up to 2 hours of capacity. An optional external battery (the ResMed Power Station II) provides up to 8 hours of additional capacity. Stellar 150 includes iVAPS,* which targets and maintains alveolar ventilation by automatically adjusting pressure support, even as the patient's respiratory needs change.



	US	Canada
Stellar 150 with iVAPS	24161	24161
Stellar 150 with iVAPS and optional AutoEPAP		24171

* iVAPS is indicated for patients weighing more than 66 lbs (30 kg)

Accessories

Stellar



24161 Stellar 150 v1.51




24927 H4i™ heated humidifier for Stellar




14922 Standard tube
36810 SlimLine™ tube



24977 Mobility bag

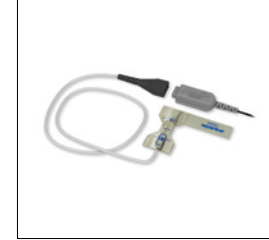

22305 Oximeter cable
1431002 One sensor & 25 tapes

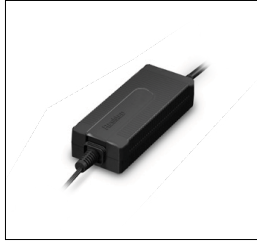

24922 Stellar DC/DC converter



24991 ResMed leak valve
24976 ResMed leak port



27903 Remote Alarm*
27904 Remote Alarm cable 2m
27905 Remote Alarm cable 5m
27906 Remote Alarm cable 10m
27907 Remote Alarm cable 20m
27912 Nurse Call cable 30m

* Cables not included

Ventilators



HCPCS codes for ResMed's ventilation devices

CODE	PRODUCT
------	---------

E0465	Home ventilator: Invasive
-------	---------------------------

E0466	Home ventilator: Non-invasive
-------	-------------------------------

Ventilators

Astral series



Astral 100

ResMed Astral™ is a portable life support ventilator offering a range of treatment options for adult and pediatric patients in the home and healthcare environments. It is indicated for patients more than 11 lbs (5 kg).*

- Through the ResMed Connectivity Module (RCM), Astral can connect to AirView™, a ResMed patient monitoring tool, giving you access to key patient data¹
- Offers an impressive 8 hours of backup internal battery power² for peace of mind
- With Big Button technology and 5 simple steps to setup, Astral is easy to use for clinicians and caregivers
- Responsive algorithms deliver adaptive, personalized ventilation as the disease progresses

* iVAPS is indicated for patients weighing more than 66 lbs (30 kg)

¹ When Astral is paired with the ResMed Connectivity Module (RCM), wireless connectivity is enabled, allowing key patient therapy data to be transmitted directly from Astral to the ResMed secure, cloud-based management system, AirView.

² An internal battery provides up to 8 hours of extra backup power if needed. A maximum of two external batteries have the capacity to deliver up to 16 hours of continuous, portable power when connected to Astral.

Astral 100

US Canada
27001 27001



Astral 150

ResMed Astral™ is a portable life support ventilator offering a range of treatment options for adult and pediatric patients in the home and healthcare environments. It is indicated for patients more than 11 lbs (5 kg).*

- Comes with all the standard features of Astral 100, plus FiO₂ monitoring, double limb circuit capabilities and additional program setting capabilities
- Through the ResMed Connectivity Module (RCM), Astral can connect to AirView™, a ResMed patient monitoring tool, giving you access to key patient data¹
- Offers an impressive 8 hours of backup internal battery power² for peace of mind
- With Big Button technology and 5 simple steps to setup, Astral is easy to utilize for clinicians and caregivers
- Responsive algorithms deliver adaptive, personalized ventilation as the disease progresses

Astral 150

US Canada
27003 27003

Accessories

Astral



27919 Astral external battery



27949 Astral SlimFit bag



27201 ResMed Connectivity Module



279001 Device trolley



27955 Mouthpiece ventilation circuit support arm



27915 Table stand



27917 Mobility bag

Portable Oxygen

Mobi

ResMed Mobi™ is indicated for patients diagnosed as requiring supplemental oxygen. It is a portable oxygen concentrator that can be used in the home, institution or travel environment.

- Pulse-Wave™ technology efficiently delivers the selected dose of oxygen that is congruent with the patient’s breath rate
- Ideal balance of oxygen, weight and battery life
- Designed to go wherever your patients go at only 5.5 lbs
- One of the longest-lasting batteries* of any portable oxygen concentrator at 9 hours on setting 2
- External battery option provides up to 3 hours of additional battery life on setting 2
- Simple, easy-to-use design
- Clear alarms and on-screen notifications for low oxygen purity, low battery, no breath detected and high temperature
- Commences automatic dosing of oxygen if breathing cannot be detected
- Stylish carry bag conceals tubing and batteries



Portable oxygen concentrator	E1392
Mobi set up pack	35001



35051 Mobi carry bag†



35050 Mobi backpack



35053 AC/DC charger†
XYC104 DC/DC charger†
35052 Mobi external battery



35055 Mobi mounting straps



XYC262 Cannula, adult nasal

* Based on published reports of battery performance
† Item can be purchased individually




Business efficiencies

Since 2014, ResMed Air Solutions has been delivering on a promise to help streamline sleep therapy for both patients and providers. In 2016, the ResMed Connectivity Module (RCM) was introduced to provide remote access to key ventilation data within our cloud-based patient management system, AirView.

"With the Remote Assist feature in AirView, I can now easily resolve 99.9% of patient callbacks from my office. It's a huge time-saver."

— Quinlan's Pharmacy and Medical Supply
Wayland, NY

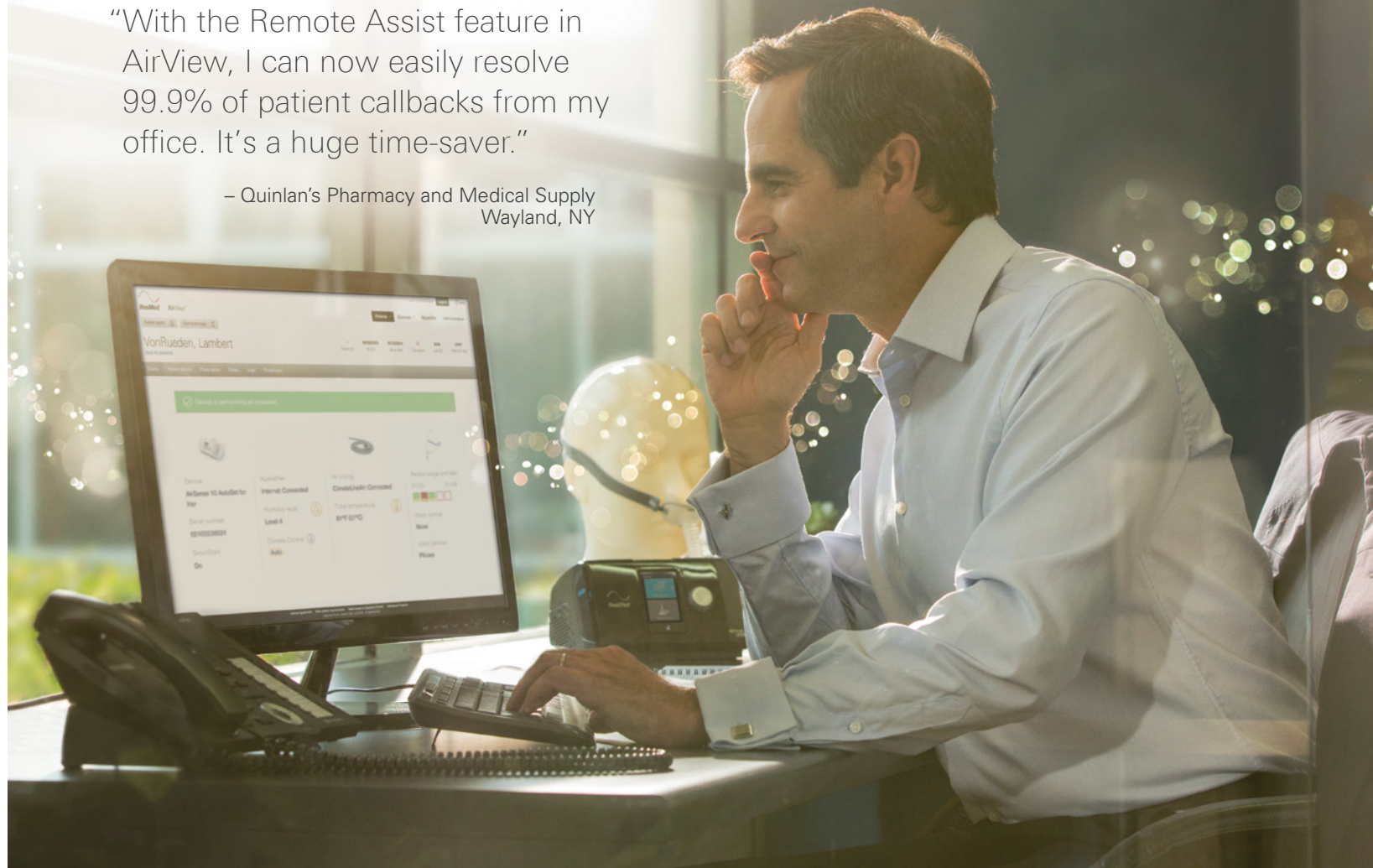

AirView™

Streamlined patient management

With AirView, you can manage patients, simplify workflows and collaborate more efficiently.

Remote Assist


The Remote Assist feature provides access to information from your patients' devices, so you can answer questions and troubleshoot issues without leaving your desk.

Device status

By displaying a green bar when performing properly and a red bar when there's something wrong, you automatically know whether the issue can be resolved over the phone.

Ambient humidity

Quickly view the patient's room humidity level and device humidity settings to recommend changes.

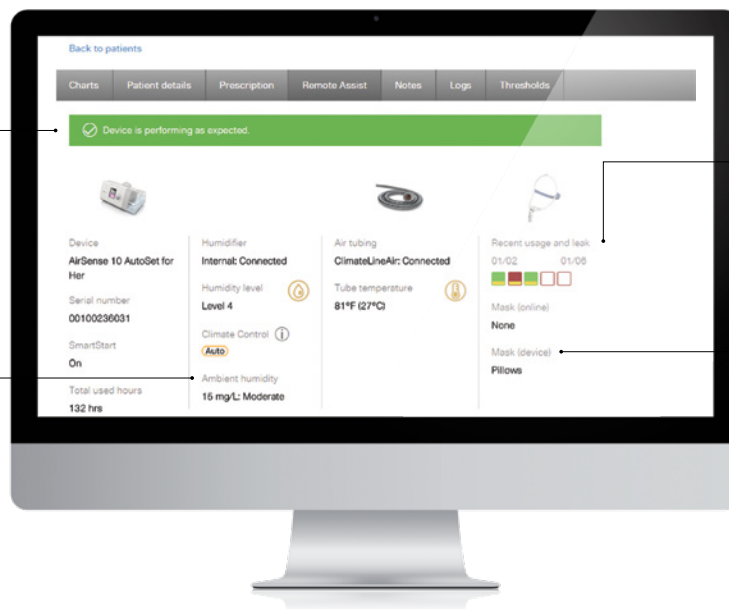

Last 5 days of usage/leak data

Understand how well the patient has been sleeping and how well their mask is fitting.

Mask type

Troubleshoot comfort issues or recommend a different mask type based on patient feedback.

With Remote Assist, you get the complete picture and can help resolve patient callbacks in minutes.

AirView™

Streamlined patient management

Daily compliance status icons (located in “Wireless patients” dashboard)

Compliant	Non-compliant		Compliant	Non-compliant	
					High AHI, high leak
		High leak			AHI, leak data not supported
		Leak data not supported			Device not used
		High AHI			Monitoring – no data received yet
		AHI data not supported			Not monitoring – no data available
		High AHI, leak data not supported			Device has fault – no data available
		High leak, AHI data not supported			

Icon legend: Top – AHI / Bottom – Leak: below threshold (green), above threshold (yellow), data not supported (white)
Overall color – Usage: above threshold (green) – Compliant, below threshold (red) – Non-compliant



URL: AirView.ResMed.com
Training and education: ResMedEducation.com

My username

ResMed tech support

Phone 800-424-0737 option 2, then option 1

Hours 5:30am - 5:00pm PT

Email airviewsupport@resmed.com

AirView for Ventilation

Streamlined patient management

Improve business efficiencies and take patient care to the next level

ResMed AirView™ for Ventilation* has been completely reimagined with new features – like management by exception and interactive clinical reporting – to help transform the remote management of ventilation patients. Combining clinical depth with clinical efficiency, AirView for Ventilation offers a simplified user experience and features designed to help you:


- **Confidently manage more respiratory patients with greater efficiency.** Management by exception allows you to spend less time identifying which patients require immediate attention. And built-in best practices will allow you to standardize care by giving staff an easy, consistent way to monitor and manage patients.
- **Enhance quality of care** by being alerted to the early signs of disease progression so you can proactively intervene.
- **Strengthen the relationship between physician and HME.** With the tools to efficiently and proactively deliver advanced care, you can differentiate your business, strengthen your relationship with referring providers and confidently take on additional respiratory patients.



*Available with ResMed Astral & Stellar devices through the ResMed Connectivity Module (RCM)

AirView for Ventilation

Streamlined patient management

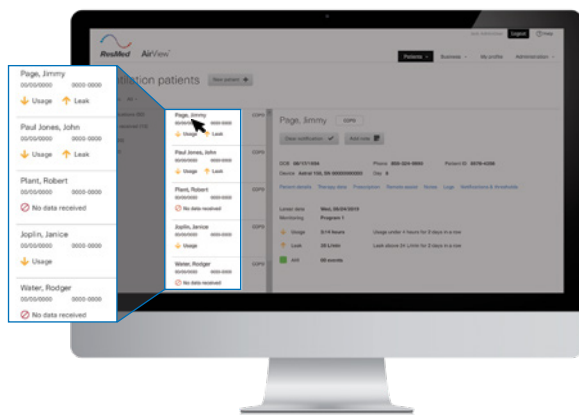

Management by exception

Using built-in best practices, management by exception improves the efficiency of the care team by automatically detecting and notifying clinicians of patient issues. With adjustable rules and thresholds that can be set by patient type, it acts as a member of the care team, bringing to the forefront patients who may require attention. Management by exception allows HMEs to handle higher patient volumes without sacrificing quality of care and intervene at the earliest signs of trouble to help address disease progression and prevent hospitalizations.




Interactive reporting

Interactive and dynamic reporting provides a holistic view of patient therapy with on-screen interactive tools designed to help visualize and detect variations in therapy trends across all programs and all settings. With the ability to zoom in and out to view clinical detail and hover over graphs for precise therapy values, interactive reporting is designed to help HMEs quickly identify variations in therapy trends.

myAir™

Therapy self-management

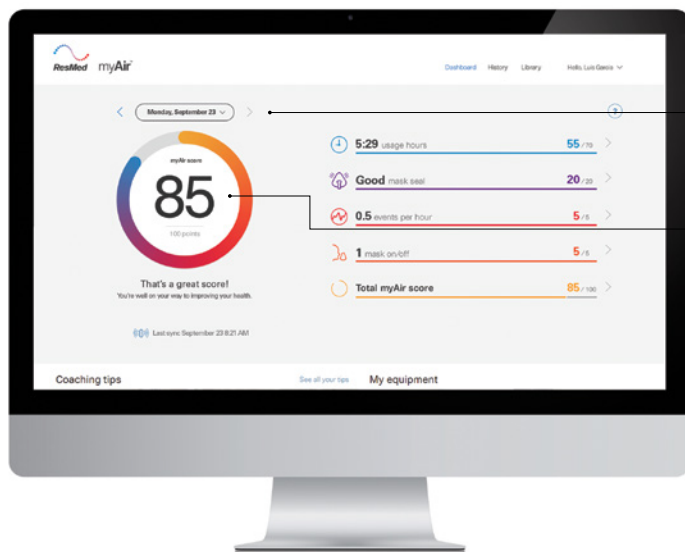

Sleep days

myAir scores are available for the past 30 nights.

myAir Score

The myAir score provides a snapshot of how well therapy went the previous night based on four sleep factors: usage hours (70 points); mask leak, events per hour and mask on/off events (30 points).


ResMed myAir allows your patients to track their nightly sleep data. It also includes a personalized coaching program to encourage, motivate and guide them to stay engaged with therapy and remain compliant over the long term. It provides instructional videos on a variety of topics, including mask fitting.

Things to remember:

- We've made it easy for you to help your patients register for myAir. Now, you can send an email directly from AirView, U-Sleep or Brightree™ inviting your patients to join myAir.
- Serial Number, Device Number, and Mask are required to register for myAir
- Patients don't have to be in AirView to register – any patient with an AirSense 10 or AirCurve 10 device is eligible
- Data is available following the first night of use after registering for myAir
- No downloads or pairing required

myAir™

Therapy self-management

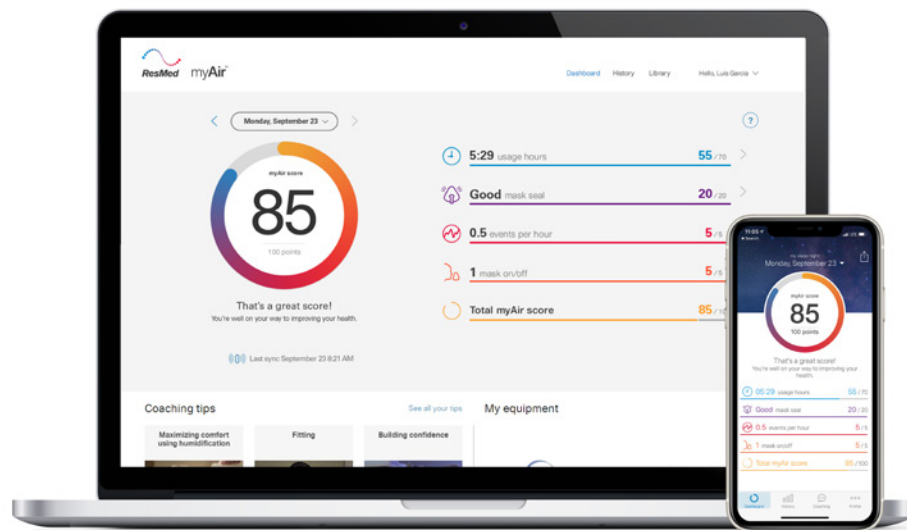

Personalized coaching

Patients can choose to receive emails and/or text messages with tips on how to improve their therapy experience. A video link is often included to help them stay on track with their:

- Usage hours
- Events per hour (AHI)
- Mask leak

Patients can register at **myAir.ResMed.com** or you can send them an invitation directly from AirView, U-Sleep or Brightree.

Patients can download the myAir for Air10™ by ResMed app* for Apple® iPhone® on the App Store® or for Android™ on the Google Play™ store. Or patients can register online at **ResMed.com/myAir**.



*The myAir app is available in English only and in the US, Canada, Australia and New Zealand only.

Resupply solutions

The importance of a resupply program

An optimized resupply program can help clinical outcomes and increase revenue. Many sleep apnea patients may not realize the importance of regular mask and supply replacement to ensure a comfortable fit and maximum effectiveness. When CPAP equipment is replaced regularly and working optimally, your patients can tell the difference.

How does a resupply program affect compliance?

60% of patients who failed on therapy cited comfort as the primary reason¹

More than 80% of patients who frequently resupply stated they use their therapy 7 nights a week¹


What is the cost of not having a resupply program?



More than 65% of patients are never given the choice to enroll in a resupply program¹

Nearly 70% of those patients would have liked the choice to enroll in a resupply program¹

Resupply materials available

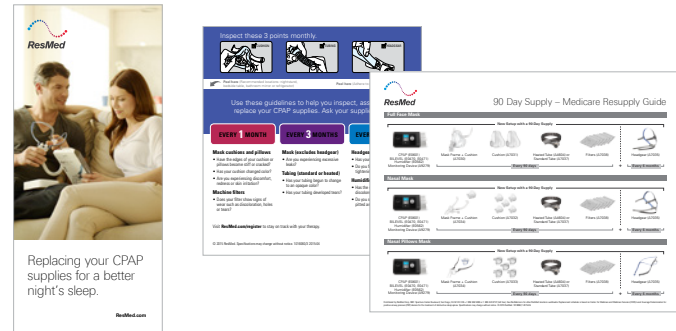

Stay on Track – Replenishment Patient Brochure (PN 1016061)

Education on the importance of replacement and what to look for

Stay on Track – Replenishment Decal (PN 1016060)

A 2-part decal patients can use as a replacement reminder

90-day Supply Guide (PN 1018890)

Shows which supplies a patient is eligible for at time of setup, every 90 days and every 6 months, according to Medicare guidelines


¹ Survey conducted for ResMed by third party market research firm (June 2013 & April 2015). Frequent resuppliers are defined as: those who receive replacement masks or parts for their mask at least every 3 months.

ResMed ReSupply™

Automated resupply made easy

ResMed ReSupply is an easy and efficient way to manage a comprehensive resupply program, from resupply enrollment through long-term retention. The program features a user-friendly platform for equipment providers, flexible communication options and an intuitive patient portal that allows patients to complete their own wellness checkup and manage their resupply needs.


How it works:

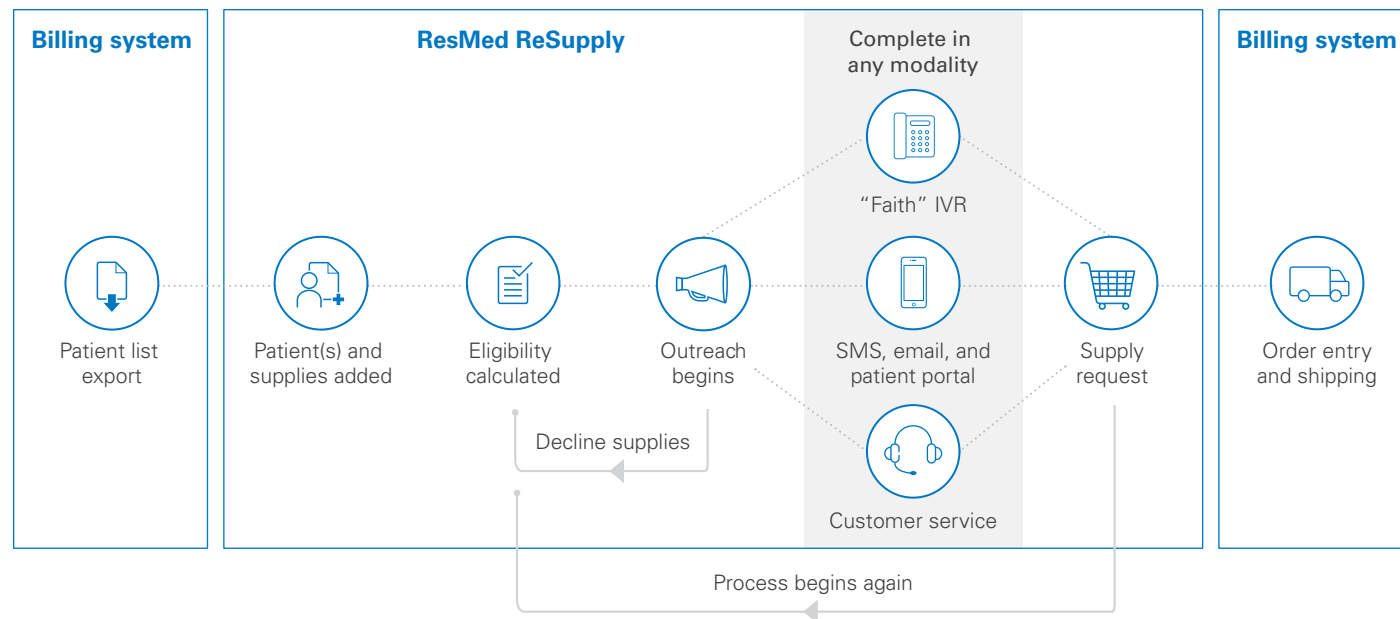

Replacement schedule is based on Centers for Medicare and Medicaid Services (CMS) Local Coverage Determination for positive airway pressure (PAP) devices for the treatment of obstructive sleep apnea but can be customized.


Inspection points

What should a patient look for?

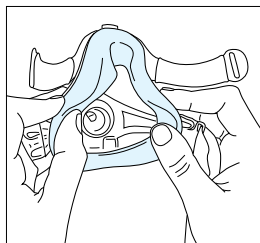


Cushions/pillows

- Deterioration over time through regular wear and tear, which can cause discomfort, compromise fit and lead to air leaks

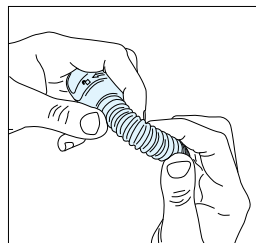


Tubing

- May develop small holes or tears, which can cause air leaks and jeopardize the delivery of prescribed therapy settings

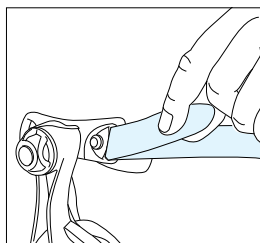

Headgear/chin straps

- May become stretched and lose their elasticity, leading to overtightening and discomfort
- Can also trap bacteria from sweat and moisture

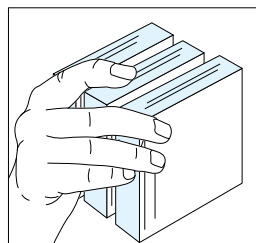

Filters

- Wear out and can become clogged, no longer protecting the patient from airborne particles, dust, pet fur and bacteria

For additional information, refer to the device's user guide.

Recommended setup and resupply schedules*

	Initial setup (1st month supplies)			Compliance follow-up (2nd and 3rd month supplies)			Quarterly outreach contact ("Every 3 months" supplies)		
Masks	Provide the patient with a 30-day supply of cushions and a 90-day supply of filters, headgear and tubing.			If patients are compliant and satisfied with their mask, then offer them two additional months of cushions.			Contact your patient for regular follow-up and any needed replacement supplies.		
	QTY	ITEM	HCPCS	QTY	ITEM	HCPCS	QTY	ITEM	HCPCS
Full face masks	1	Mask	A7030				1	Mask	A7030
	1	Headgear	A7035				1	Tubing	A7037 or A4604
	1	Tubing	A7037 or A4604				6	Filters	A7038
	6	Filters	A7038	2	Cushions	A7031	2	Cushions	A7031
Nasal masks	1	Mask	A7034				1	Mask	A7034
	1	Headgear	A7035				1	Tubing	A7037 or A4604
	1	Tubing	A7037 or A4604				6	Filters	A7038
	6	Filters	A7038				5	Cushions	A7032
	1	Cushion	A7032	4	Cushions	A7032			
Nasal pillows masks	1	Mask	A7034				1	Mask	A7034
	1	Headgear	A7035				1	Tubing	A7037 or A4604
	1	Tubing	A7037 or A4604				6	Filters	A7038
	6	Filters	A7038				5	Pillows	A7033
	1	Pillow	A7033	4	Pillows	A7033			

Replacement schedule is based on Centers for Medicare and Medicaid Services (CMS) positive airway pressure and respiratory assist device Local Coverage Determinations.

*Replacement cycles may vary depending on your patients' insurance requirements




The ResMed Online Store

Save time and money while tracking down product info, pricing and shipments with our complete e-business solution. Log in or sign up for free at **OnlineStore.ResMed.com**.

You can log onto ROS anytime, anywhere and:

- Place orders with a single click
- Find product info, real-time pricing, inventory status, part numbers and reimbursement codes
- Check your product's warranty
- Track shipments
- Access digital packing slips
- Initiate a return merchandise authorization (RMA) in just minutes
- View and pay invoices

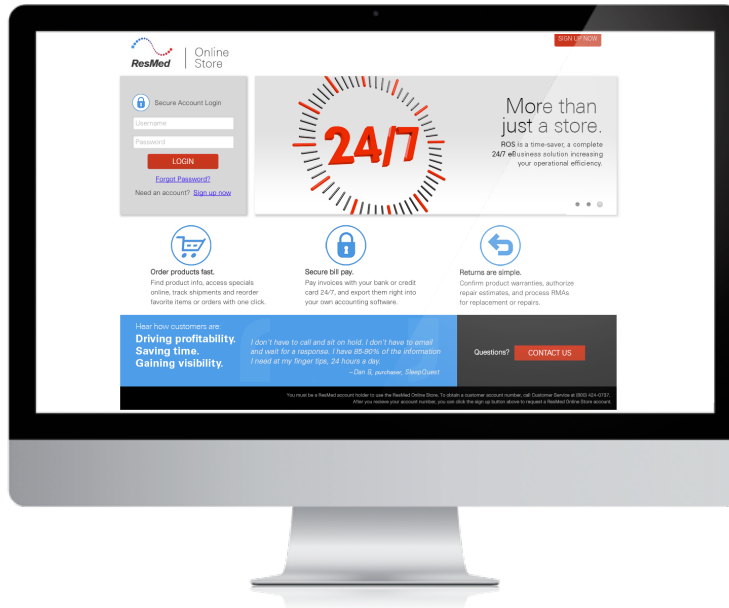

ROS username

ResMed account#



Patient comfort

We offer a wide range of effective and comfortable masks to fit all face types – men and women, adults and children.

An easier way to help patients choose the right mask

Helping your patients choose the right mask at setup is one of the most important parts of CPAP therapy success. ResMed mask categories make it easy to understand key benefits and can serve as a conversation starter between you and your patients. Every ResMed mask is designed to deliver effective therapy and a comfortable night's sleep, but each category offers something a little different.

Minimalist category

Small. Light. Simple.



AirFit™ N30



AirFit™ F30



AirFit™ P10

Freedom category

Sleep in any position.



AirFit™ F30i



AirFit™ N30i



AirFit™ P30i

Ultra Soft category

Patented memory foam offers unique comfort.



AirTouch™ F20



AirTouch™ N20

Universal Fit category

Classic designs that fit nearly every face.



AirFit™ F20



AirFit™ N20

Full face masks

MASK LINEUP



AirFit™ F30i



AirFit™ F30



AirTouch™ F20
AirTouch™ F20 for Her



AirFit™ F20
AirFit™ F20 for Her

Minimalist category

Freedom category

Ultra Soft category

Universal Fit category

HCPSC codes for ResMed full face masks* (excluding Mirage Liberty)

CODE	PRODUCT	REPLACEMENT SCHEDULE
A7030	Full face mask system	1 per 3 months
A7031	Full face mask cushion – replacement	1 per month
A7035	Headgear	1 per 6 months

HCPSC codes for ResMed Mirage Liberty*

CODE	PRODUCT	REPLACEMENT SCHEDULE
A7027	Oral/nasal mask system	1 per 3 months
A7028	Oral cushion replacement for oral/nasal mask	2 per month
A7029	Nasal pillow replacement for oral/nasal mask	2 per month
A7035	Headgear	1 per 6 months

OTHER MASKS INCLUDE:

AirFit™ F10/AirFit F10 for Her
Quattro™ Air/Quattro Air for Her
Mirage Quattro™
Mirage Liberty™
Quattro™ FX/Quattro FX for Her

For all legacy mask information, visit [ResMed.com/MaskGuide](https://www.resmed.com/MaskGuide).

* The Medicare schedule states how often accessories can be replaced if medically necessary. Check with all payers directly to verify current coverage requirements and policies.

AirFit F30i

Full face mask

Quick fitting guide

ResMed AirFit™ F30i is a full face mask with a top-of-head tube connection designed to deliver a reliable seal with the freedom to sleep in any position. AirFit F30i was preferred by 3 out of 4 CPAP users to DreamWear™ full face for its overall performance in a direct comparison study.*

* Calculations based on CPAP users in the study who had a mask preference

1 ResMed guided external clinical study of 48 evaluable current full face CPAP mask users from Oct. 24 to Nov. 20, 2018, in Sydney, AU and San Diego, CA, who trialed AirFit F30i and Phillips™ Respironics DreamWear full face mask at home.



- With one hand holding the bottom of the cushion and the other holding the conduit frame, position the mask in front of the face.



- Place the cushion under the nose and against the face. With the ResMed logo and headgear facing up, pull the headgear over the head.



- Bring the lower headgear straps under the ears and attach the magnetic clips to the frame connectors.



- Undo the fastening tabs on the upper headgear straps and pull evenly.



- Repeat with the lower headgear straps.



- Squeeze the side buttons on the elbow and detach it from the frame. Connect the air tube from the device to the elbow, then reattach it to the frame.

myAir™ Compatible with myAir™ patient engagement program when paired with most Air10™ devices.

Refer to the user guide for further instructions and complete labeling information, including indications, contraindications, warnings and precautions.

AirFit F30i

Full face mask

Systems and components

AirFit F30i

Complete mask system

(Includes mask frame, full face cushion and headgear)

Small system

Small **63330**

Standard system

Small **63331**

Medium **63332**

Wide **63333**

Large system

Medium **63334**

Frame system

(includes mask frame and full face cushion)

Small frame

Small **63355**

Medium **63356**

Small-Wide **63357**

Wide **63358**

Standard frame

Small **63359**

Medium **63360**

Small-Wide **63361**

Wide **63362**

Large frame

Small **63363**

Medium **63364**

Small-Wide **63365**

Wide **63366**

Complete mask system

63330

63331

63332

63333

63334

Elbow

Conduit frame

63368 (SML)

63369 (STD)

63370 (LGE)

Full face cushion

63350 (S)

63351 (M)

63352 (SW)

63353 (W)

Headgear

63372 (STD)

AirFit F30

Full face mask

Quick fitting guide

ResMed AirFit™ F30 is a smaller style full face mask. It's designed to sit low on – and cover less of the face to prevent top-of-the-nose red marks. Patients will feel more confident using therapy around loved ones with a simple, yet effective minimal-contact mask that gives them a modern look.

We've also included the QuietAir diffuser vent used with the F20 series.



- Place the UltraCompact cushion under the patient's nose and ensure it sits comfortably on the patient's face. With the ResMed logo on the headgear facing up, pull the headgear over the patient's head.



- Bring the lower headgear straps under the patient's ears and attach the magnetic headgear clips to the frame.



- Undo the fastening tabs on the upper headgear straps and pull evenly. Repeat with the lower headgear straps.



- If the headgear is too close to the patient's ears, adjust the headgear strap on the crown of the head.



- Connect the air tubing from the CPAP machine to the elbow.



- Reattach the elbow to the frame. The mask is now ready to use.

Compatible with AirMini™ when used with the F20/AirMini connector elbow.



myAir™ Compatible with myAir™ patient engagement program when paired with most Air10™ devices.

Refer to the user guide for further instructions and complete labeling information including indications, contraindications, warnings and precautions.

AirFit F30

Full face mask

Systems and components


AirFit F30 Complete mask system

Small **64100**

Includes frame, small cushion, elbow, soft sleeves and standard blue headgear with headgear clips

Medium **64101**

Includes frame, medium cushion, elbow, soft sleeves and standard blue headgear with headgear clips

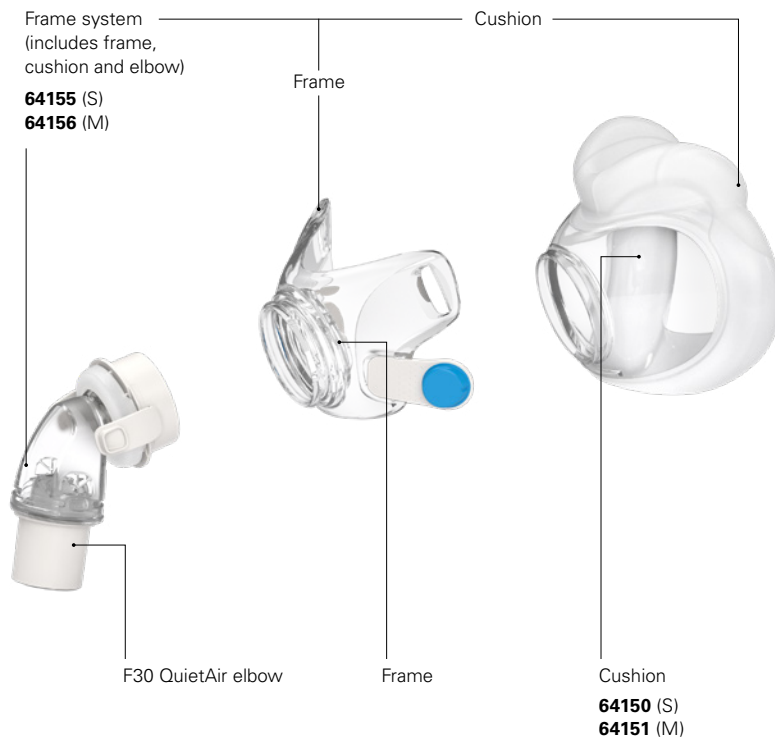

AirFit F30 headgear: Standard
64161



AirTouch F20/AirTouch F20 for Her

Full face mask

Quick fitting guide

When it comes to effective sleep therapy, comfort is key. That's why the AirTouch F20 features the ResMed UltraSoft™ memory foam cushion, a breakthrough in CPAP mask comfort technology. This innovative new cushion conforms to the unique contours of each face for a truly personalized fit. It's designed to deliver a light, breathable seal across a range of therapy pressures and reduce facial markings and potential skin irritation. Featuring a modular frame, you can conveniently switch between the AirTouch F20 UltraSoft™ and AirFit F20 InfinitySeal™ cushions until you find the best option for each patient.

We've also pioneered QuietAir technology – a diffuser vent that makes ResMed full face masks 89% quieter and 70% gentler!^{1,2}

¹ ResMed internal verification testing conducted from July–August of 2017.

² Measures as the reduction in peak velocity when compared to the standard F20 elbow.

Compatible with AirMini™ when used with the F20/AirMini connector elbow.



myAir™ Compatible with myAir™ patient engagement program when paired with most Air10™ devices.



- Start by ensuring the magnets on the lower straps are unclipped.
- If they are clipped, detach them by twisting and pulling them away from the frame.



- Pull the headgear over the patient's head making sure the logo is facing outwards.



- Position the lower straps so they both sit below the patient's ears.
- Bring each of the magnets up to meet their corresponding clips on the frame.



- Adjust the fastening tabs on the upper straps until the mask feels even and stable.
- Reattach the fastening tabs.



- Repeat step 4 with the lower straps.
- Reattach the fastening tabs.



- Connect the device's air tube to the elbow, then click the elbow into the front of the mask.
- Start the device.

Refer to the user guide for further instructions and complete labeling information including indications, contraindications, warnings and precautions.

AirTouch F20/AirTouch F20 for Her

Full face mask

Systems and components

AirTouch F20 Complete mask system

Small **63000**

Includes frame, small cushion, elbow, soft sleeves and standard blue headgear with headgear clips

Medium **63001**

Includes frame, medium cushion, elbow, soft sleeves and standard blue headgear with headgear clips

Large **63002**

Includes frame, large cushion, elbow, soft sleeves and standard blue headgear with headgear clips

AirTouch F20 for Her Complete mask system

Small **63003**

Includes frame, small cushion, elbow, soft sleeves and small lavender headgear with headgear clips

Medium **63004**

Includes frame, medium cushion, elbow, soft sleeves and small lavender headgear with headgear clips

Frame system (includes frame, soft sleeves, cushion and elbow)

63021 (S)

63022 (M)

63023 (L)

63024 for Her (S)

63025 for Her (M)

Frame

Cushion

Elbow

Frame

AirTouch F20

AirTouch F20 for Her

Cushion

63028 (S)

63029 (M)

63030 (L)



AirTouch F20 headgear with headgear clips

63470 (S)

63471 (Std)

63472 (L)



AirTouch F20 for Her headgear with headgear clips

63473 (S)



AirFit F20/AirFit F20 for Her

Full face mask

Quick fitting guide

Finding the right mask fit for a broad range of patients can be challenging. That's why we created the AirFit™ 20 series – masks that are specially-designed to fit all patients, regardless of facial profile. The AirFit F20 is designed to perform across a wide range of therapy pressures, and fit a remarkable 97% of patients in a recent international fitting study.¹⁻³

We've also pioneered QuietAir technology – a diffuser vent that makes ResMed full face masks 89% quieter and 70% gentler!^{4,5}

¹ ResMed AirFit F20 internal Australian fitting study of existing CPAP patients conducted between March–April 2016.

² ResMed AirFit F20 internal USA fitting study of existing CPAP patients conducted April 2016.

³ ResMed AirFit F20 internal EU and APAC fitting study of existing CPAP patients conducted June 2016.

⁴ ResMed internal verification testing conducted from July–August of 2017.

⁵ Measures as the reduction in peak velocity when compared to the standard F20 elbow.

Compatible with AirMini™ when used with the F20/AirMini connector elbow.



myAir™ Compatible with myAir™ patient engagement program when paired with most Air10™ devices.



- Start by ensuring the magnets on the lower straps are unclipped.
- If they are clipped, detach them by twisting and pulling them away from the frame.



- Hold the mask to the face, ensuring the top of the mask is aligned with the top of the nasal bridge.
- Pull the headgear over the head, making sure the logo is facing outward.



- Position the lower straps so they both sit below the ears.
- Bring each of the magnets up to meet their corresponding clips on the frame.



- Adjust the fastening tabs on the upper straps until the mask feels even and stable.
- Reattach the fastening tabs.



- Repeat step 4 with the lower straps.
- Reattach the fastening tabs.



- Connect the device's air tube to the elbow, then click the elbow into the front of the mask.
- Start the device.

Refer to the user guide for further instructions and complete labeling information including indications, contraindications, warnings and precautions.

AirFit F20/AirFit F20 for Her

Full face mask

Systems and components

AirFit F20 Complete mask system

Small **63400**

Includes frame, small cushion, elbow, soft sleeves and standard blue headgear with headgear clips

Medium **63401**

Includes frame, medium cushion, elbow, soft sleeves and standard blue headgear with headgear clips

Large **63402**

Includes frame, large cushion, elbow, soft sleeves and standard blue headgear with headgear clips


AirFit F20 for Her Complete mask system

Small **63403**

Includes frame, small cushion, elbow, soft sleeves and small lavender headgear with headgear clips

Medium **63404**

Includes frame, medium cushion, elbow, soft sleeves and small lavender headgear with headgear clips

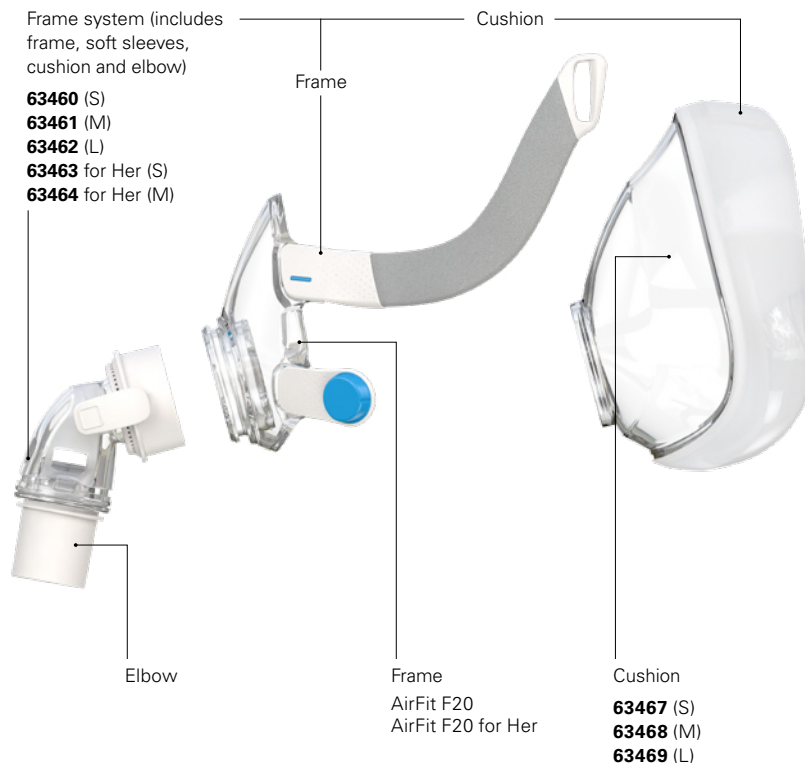

AirFit F20 headgear with headgear clips

63470 (S)
63471 (Std)
63472 (L)



AirFit F20 for Her headgear with headgear clips

63473 (S)



AirFit F10/AirFit F10 for Her

Full face mask

Systems and components

AirFit F10

Complete mask system

Small **63101**

Includes frame, small cushion, elbow with swivel, soft sleeves and standard blue headgear

Medium **63102**

Includes frame, medium cushion, elbow with swivel, soft sleeves and standard blue headgear

Large **63103**

Includes frame, large cushion, elbow with swivel, soft sleeves and standard blue headgear

AirFit F10 for Her

Complete mask system

X-Small **63139**

Includes frame, x-small cushion, elbow with swivel, soft sleeves and small pink headgear

Small **63140**

Includes frame, small cushion, elbow with swivel, soft sleeves and small pink headgear

Medium **63141**

Includes frame, medium cushion, elbow with swivel, soft sleeves and small pink headgear

Frame system (includes frame, cushion, soft sleeves and elbow with swivel)

63160 (XS)

63161 (S)

63162 (M)


63163 (L)

Frame

Cushion



Elbow with swivel

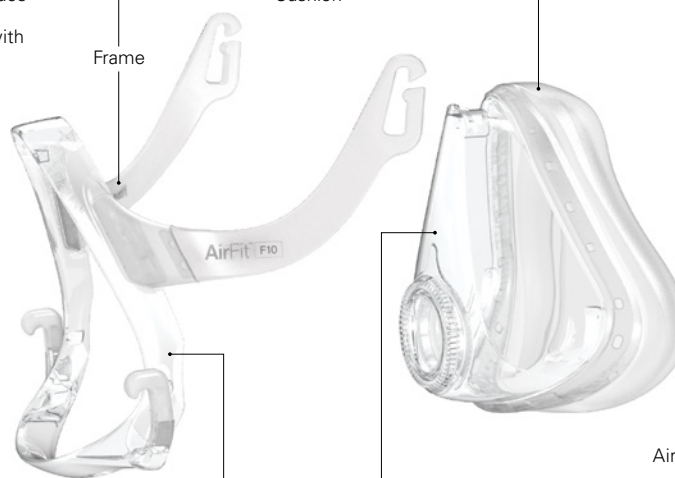

Frame

Cushion

62736 (XS)

62737 (S)

62738 (M)

62739 (L)



AirFit F10 headgear

63165 (S)

63164 (Std)

63166 (L)



AirFit F10 for Her headgear

63167 (S)



Quattro Air/Quattro Air for Her

Full face mask

Systems and components

Quattro Air Complete mask system

Small **62701**

Includes frame, small cushion,
elbow with swivel and standard
blue headgear

Medium **62702**

Includes frame, medium cushion,
elbow with swivel and standard
blue headgear

Large **62703**

Includes frame, large cushion,
elbow with swivel and standard
blue headgear

Quattro Air for Her Complete mask system

X-Small **62740**


Includes frame, x-small cushion,
elbow with swivel and small
pink headgear

Small **62741**

Includes frame, small cushion,
elbow with swivel and small
pink headgear

Medium **62742**

Includes frame, medium cushion,
elbow with swivel and small
pink headgear

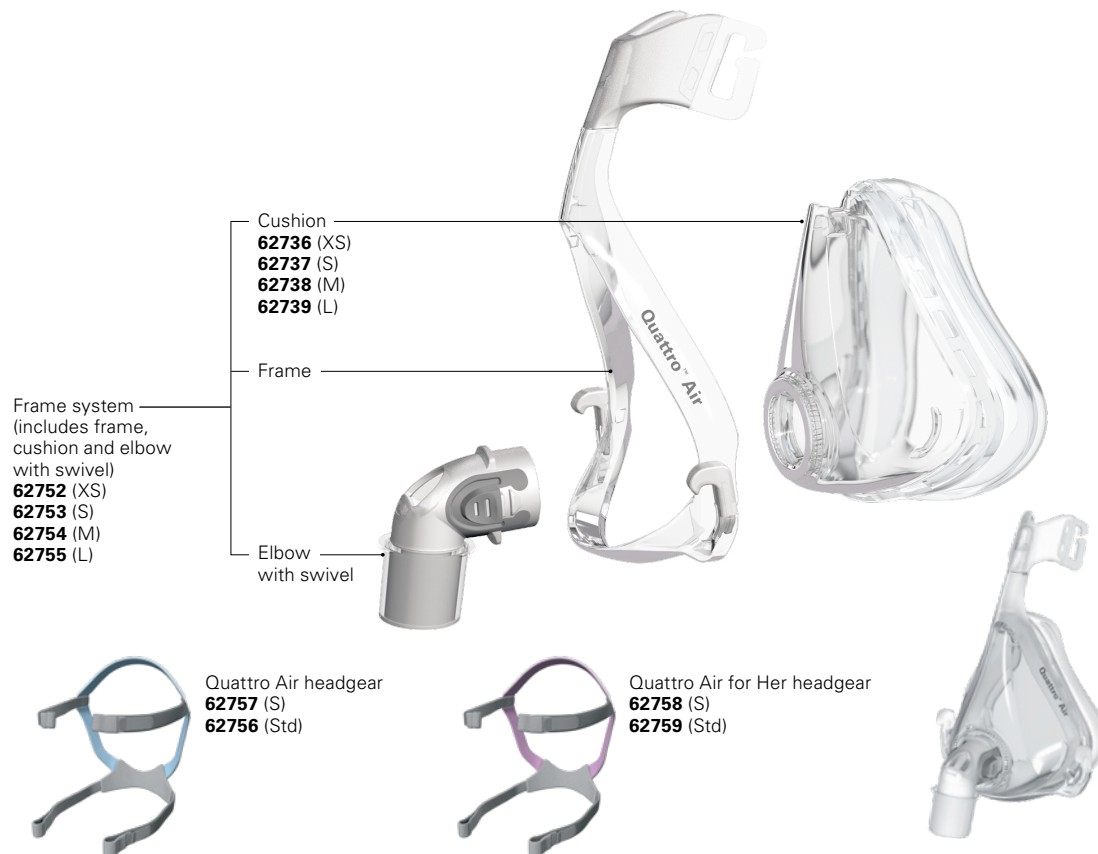

Mirage Quattro

Full face mask

Systems and components

Complete mask system¹

Extra Small	61200
Small	61201
Medium	61202
Large	61203

Frame system²

(no headgear)

Extra Small	61260
Small	61261
Medium	61262
Large	61263

Frame assembly²


(no headgear, no cushion)

Extra Small	61264
Small	61265
Medium	61266
Large	61267

Headgear

(no clips)

Small	16118
Standard	16117
Large	16119

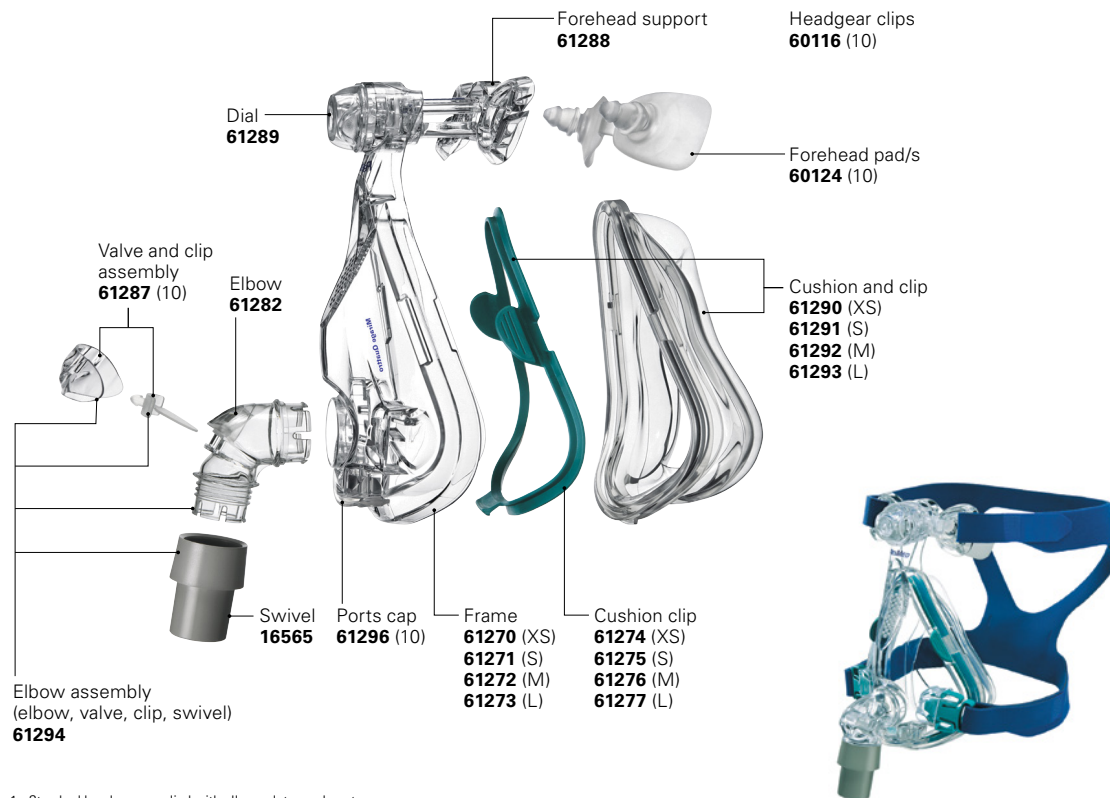

1 Standard headgear supplied with all complete mask systems

2 Includes headgear clips

Mirage Liberty

Full face mask

Systems and components

Complete mask system¹

Small **61300**
Large **61301**

Frame system²

(no headgear)

Cushion Pillows

Small	Small	61318
Small	Medium	61319
Small	Large	61320
Large	Small	61321
Large	Medium	61322
Large	Large	61323

Frame assembly²


(no headgear, no cushion,
no pillows)

Small **61325**
Large **61326**

Headgear

(no clips)

Small **61348**
Standard **61349**

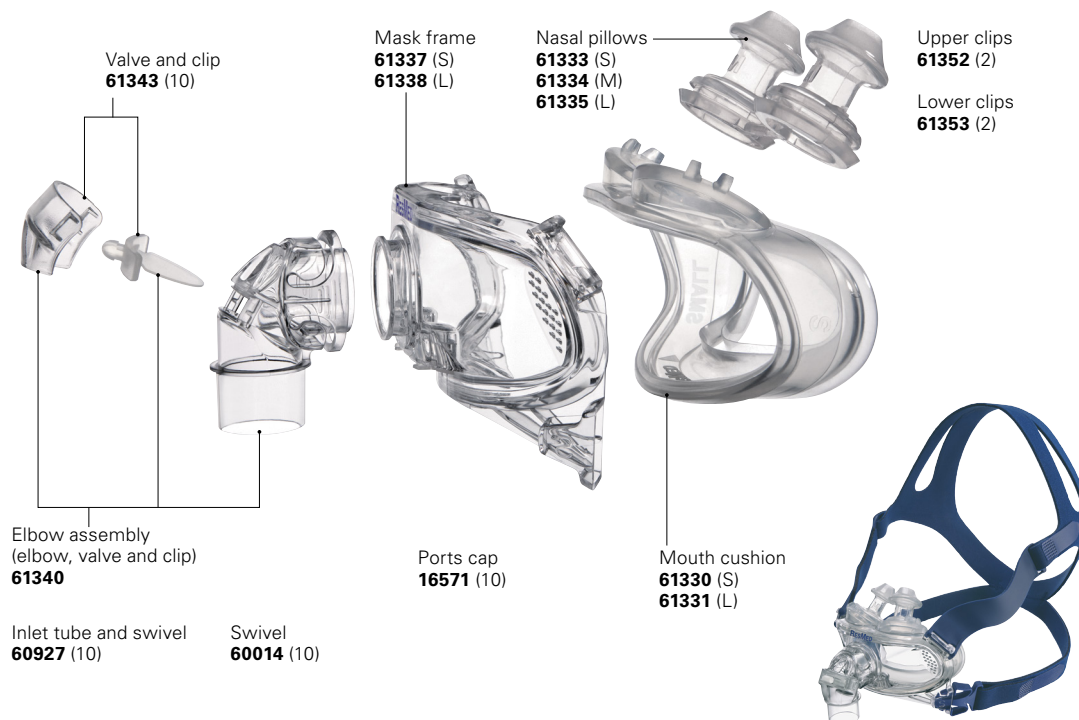

¹ Standard headgear supplied with all complete mask systems

² Includes headgear clips

Quattro FX/Quattro FX for Her

Full face mask

Systems and components

Quattro FX Complete mask system¹

Small	61700
Medium	61701
Large	61702

Frame system²

(headgear clips but no headgear)

Small	61724
Medium	61725
Large	61726

Quattro FX for Her Complete mask system



(standard headgear supplied with Quattro FX for Her complete mask systems)

Small	62501
Medium	62502

Frame system²

(headgear clips but no headgear)

Small	62513
Medium	62514

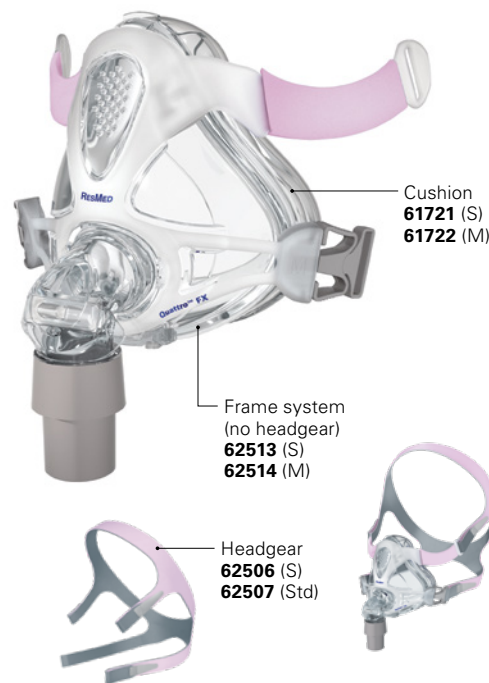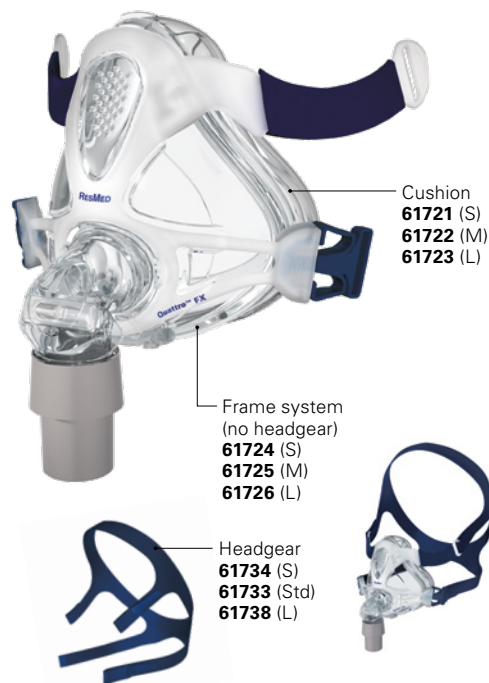

¹ Standard headgear supplied with all complete mask systems

² Includes headgear clips

Nasal masks

MASK LINEUP



AirFit™ N30

Minimalist category



AirFit™ N30i

Freedom category



AirTouch™ N20
AirTouch™ N20 for Her

Ultra Soft category



AirFit™ N20
AirFit™ N20 for Her

Universal Fit category

OTHER MASKS INCLUDE:

AirFit™ N10/AirFit N10 for Her

Mirage™ FX /Mirage FX for Her

Mirage Activa™ LT

Swift™ FX Nano/ Swift FX Nano for Her

For all legacy mask information, visit [ResMed.com/MaskGuide](https://www.resmed.com/MaskGuide).

HCPCS codes for ResMed's nasal masks*

CODE	PRODUCT	REPLACEMENT SCHEDULE
A7034	Nasal mask system	1 per 3 months
A7032	Nasal cushion replacement	2 per month
A7035	Headgear	1 per 6 months

* The Medicare schedule states how often accessories can be replaced if medically necessary. Check with all payers directly to verify current coverage requirements and policies.

AirFit N30

Nasal cradle mask

Quick fitting guide

ResMed AirFit™ N30 features a nasal cradle cushion with a front-end tube connection. It's easy to use and well-balanced, with unique features that maximize comfort and performance while minimizing facial contact. And with a small, simple and discreet design, AirFit N30 is designed to make adjusting to and staying compliant with therapy easier than ever.

Refer to the user guide for further instructions and complete labeling information including indications, contraindications, warnings and precautions.

myAir™ Compatible with myAir™ patient engagement program when paired with most Air10™ devices.



- While holding the nasal cradle cushion away from the nose, make sure the cushion size (SML, MED or SML-WIDE), and left (L) and right (R) indicators are facing towards you. Place the nasal cradle cushion under the nose.



- Pull the QuickFit™ headgear over the head by holding the bottom band and guiding it around the back of the head. The top band should sit comfortably on the crown of the head.



- Briefly pull the nasal cradle cushion away from the nose, and then place it back under the nose to ensure the cushion fits well.



- To achieve a comfortable fit, spread the split bands on the QuickFit headgear apart to loosen the mask or draw them closer together to tighten the mask.



- To fix a mask leak, adjust the QuickFit headgear by pulling it through the single side-clip. Adjust it just enough for a comfortable seal. Do not overtighten. You may need to adjust the headgear over time.



- The mask is now fitted and ready to use.

AirFit N30

Nasal cradle mask

Systems and components

AirFit N30 Complete mask system


(Includes mask frame, nasal cradle cushion and headgear)

Small **64222**
Medium **64223**
Small-Wide **64224**

Frame system

(includes mask frame and nasal cradle cushion)

Small **64210**
Medium **64211**
Small-Wide **64212**

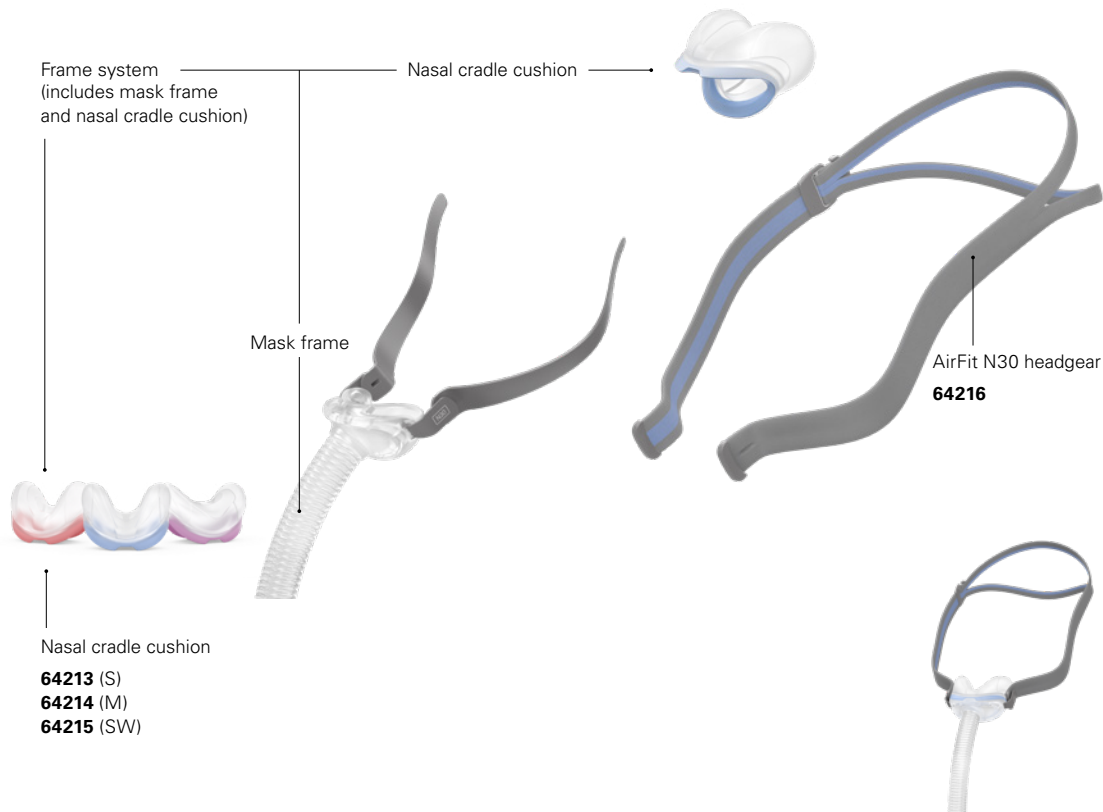

AirFit N30i

Nasal mask

Quick fitting guide

ResMed AirFit N30i features a top-of-the-head tube design and nasal cradle cushion that allows patients the freedom to sleep in the position of their choice. The SpringFit frame is designed to adapt to each patient's face, while only having two frame sizes helps streamline fitting time and SKU management.




- Pull the frame and headgear over your head until the elbow is positioned on top and the headgear is behind your head.



- Place the cushion under your nose.



- Tighten or loosen the headgear until the cushion fits comfortably under your nose.

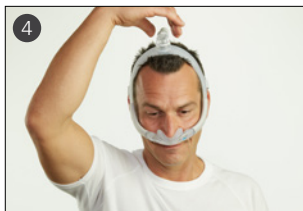

- Squeeze the side buttons on the elbow and detach from the frame.



- Connect the air tubing from your device to the elbow.



- Reattach the elbow to the frame. The mask is now ready to use.

Refer to the user guide for further instructions and complete labeling information including indications, contraindications, warnings and precautions.

myAir™ Compatible with myAir™ patient engagement program when paired with most Air10™ devices.

AirFit N30i

Nasal mask

Systems and components

AirFit N30i

Complete mask system

(Includes frame system, headgear and 2 more cushions sizes)

Standard **63800**
Small **63801**

Frame system

(frame, cushion, elbow)

Standard frame

Small **63802**
Medium **63804**
Wide **63806**
Small-Wide **63808**

Small frame

Small **63803**
Medium **63805**
Wide **63807**
Small-Wide **63809**

Frame system
(includes frame,
cushion and elbow)

Elbow

Frame

Cushion

63812 (W)
63810 (SW)
63811 (M)
63813 (S)

Headgear

63814

AirTouch N20/AirTouch N20 for Her

Nasal mask

Quick fitting guide

When it comes to effective sleep therapy, comfort is key. That's why ResMed AirTouch™ N20 nasal mask features our UltraSoft™ memory foam cushion, making it our softest nasal mask ever. AirFit N20 provides patients with a unique comfort experience with a remarkable 9 out of 10 CPAP users who tried it found the memory foam cushion to be comfortable.*

* CPAP users rated mask criteria on a Likert Scale of 0-10. A 'comfortable mask' is equal to a score of 8 or higher.

1 ResMed guided external clinical study of 25 evaluable, existing CPAP mask users from Mar. 15 to Apr. 5, 2018, in Bella Vista, AU who trialed ResMed AirTouch N20 and ResMed AirFit N20 at home for 7 nights each in place of their prescribed nasal mask.

myAir™ Compatible with myAir™ patient engagement program when paired with most Air10™ devices.



- Twist and pull both magnetic clips away from the frame.



- With both lower headgear straps released, hold the mask against the patient's face.
- Pull the headgear over the patient's head, ensuring that the ResMed logo on the headgear is facing up.



- Bring the lower straps below the patient's ears and attach the magnetic clip to the frame.



- Undo the fastening tabs on the upper headgear straps.
- Adjust the straps evenly until the mask sits comfortably just over the nose.



- Undo the fastening tabs on the lower headgear straps.
- Adjust the straps evenly until the mask is stable.
- Reattach the fastening tabs.



- Connect the air tube on the device to the short tube.
- Attach the elbow to the mask. First, press the elbow's side buttons, and then push the elbow into the mask, ensuring it clicks in on both sides.
- Start the device.

Refer to the user guide for further instructions and complete labeling information including indications, contraindications, warnings and precautions.

AirTouch N20/AirTouch N20 for Her

Nasal mask

Systems and components

AirTouch N20 Complete mask system

Small **63903**

Includes frame, small cushion, elbow with short tube and standard blue headgear with headgear clips

Medium **63901**

Includes frame, medium cushion, elbow with short tube and standard blue headgear with headgear clips


Large **63902**

Includes frame, large cushion, elbow with short tube and standard blue headgear with headgear clips

AirTouch N20 for Her Complete mask system

Small **63900**

Includes frame, small cushion, elbow with short tube and small lavender headgear with headgear clips

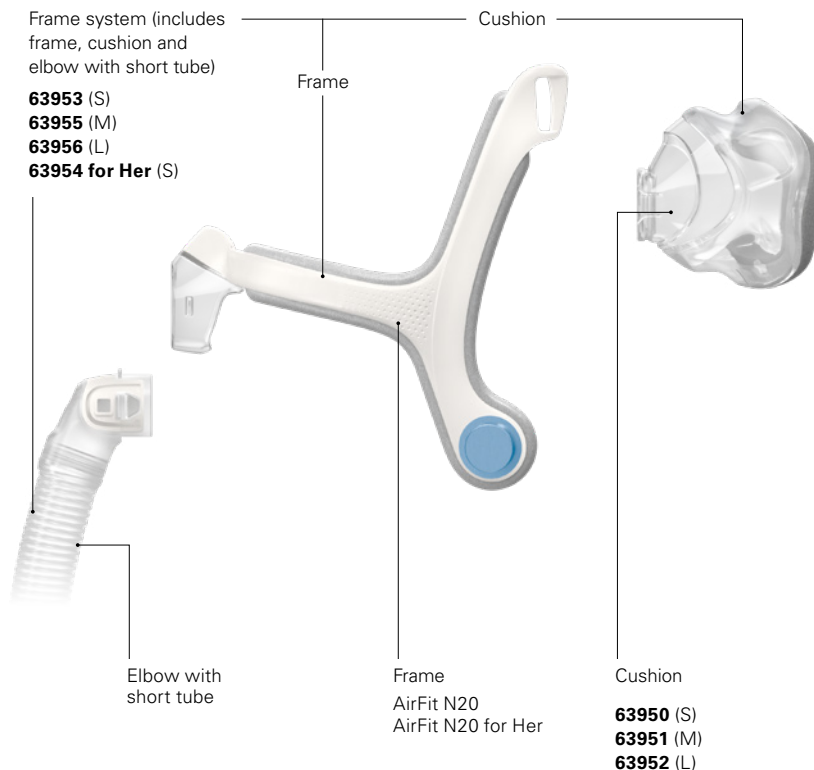

AirFit N20 headgear with headgear clips

63560 (S)

63561 (Std)

63562 (L)



AirFit N20 for Her headgear with headgear clips

63558 (S)



AirFit N20/AirFit N20 for Her

Nasal mask

Quick fitting guide

Everyone's face is different. But with our AirFit™ 20 series, fitting a broad range of patients has never been easier. That's because each mask has been specially designed to fit all patients regardless of facial profile, taking the constant guesswork out of mask selection. And with a recent international study showing the AirFit N20 fits 99% of patients,¹ you can trust that it will perform for your diverse group of patients.

¹ ResMed AirFit N20 internal international fitting study of new and existing CPAP patients, conducted Nov 2015.

Compatible with AirMini™ when used with the N20/AirMini connector elbow.

myAir™ Compatible with myAir™ patient engagement program when paired with most Air10™ devices.



- Start by ensuring the magnets on the lower straps are unclipped.
- If they are clipped, detach them by twisting and pulling them away from the frame.



- Pull the headgear over the head, making sure the logo is facing outward.



- Position the lower straps so they both sit below the ears.
- Bring each of the magnets up to meet their corresponding clips on the frame.



- Adjust the fastening tabs on the upper straps until the mask feels even and stable.
- Reattach the fastening tabs.



- Repeat step 4 with the lower straps.
- Reattach the fastening tabs.



- Connect the device's air tube to the elbow, then click the elbow into the front of the mask.
- Start the device.

Refer to the user guide for further instructions and complete labeling information including indications, contraindications, warnings and precautions.

AirFit N20/AirFit N20 for Her

Nasal mask

Systems and components

AirFit N20 Complete mask system

Small **63503**

Includes frame, small cushion, elbow with short tube and standard blue headgear with headgear clips

Medium **63501**

Includes frame, medium cushion, elbow with short tube and standard blue headgear with headgear clips

Large **63502**

Includes frame, large cushion, elbow with short tube and standard blue headgear with headgear clips

AirFit N20 for Her Complete mask system

Small **63500**

Includes frame, small cushion, elbow with short tube and small lavender headgear with headgear clips

Frame system (includes frame, cushion and elbow with short tube)

63553 (S)


63555 (M)

63556 (L)

63554 for Her (S)

Frame

Cushion

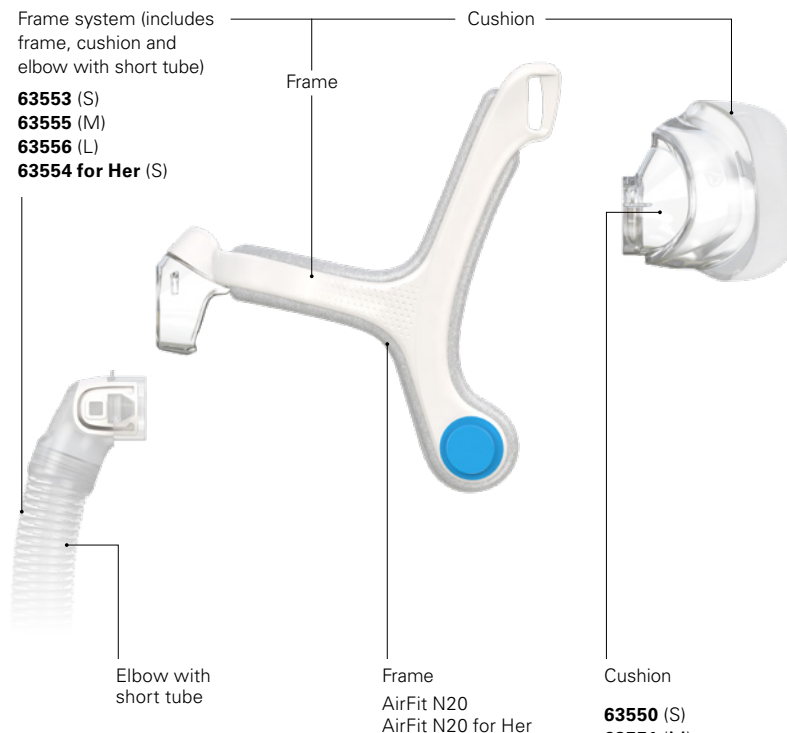

Elbow with short tube

Frame
AirFit N20
AirFit N20 for Her

Cushion
63550 (S)
63551 (M)
63552 (L)



AirFit N20 headgear with headgear clips

63560 (S)

63561 (Std)

63562 (L)



AirFit N20 for Her headgear with headgear clips

63558 (S)



AirFit N10/AirFit N10 for Her

Nasal mask

Systems and components

AirFit N10 (Small) Complete mask system

(Complete system includes frame, small cushion, standard blue headgear, soft sleeves and headgear clips)

63229

AirFit N10 (Standard) Complete mask system

(Complete system includes frame, standard cushion, standard blue headgear, soft sleeves and headgear clips)

63200

AirFit N10 (Wide) Complete mask system

(Complete system includes frame, wide cushion, standard blue headgear, soft sleeves and headgear clips)

63202

AirFit N10 for Her (Small) Complete mask system

(Complete system includes frame, small cushion, small pink headgear, soft sleeves and headgear clips)

63201

Frame system
(includes frame,
cushion, soft sleeves
and headgear clips)

63231 (S)

63230 (Std)

63232 (W)

Cushion

Frame



AirFit N10 headgear

63262 (S)

63260 (Std)



AirFit N10 for Her headgear


63261 (S)

Cushion

63241 (S)

63240 (Std)

63242 (W)

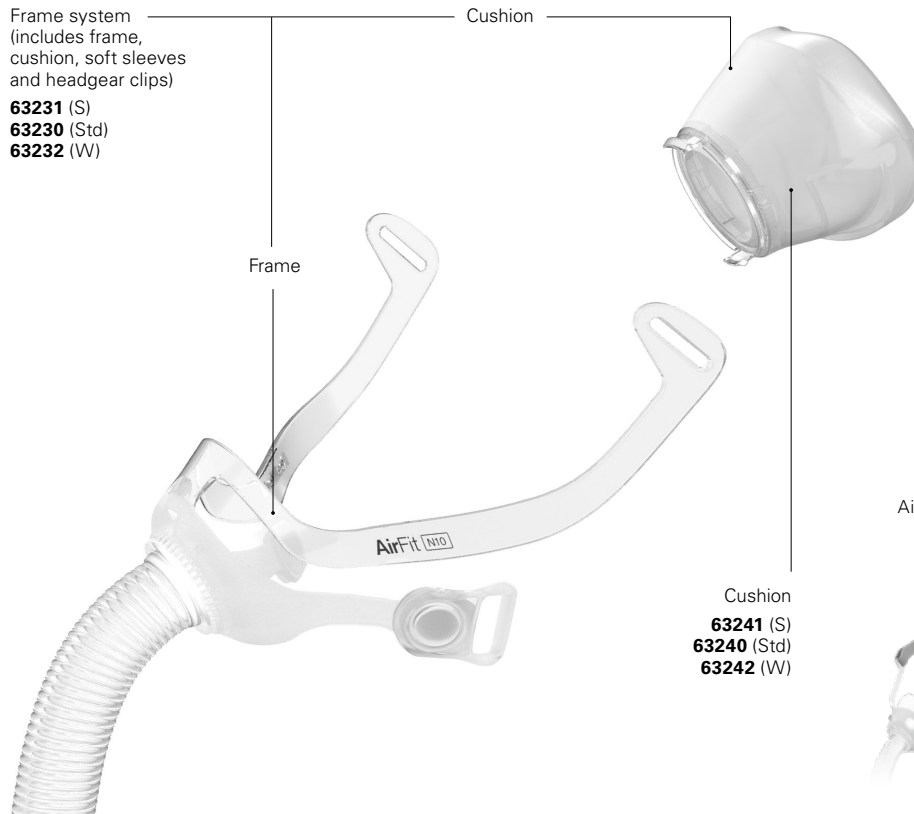

Mirage FX/Mirage FX for Her

Nasal mask

Systems and components

Mirage FX Complete mask system

Mirage FX **62103**
Mirage FX (W) **62118**

Frame system (no headgear)

Mirage FX **62112**
Mirage FX (W) **62126**

Headgear

Mirage FX **62110**
Mirage FX (S) **62138**

Mirage FX for Her Complete mask system




Mirage FX for Her **62128**
Mirage FX
for Her (S) **62109**

Frame system (no headgear)

Mirage FX for Her **62112**
Mirage FX
for Her (S) **62137**

Headgear

Mirage FX for Her **62129**

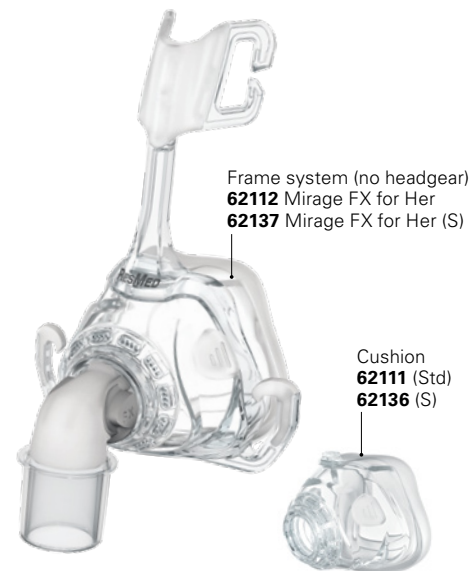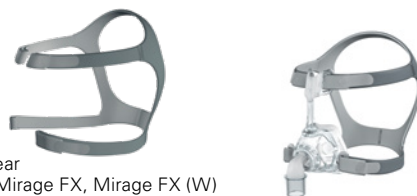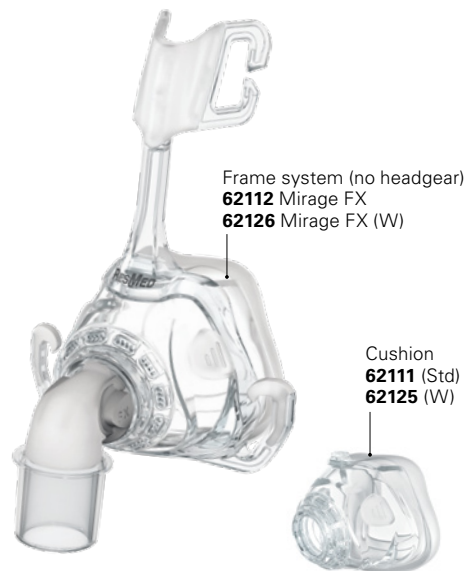

Mirage Activa LT

Nasal mask

Systems and components

Complete mask system¹

Small	60182
Medium	60148
Large	60149
Large Wide	60150

Frame system²

(no headgear)

Small	60195
Medium	60172
Large	60173
Large Wide	60174

Frame assembly²⁻³

(no headgear, no cushion)

Small	60196
Standard	60175


Headgear

(no clips)

Small	16118
Standard	16117
Large	16119

Elbow assembly

(elbow, swivel)

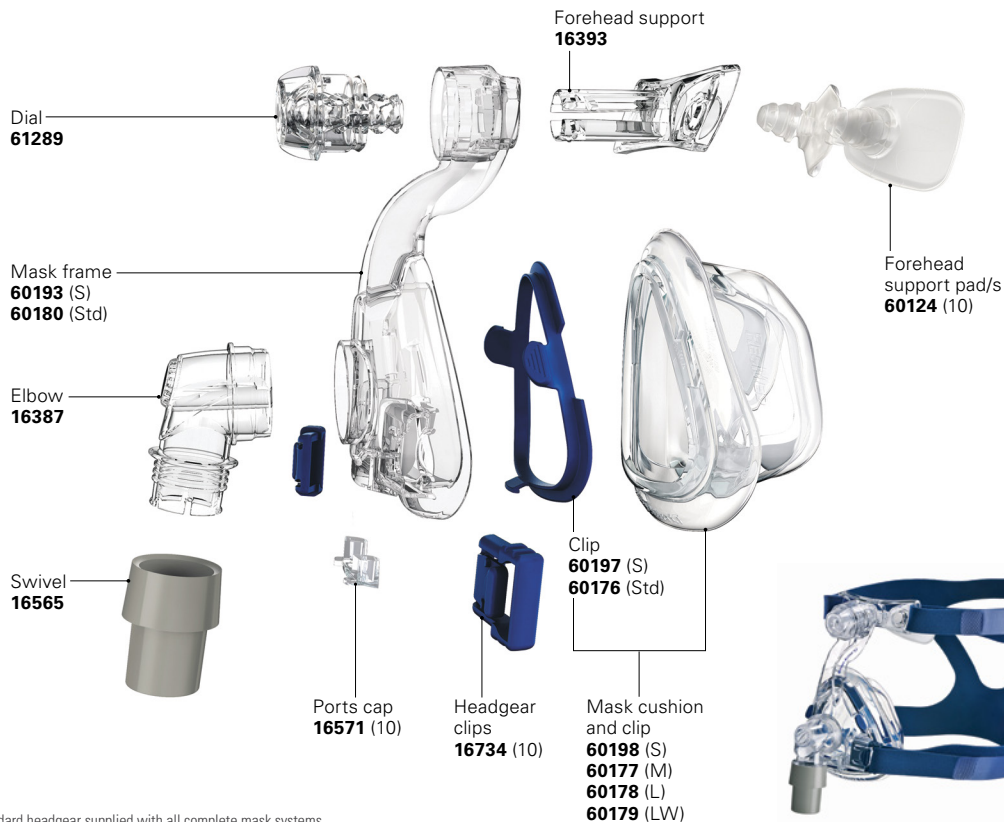

1 Standard headgear supplied with all complete mask systems

2 Includes headgear clips

3 Small cushion compatible with Small frame only

Swift FX Nano/Swift FX Nano for Her

Nasal mask

Systems and components

Swift FX Nano (Standard) Complete mask system

(Complete system includes Standard cushion, Swift FX Nano headgear assembly, plus gray backstrap, soft wraps and frame assembly)

62200

Swift FX Nano (Wide) Complete mask system


(Complete system includes Wide cushion, Swift FX Nano headgear assembly, plus gray backstrap, soft wraps and frame assembly)

62251

Swift FX Nano for Her (Small) Complete mask system

(Complete system includes Small cushion, Swift FX Nano For Her headgear assembly, plus pink backstrap, soft wraps and frame assembly)

62201

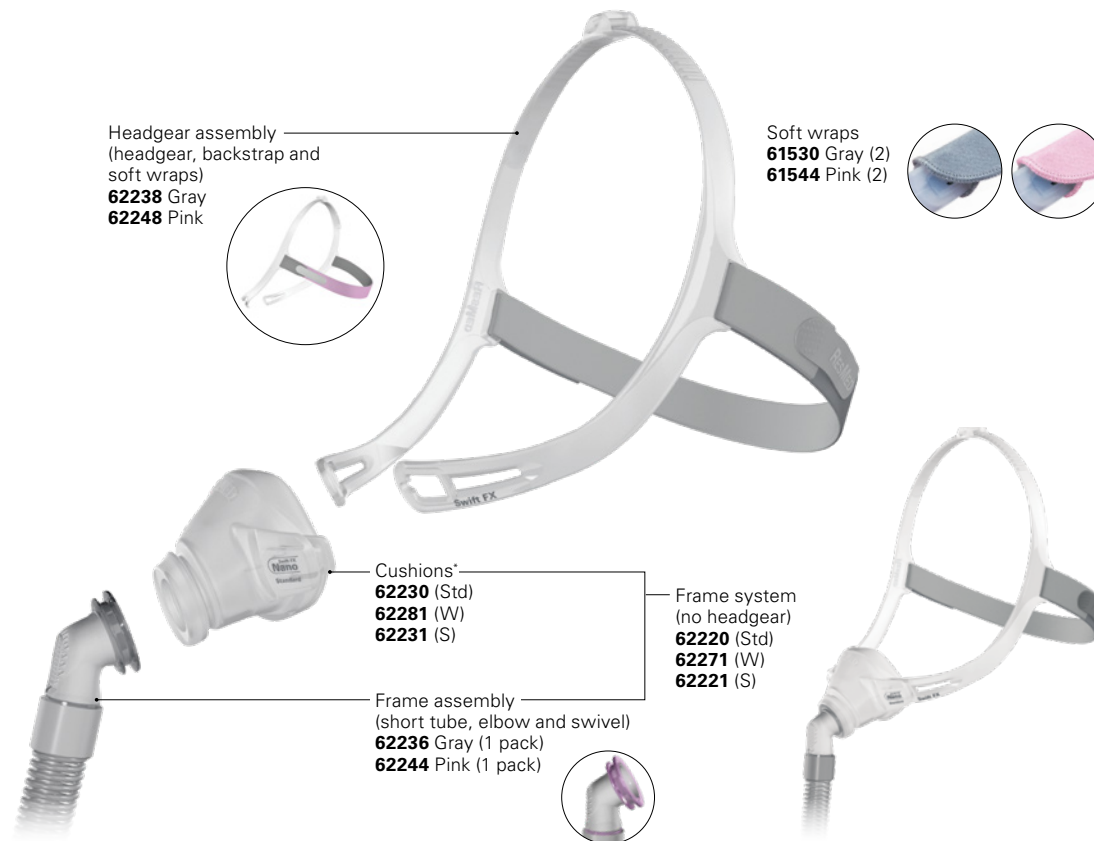

* Note: Nano Cushions are coded by PDAC as HCPCS A7034

Nasal pillows masks

MASK LINEUP



AirFit™ P10
AirFit™ P10 for Her

Minimalist category



AirFit™ P30i

Freedom category

OTHER MASKS INCLUDE:

Swift™ FX/Swift™ FX for Her

Swift™ FX Bella Gray/Swift™ FX Bella

For all legacy mask information, visit [ResMed.com/MaskGuide](https://www.resmed.com/MaskGuide).

HCPCS codes for ResMed's nasal pillows*

CODE	PRODUCT	REPLACEMENT SCHEDULE
A7034	Nasal mask system	1 per 3 months
A7033	Nasal pillow replacement	2 pairs per month
A7035	Headgear	1 per 6 months

* The Medicare schedule states how often accessories can be replaced if medically necessary. Check with all payers directly to verify current coverage requirements and policies.

AirFit P30i

Nasal pillows mask

Quick fitting guide

A new tube-up option

ResMed AirFit™ P30i offers a nasal pillows cushion option for patients who prefer a top-of-the-head tube design. Just like AirFit N30i, it features a SpringFit frame designed to adapt to each patient's face and only two frame sizes, which streamlines fitting and SKU management.

Patients can also easily disconnect from therapy if they have to get up at night using the quick-release elbow, a feature that ensures parts don't get stuck or lost.



- Pull the frame and headgear over the head until the elbow is positioned on top and the headgear is behind the head.



- Place the cushion under the nose and align the pillows into the nostrils.



- Tighten or loosen the headgear until the cushion fits comfortably under the nose.



- Squeeze the side buttons on the elbow and detach from the frame. Connect the air tubing from the device to the elbow, then reattach to the frame.

Refer to the user guide for further instructions and complete labeling information including indications, contraindications, warnings and precautions.

myAir™ Compatible with myAir™ patient engagement program when paired with most Air10™ devices.

1 ResMed guided external clinical study of 23 current nasal pillows mask users from May to June, 2018 in Sydney, Australia who trialed ResMed AirFit P30i and Phillips Respironics DreamWear gel pillows in randomized order. One participant withdrew consent after study completion (n=22)

AirFit P30i

Nasal pillows mask

Systems and components

AirFit P30i Complete mask system

(frame system, headgear and additional cushion sizes)

Standard **63850**
Small **63851**

Frame system

(frame, cushion, elbow)

Standard frame
Small **63856**
Medium **63854**
Large **63852**

Small frame
Small **63857**
Medium **63855**
Large **63853**

Frame system
(includes frame,
cushion and elbow)

Elbow

Frame

Cushion

63861 (S)
63862 (M)
63863 (L)

Headgear
63814

AirFit P10/AirFit P10 for Her

Nasal pillows mask

Quick fitting guide

A quieter mask to fit your sleep

- The AirFit P10 offers the reliability and performance of the Swift FX in a design that's approximately 50% lighter and 50% quieter.¹
- A study found that patients slept more than 40 minutes longer per night when using the AirFit P10.²
- The QuietAir™ woven-mesh venting gently diffuses exhaled air, so it's barely noticeable to patients.

Light and comfortable

- Comfortable and unobtrusive, the AirFit P10 is 50% lighter than Swift FX.¹
- Lightweight flexible tube stretches and bends to allow for flexibility without sacrificing stability.
- Dual-wall nasal pillows with trampoline action and adjustable, split-strap headgear let patients move around freely without compromising seal.

Easy to use

- The QuickFit™ elastic headgear is auto-adjusting and simple to fit, and the mask easily slips on and off without any extra steps.
- Key-clip nasal pillows seal on contact with a simple one-click action that makes assembling and cleaning the mask fast and easy.
- Each of the four pillows sizes is a different color for easy identification: pink (extra small), clear (small), gray (medium) and blue (large).

Compatible with the AirMini™

myAir™ Compatible with myAir™ patient engagement program when paired with most Air10™ devices.




- Hold the pillows at your nose, ensuring that the pillow size indicator is on the bottom and you can see the left (L) and right (R) orientation indicators.



- Stretch the lower strap around the back of your head.
- The top strap should sit comfortably on the top of your head.



- Pull the headgear over your head by holding the lower strap.
- If adjustment is needed, simply spread the upper and lower straps apart to loosen, or draw them closer together to tighten.

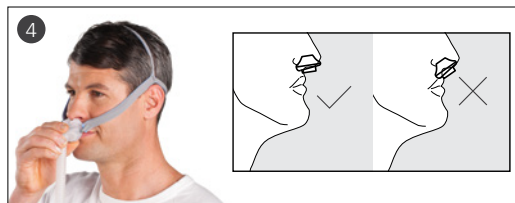

- Move the mask away from your face, and then place it back again.
- Make sure the nasal pillows fit securely into your nostrils with both pillows positioned upright.



- Finally, connect the end of the mask's tube into your device's air tubing.
- Turn on your device – your mask is now ready to use.

¹ ResMed Ref: D000-0201; eLB0798-206.

² ResMed Clinical Study Comparing AirFit P10 to Swift FX; FECS3 Ref: D629-117 (page 11), ResMed objective data (pressure, leak, usage, AHI). Refer to the user guide for further instructions and complete labeling information including indications, contraindications, warnings and precautions.

AirFit P10/AirFit P10 for Her

Nasal pillows mask

Systems and components

AirFit P10 Complete mask system


(Complete system includes frame, pillows and blue QuickFit elastic headgear. S, M and L pillows included.)

62900

AirFit P10 for Her Complete mask system

(Complete system includes frame, pillows and pink QuickFit elastic headgear. XS, S and M pillows included.)

62910

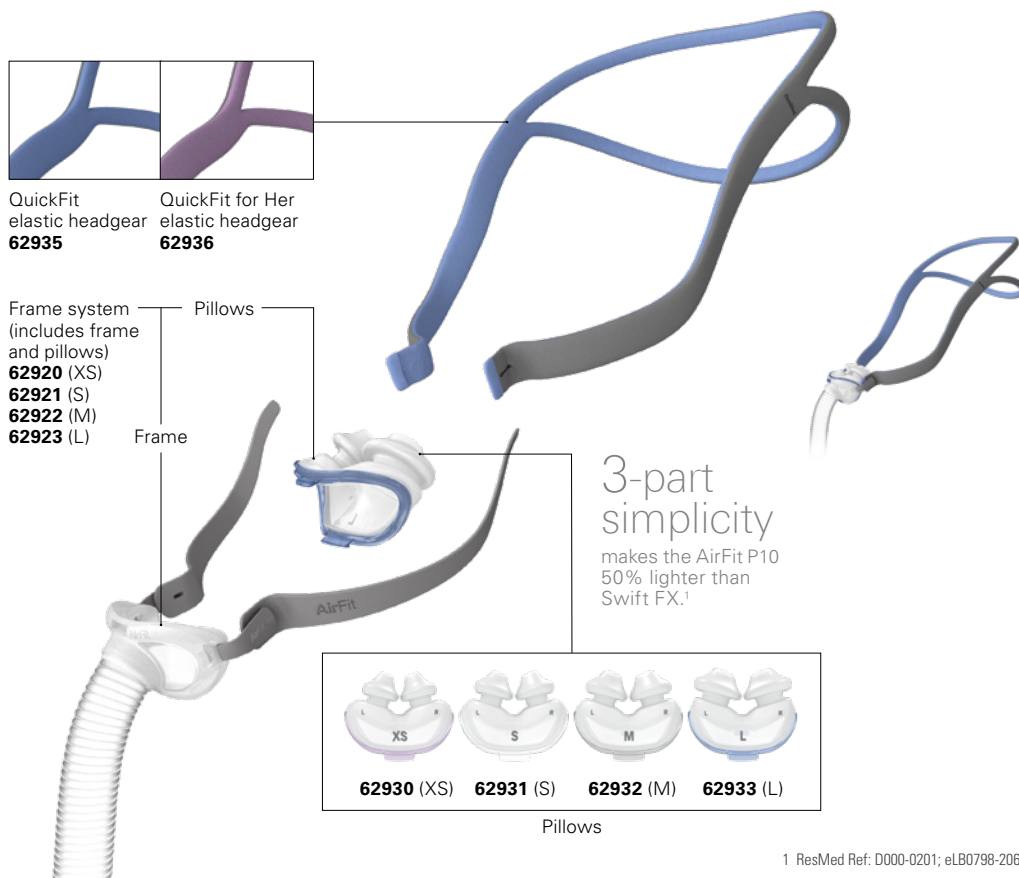

¹ ResMed Ref: D000-0201; eLB0798-206.

Swift FX/Swift FX for Her

Nasal pillows mask

Systems and components

Swift FX

Complete mask system

(Includes frame, headgear, soft wraps, backstrap and Small, Medium, and Large pillows.)

61500

Frame system

(no headgear)

Extra Small¹ **61510**

Small **61511**

Medium **61512**

Large **61513**

Headgear assembly

(headgear, backstrap and soft wraps)

61529

Swift FX for Her

Complete mask system

(Extra Small, Small and Medium pillows, pink backstrap and pink soft wraps)

61540

Frame system

(no headgear)

Extra Small **61510**

Small **61511**

Medium **61512**

Large² **61513**


Headgear assembly

(headgear, pink and grey backstraps and pink soft wraps)

61543



- 1 Extra Small (XS) pillows sold separately
2 Large (L) pillows sold separately

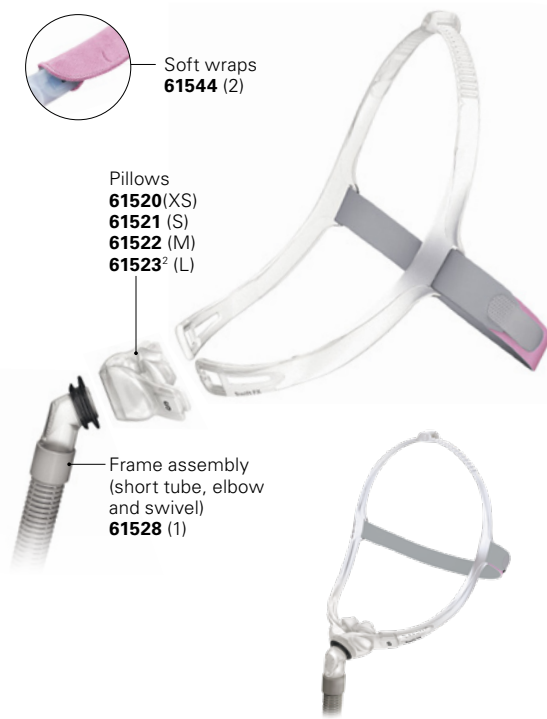

Swift FX Bella Gray/Swift FX Bella

Nasal pillows mask

Systems and components

Swift FX Bella Gray Complete mask system



(Swift FX Bella Gray complete system includes Bella Gray loops attached to the pillows frame, Swift FX headgear with blue backstrap and gray soft wraps. Also included are Small [S], Medium [M] and Large [L] pillows.¹)

61568

Swift FX Bella Complete mask system

(Swift FX Bella complete system includes Bella loops attached to the pillows frame, Swift FX for Her headgear with pink backstrap and pink soft wraps. Also included are Extra Small [XS], Small [S] and Medium [M] pillows.²)

61560

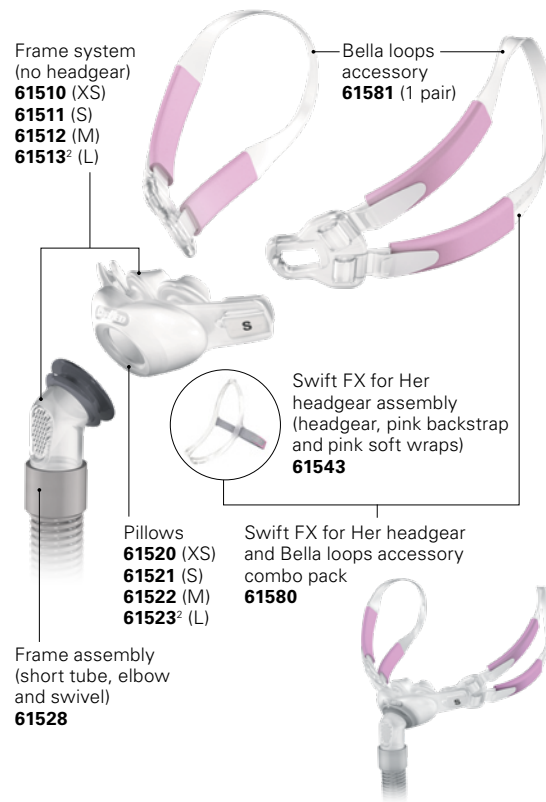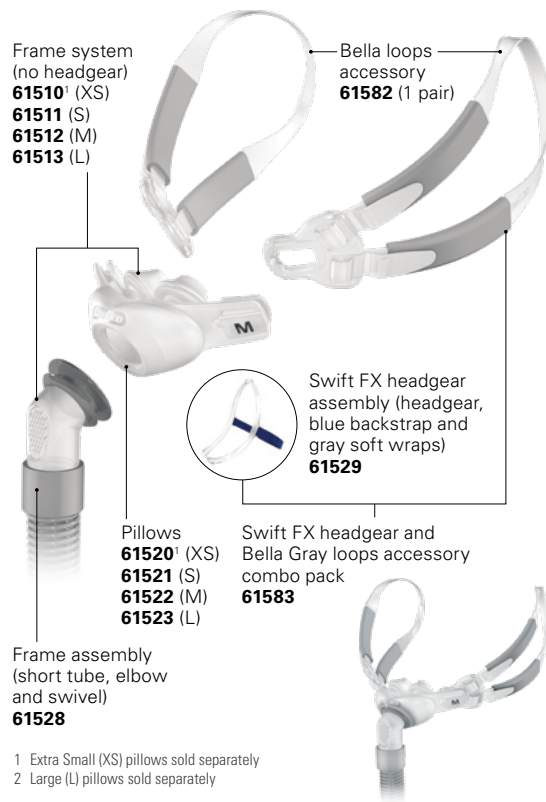

Pediatric nasal masks

MASK LINEUP



HCPCS codes for ResMed's pediatric masks*

CODE	PRODUCT	REPLACEMENT SCHEDULE
A7034	Nasal mask system	1 per 3 months
A7032	Nasal cushion replacement	2 per month
A7035	Headgear	1 per 6 months

* The Medicare schedule states how often accessories can be replaced if medically necessary. Check with all payers directly to verify current coverage requirements and policies.

Pixi

Pediatric nasal mask

Quick fitting guide

Three points of adjustment

- The mask can be adjusted at the top, side and bottom straps, providing a greater range of adjustment to cater to children of various ages.

Out of sight

- The headgear is designed to sit away from the child's eyes and ears to reduce irritation and obstruction.

Clear line of sight

- Not only can the child see clearly while wearing the mask, but caregivers are also reassured and comforted by being able to see their child's eyes and more of their face.

Fitting tip

- When fully fitted, ensure that the top strap is not overtightened and is aligned with the top of the ears; the middle strap should sit across the back of the head aligned with the bottom of the ears.



Refer to the user guide for further instructions and complete labeling information including indications, contraindications, warnings and precautions.

Pixi hint: If fitting the mask for the first time, we recommend the child be sitting in an upright position to ensure correct fit.




- Hold the mask gently against the child's face, checking that the nose is fully enclosed by the cushion and the mouth is clear.



- Pull the headgear over the head and move the middle strap down until it sits across the middle of the back of the head.



- Keeping your hand on the frame to maintain correct headgear position, bring the lower headgear strap behind the back of the neck.

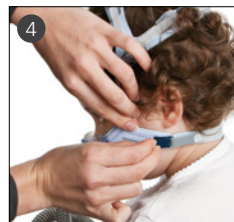

- Attach the lower headgear strap by holding the tab and hooking the quick-release clip onto the mask.



- If necessary, adjust the straps by removing the Velcro and gently tightening the headgear straps until comfortable. Do not overtighten.



- Connect one end of the air tubing to the swivel and the other end to the airflow device.
- Turn on the airflow device.

- For optimal results, adjust the straps starting at the bottom strap and work upward.


Pixi

Pediatric nasal mask

Systems and components

Complete mask system
61030

Frame system **61035**
(no headgear)

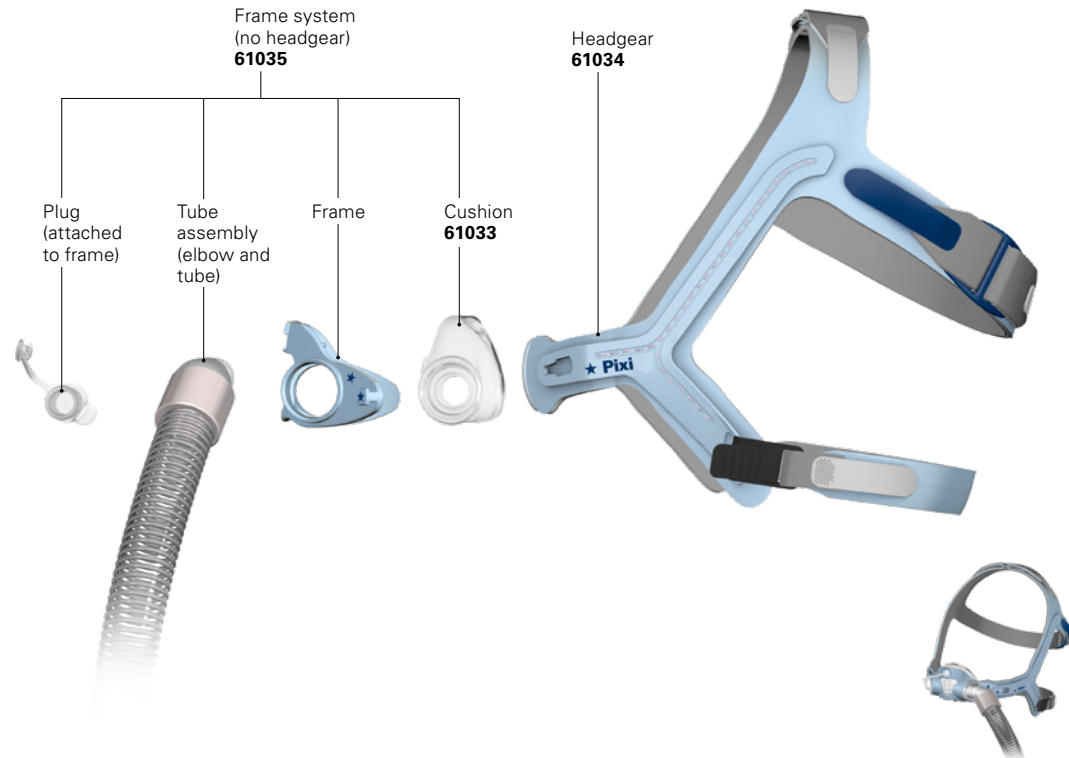

Pressure vs. flow rates for ResMed masks

These tables show the airflow rates through mask vents at varying mask pressures.
This information does not replace the relevant device user guide.

Full Face Masks

Pressure (cm H ₂ O)	AirTouch™ F20/ AirTouch™ F20 for Her	AirFit™ F20/AirFit™ F20 for Her	AirFit™ F10/AirFit™ F10 for Her	Quattro™ Air/Quattro™ Air for Her	AirFit™ F30i	AirFit™ F30	Mirage Quattro™	Quattro™ FX	Quattro™ FX for Her	Ultra Mirage	Mirage Liberty™
3	19	-	19	-	-	-	-	-	-	-	-
4	22	22	22	22	22	22	22	22	22	22	22
8	32	-	32	32	32	32	32	32	32	32	32
10	37	-	37	-	-	-	-	-	-	-	-
11	-	38	-	-	-	-	-	-	-	-	-
12	41	-	41	41	41	41	41	41	41	41	41
16	48	-	48	48	48	48	48	48	48	48	48
17	-	48	-	-	-	-	-	-	-	-	-
20	54	-	54	54	54	54	54	54	54	54	54
24	60	59	60	60	60	60	60	60	60	60	60
28	66	-	66	66	66	66	66	66	66	66	66
30	69	67	69	-	69	69	69	69	69	69	69
32	72	-	72	-	72	72	72	72	72	72	72
36	77	-	77	-	77	77	77	77	77	77	77
40	82	-	82	-	82	82	82	82	82	82	82

Nasal Masks

Pressure (cm H ₂ O)	AirFit™ N30	AirFit™ N30i	AirTouch™ N20 AirFit™ N20/AirFit™ N20 for Her/ AirFit™ N20 Classic	AirFit™ N10/AirFit™ N10 for Her	Mirage™ FX/Mirage FX for Her	Mirage™ SoftGel	Mirage Activa™ LT	Mirage Micro™	Mirage Activa™	Mirage Kidsta™	Swift™ FX Nano	Swift™ FX Nano for Her	Ultra Mirage™ II
4	20	20	20	20	19	-	-	-	-	21	20	18	-
7	-	-	27	-	-	-	-	-	-	-	-	-	-
8	29	-	-	29	28	-	-	-	-	30	29	29	-
9	-	31	-	-	-	-	-	-	-	-	-	-	-
11	-	-	35	35	-	-	-	-	-	-	-	-	-
12	37	-	-	37	34	-	-	-	-	36	37	38	-
14	-	-	40	-	-	-	-	-	-	-	-	-	-
15	-	41	-	-	-	-	-	-	-	-	-	-	-
16	43	-	-	43	40	-	-	-	-	42	43	46	-
17	-	-	44	44	-	-	-	-	-	-	-	-	-
20	49	49	49	49	45	-	-	-	-	47	49	53	-
22.5	-	-	53	53	-	-	-	-	-	-	52	-	-
24	-	-	54	54	-	-	-	-	-	-	-	-	-
25	-	55	-	-	-	-	-	-	-	-	-	-	-
27	-	-	58	-	-	-	-	-	-	-	-	-	-
30	-	-	61	61	-	-	-	-	-	-	61	-	-

Nasal Pillows

Pressure (cm H ₂ O)	AirFit™ P30i	AirFit™ P10/AirFit™ P10 for Her	Swift™ FX/Swift™ FX for Her	Swift™ FX Bella Gray	Swift™ FX Bella	Mirage™ Swift™ II	Swift™ LT/Swift™ LT for Her
4	20	20	20	20	20	20	20
8	-	29	29	29	29	29	29
9	31	-	-	-	-	-	-
12	-	37	37	37	37	37	37
15	41	-	-	-	-	-	-
16	-	43	43	43	43	43	43
20	49	49	49	49	49	49	49
25	55	-	-	-	-	-	-

Pediatric

Pressure (cm H ₂ O)	Pixi™
3	18
4	-
7	28
8	-
10	-
12	36
16	43
20	49

Cleaning the mask at home

Daily/after each use

1. Disassemble the mask according to the disassembly instructions.
2. For AirTouch: If required, wipe your cushion with a CPAP wipe or equivalent and then allow to dry. Important – do not get the cushion wet.
3. Rinse the frame and elbow under running water. Clean with a soft brush until dirt is removed. Pay close attention to the vents.
4. Soak the components (excluding AirTouch cushion) in warm water with a mild soap for up to ten minutes.
5. Shake the components in the water.
6. Brush the moving parts of the elbow and around the vent holes.
7. Brush the areas of the frame where the arms connect, and inside and outside the frame where the elbow and cushion connect.
8. Rinse the components under running water (excluding AirTouch cushion).
9. For the frame sleeves, squeeze with a clean towel to remove excess water. For the AirFit P30i mask, shake the QuietAir vent to remove any excess water.
10. Leave the components to air dry out of direct sunlight. Make sure to squeeze the arms of the frame to ensure that excess water is removed.

Weekly

1. Disassemble the mask. If your mask has magnets, they can remain attached to the headgear during cleaning.
2. Handwash the headgear in warm water with mild liquid detergent.
3. Rinse the headgear under running water. Inspect to ensure the headgear is clean and detergent free. Wash and rinse again, if necessary.
4. Squeeze the headgear to remove excess water.
5. Leave the headgear to air dry out of direct sunlight.



WARNING

- Do not use aromatic-based solutions or scented oils (e.g., eucalyptus or essential oils), bleach, alcohol or products that smell strongly (e.g., citrus) to clean any of the mask components. Residual vapors from these solutions can be inhaled if not rinsed thoroughly. They may also damage the mask, causing cracks.
- Do not iron the headgear as the material is heat-sensitive and will be damaged.

For additional information, refer to the device's user guide.

Mask/device compatibility list

This document specifies the masks compatible with ResMed's positive airway pressure devices and the mask setting, if one is required. For full details regarding the correct use of these devices, please refer to the specific User Guides. **Note:** *Not all products are available in all regions.*

Full Face Masks

Device	AcuCare F1-4	AirFit™ F30i ¹	AirFit™ F30	AirTouch™ F20 / AirTouch™ F20 for Her AirFit™ F20 / AirFit™ F20 for Her	AirFit™ F10 / AirFit™ F10 for Her	Mirage™	Mirage Series 2	Ultra Mirage™	Mirage Quattro™	Quattro™ FX / Quattro™ FX for Her Quattro™ Air / Quattro™ Air for Her	Mirage Liberty™	Hospital Full Face
AirMini™	X	X	✓	✓	X	X	X	X	X	X	X	X
AirStart™ 10 / AirSense™ 10	N/V	Full Face	Full Face	Full Face	Full Face	N/V	N/V	Full Face	Full Face	Full Face	Full Face	Full Face
AirCurve™ 10 VAuto / AirCurve™ 10 S / AirCurve™ 10 ST	N/V	Full Face	Full Face	Full Face	Full Face	N/V	N/V	Full Face	Full Face	Full Face	Full Face	Full Face
AirCurve™ 10 CS Pacewave	N/V	Full Face	Full Face	Full Face	Full Face	N/V	N/V	Full Face	Full Face	Full Face	Full Face	Full Face
AirCurve™ 10 CS-A Pacewave	N/V	N/V	Full Face	Full Face	Full Face	N/V	N/V	Full Face	Full Face	Full Face	Full Face	Full Face
Lumis™	N/V	N/V	N/V	Full Face	Full Face	N/V	N/V	Full Face	Full Face	Full Face	Full Face	Full Face
S9 AutoSet™ / S9 Elite™	N/V	Full Face	N/V	Full Face	Full Face	N/V	N/V	Full Face	Full Face	Full Face	Full Face	Full Face
S9 Escape™	N/V	Full Face	N/V	Full Face	Full Face	N/V	N/V	Full Face	Full Face	Full Face	Full Face	Full Face
AutoSet CS™ 2	N/V	N/V	N/V	N/V	FULL FACE	FULL FACE	FULL FACE	FULL FACE	FULL FACE	FULL FACE	FULL FACE ¹	FULL FACE
AutoSet CS™ (S9)	N/V	Full Face	N/V	Full Face	Full Face	N/V	N/V	Full Face	Full Face	Full Face	Full Face	Full Face
AutoSet CS™ (S9 PaceWave)	N/V	Full Face	N/V	Full Face	Full Face	N/V	N/V	Full Face	Full Face	Full Face	Full Face	Full Face
AutoSet CS-A™ (S9 PaceWave)	N/V	N/V	N/V	Full Face	Full Face	N/V	N/V	Full Face	Full Face	Full Face	Full Face	Full Face
C-Series Tango™	N/V	N/V	N/V	N/V	No setting required. Consult your clinician before changing to a different mask type.							
S7™ Lightweight	N/V	N/V	N/V	N/V	No setting required. Consult your clinician before changing to a different mask type.							
S7 Elite™ / S7 AutoSet Spirit	N/V	N/V	N/V	N/V	MIR FULL	MIR FULL	MIR FULL	MIR FULL	MIR FULL	MIR FULL	MIR FULL	MIR FULL
S7 AutoSet Respond™	N/V	N/V	N/V	N/V	MIRAGE	STANDARD	ULTRA	MIRAGE	MIRAGE	MIRAGE	MIR FULL	MIR FULL
S8™ Auto 25	N/V	N/V	N/V	MIR FULL	MIR FULL	MIR FULL	MIR FULL	MIR FULL	MIR FULL	MIR FULL	MIR FULL	MIR FULL
S8 AutoSet Spirit™	N/V	N/V	N/V	MIR FULL	MIR FULL	MIR FULL	MIR FULL	MIR FULL	MIR FULL	MIR FULL	MIR FULL	MIR FULL
S8 AutoSet Vantage™	N/V	N/V	N/V	MIR FULL	MIR FULL	MIR FULL	MIR FULL	MIR FULL	MIR FULL	MIR FULL	MIR FULL	MIR FULL
S8 Escape™	N/V	N/V	N/V	N/V	Requires pressure calibration at 12 cm H2O. Consult your clinician and check the user guide.							
S8 Escape™ II	N/V	N/V	N/V		Requires pressure calibration at 12 cm H2O. Consult your clinician and check the user guide.							
S8 Elite™ / S8 Elite™ II	N/V	N/V	N/V	MIR FULL	MIR FULL	MIR FULL	MIR FULL	MIR FULL	MIR FULL	MIR FULL	MIR FULL	MIR FULL
S8 AutoSet Spirit™ II	N/V	N/V	N/V	MIR FULL	MIR FULL	MIR FULL	MIR FULL	MIR FULL	MIR FULL	MIR FULL	MIR FULL	MIR FULL
S8 Compact™	N/V	N/V	N/V	N/V	No setting required. Consult your clinician before changing to a different mask type.							
VPAP™ III / VPAP™ III ST	N/V	N/V	N/V	N/V	MIR FULL	MIR FULL	MIR FULL	MIR FULL	MIR FULL	MIR FULL	MIR FULL	MIR FULL
VPAP™ III ST-A	N/V	N/V	N/V	N/V	MIR FULL	MIR FULL	MIR FULL	MIR FULL	MIR FULL	MIR FULL	MIR FULL	MIR FULL
VPAP™ IV / VPAP™ IV ST	N/V	N/V	N/V	MIR FULL	MIR FULL	MIR FULL	MIR FULL	MIR FULL	MIR FULL	MIR FULL	MIR FULL	MIR FULL
VPAP™ Adapt SV	N/V	N/V	N/V	N/V	FULL FACE	FULL FACE	FULL FACE	FULL FACE	FULL FACE	FULL FACE	FULL FACE ¹	FULL FACE
VPAP™ Adapt (S9)	N/V	Full Face	N/V	Full Face	Full Face	N/V	N/V	Full Face	Full Face	Full Face	Full Face	Full Face
VPAP™ Adapt (S9 PaceWave)	N/V	Full Face	N/V	Full Face	Full Face	N/V	N/V	Full Face	Full Face	Full Face	Full Face	Full Face
VPAP™ Malibu	N/V	N/V	N/V	N/V	FULL FACE	FULL FACE	FULL FACE	FULL FACE	FULL FACE	FULL FACE	FULL FACE	FULL FACE

Full face masks, continued

Device	AcuCare F1-4	AirFit™ F30i¹	AirFit™ F30	AirTouch™ F20 / AirTouch™ F20 for Her AirFit™ F20 / AirFit™ F20 for Her	AirFit™ F10 / AirFit™ F10 for Her	Mirage™	Mirage Series 2	Ultra Mirage™	Mirage Quattro™	Quattro™ FX / Quattro™ Air	Quattro™ FX for Her	Mirage Liberty™	Hospital Full Face
VPAP™ Auto 25	N/V	N/V	N/V	MIR FULL	MIR FULL	MIR FULL	MIR FULL	MIR FULL	MIR FULL	MIR FULL	MIR FULL	MIR FULL	MIR FULL
Auto 25 (S9)	N/V	Full Face	N/V	Full Face	Full Face	N/V	N/V	Full Face	Full Face	Full Face	Full Face	Full Face	Full Face
VPAP™ S / VPAP™ ST	N/V	N/V	N/V	MIR FULL	MIR FULL	MIR FULL	MIR FULL	MIR FULL	MIR FULL	MIR FULL	MIR FULL	MIR FULL	MIR FULL
VPAP™ S (S9) / VPAP™ ST (S9)	N/V	Full Face	N/V	Full Face	Full Face	N/V	N/V	Full Face	Full Face	Full Face	Full Face	Full Face	Full Face
VPAP™ ST-A (S9)	N/V	N/V	N/V	Full Face	Full Face	N/V	N/V	Full Face	Full Face	Full Face	Full Face	Full Face	Full Face
VPAP™ ST (S9 iVAPS) / VPAP™ ST-A (S9 iVAPS)	N/V	N/V	N/V	Full Face	Full Face	N/V	N/V	Full Face	Full Face	Full Face	Full Face	Full Face	Full Face
S9 VPAP™ Tx	N/V	Full Face²	N/V	N/V	Full Face	N/V	N/V	Full Face	Full Face	Full Face	Full Face	Full Face	Full Face
VPAP™ Tx	N/V	N/V	N/V	N/V	FULL FACE	N/V	N/V	FULL FACE	FULL FACE	FULL FACE	FULL FACE	FULL FACE	FULL FACE
Stellar™	N/V	N/V	N/V	FULL FACE	FULL FACE	N/V	N/V	FULL FACE	FULL FACE	FULL FACE	FULL FACE	FULL FACE	FULL FACE
VS Ultra™ / VS Integra™	✓	N/V	N/V	N/V	✓	N/V	N/V	✓	✓	N/V	✓	✓	✓
VS III™	N/V	N/V	N/V	N/V	✓	N/V	N/V	✓	✓	N/V	✓	✓	✓
Astral™	FULL FACE	N/V	N/V	FULL FACE	N/V	N/V	N/V	FULL FACE	FULL FACE	FULL FACE	FULL FACE	FULL FACE	FULL FACE

¹ SmartStop may not operate as described when using this mask with some CPAP or bi-level devices.

² Excluding PAC and iVAPS modes. Only compatible for pressure range 4 - 25 cm H2O.

N/V = The mask has not been validated for this device.

✓ = Compatible X = Not Compatible

Nasal Masks 1

Device	AirFit™ N30 for AirMini	AirFit™ N30i¹ / AirFit™ N30	AirFit™ N20 Classic	AirFit™ N20 / AirFit™ N20 for Her AirTouch™ N20 / AirTouch™ N20 for Her	AirFit™ N10 / AirFit™ N10 for Her	Swift™ FX Nano / Swift™ FX Nano for Her	Mirage™ FX / Mirage FX™ for Her	Mirage Activa™ LT	Mirage™ SoftGel	Hospital Nasal
AirMini™	✓	X	✓	✓	X	X	X	X	X	X
AirStart™ 10 / AirSense™ 10	X	Pillows	Pillows	Pillows	Pillows	Pillows	Nasal	Nasal	Nasal	Nasal
AirCurve™ 10 VAuto / AirCurve™ 10 S / AirCurve™ 10 ST	X	Pillows	Pillows	Pillows	Pillows	Pillows	Nasal	Nasal	Nasal	Nasal
AirCurve™ 10 CS Pacewave	X	Pillows	Pillows	Pillows	Pillows	Pillows	Nasal	Nasal	Nasal	Nasal
AirCurve™ 10 CS-A Pacewave	X	N/V	Pillows	Pillows	Pillows	Pillows	Nasal	Nasal	Nasal	Nasal
Lumis™	X	N/V	Pillows	Pillows	Pillows	Pillows	Nasal	Nasal	Nasal	Nasal
S9 AutoSet™ / S9 Elite™	X	Pillows	Pillows	Pillows	Pillows	Pillows	Nasal	Nasal	Nasal	Nasal
S9 Escape™	X	Pillows	Pillows	Pillows	Pillows	Pillows	Nasal	Nasal	Nasal	Nasal
AutoSet CS™ 2	X	N/V	N/V	N/V	N/V	N/V	VISTA¹	ACTIVA	ACTIVA	VISTA
AutoSet CS™ (S9)	X	Pillows	Pillows	Pillows	Pillows	Pillows	Nasal	Nasal	Nasal	Nasal
AutoSet CS™ (S9 PaceWave)	X	Pillows	Pillows	Pillows	Pillows	Pillows	Nasal	Nasal	Nasal	Nasal
AutoSet CS-A™ (S9 PaceWave)	X	N/V	Pillows	Pillows	Pillows	Pillows	Nasal	Nasal	Nasal	Nasal
C-Series Tango™	X	N/V	N/V	N/V	N/V	N/V	No setting required. Consult your clinician before changing to a different mask type.			
S7™ Lightweight	X	N/V	N/V	N/V	N/V	N/V	No setting required. Consult your clinician before changing to a different mask type.			
S7 Elite™ / S7 AutoSet Spirit	X	N/V	N/V	N/V	N/V	N/V	STANDARD	STANDARD	STANDARD	STANDARD
S7 AutoSet Respond™	X	N/V	N/V	N/V	N/V	N/V	STANDARD	STANDARD	STANDARD	STANDARD

Nasal masks 1, continued

Device	AirFit™ N30 for AirMini	AirFit™ N30i ¹ / AirFit™ N30	AirFit™ N20 Classic	AirFit™ N20 / AirFit™ N20 for Her / AirTouch™ N20 / AirTouch™ N20 for Her		AirFit™ N10 / AirFit™ N10 for Her	Swift™ FX Nano / Swift™ FX Nano for Her	Mirage™ FX / Mirage FX™ for Her	Mirage Activa™ LT	Mirage™ SoftGel	Hospital Nasal
S8™ Auto 25	X	N/V	N/V	SWIFT		SWIFT	SWIFT	STANDARD	ACTIVA	ACTIVA	STANDARD
S8 AutoSet Spirit™	X	N/V	N/V	SWIFT		SWIFT	SWIFT	STANDARD	ACTIVA	ACTIVA	STANDARD
S8 AutoSet Vantage™	X	N/V	N/V	SWIFT		SWIFT	SWIFT	STANDARD	ACTIVA	ACTIVA	STANDARD
S8 Escape™	X	N/V	N/V	N/V		Requires pressure calibration at 12 cm H2O. Consult your clinician and check the user guide.					
S8 Escape™ II	X	N/V	N/V	Requires pressure calibration at 12 cm H2O. Consult your clinician and check the user guide.							
S8 Elite™ / S8 Elite™ II	X	N/V	N/V	SWIFT		SWIFT	SWIFT	STANDARD	ACTIVA	ACTIVA	STANDARD
S8 AutoSet Spirit™ II	X	N/V	N/V	SWIFT		SWIFT	SWIFT	STANDARD	ACTIVA	ACTIVA	STANDARD
S8 Compact™	X	N/V	N/V	N/V		No setting required. Consult your clinician before changing to a different mask type.					
VPAP™ III / VPAP™ III ST	X	N/V	N/V	N/V		MIRAGE	MIRAGE	STANDARD	STANDARD	STANDARD	N/V
VPAP™ III ST-A	X	N/V	N/V	N/V		MIRAGE	MIRAGE	STANDARD	STANDARD	STANDARD	N/V
VPAP™ IV / VPAP™ IV ST	X	N/V	N/V	SWIFT		SWIFT	SWIFT	STANDARD	ACTIVA	ACTIVA	STANDARD
VPAP™ Adapt SV	X	N/V	N/V	N/V		N/V	N/V	VISTA ¹	ACTIVA	ACTIVA	VISTA
VPAP™ Adapt (S9)	X	Pillows	Pillows	Pillows		Pillows	Pillows	Nasal	Nasal	Nasal	Nasal
VPAP™ Adapt (S9 PaceWave)	X	Pillows	Pillows	Pillows		Pillows	Pillows	Nasal	Nasal	Nasal	Nasal
VPAP™ Malibu	X	N/V	N/V	N/V		SWIFT	SWIFT	VISTA ¹	N/V	N/V	MERIDIAN
VPAP™ Auto 25	X	N/V	N/V	SWIFT		SWIFT	SWIFT	STANDARD	ACTIVA	ACTIVA	STANDARD
Auto 25 (S9)	X	Pillows	Pillows	Pillows		Pillows	Pillows	Nasal	Nasal	Nasal	Nasal
VPAP™ S / VPAP™ ST	X	N/V	N/V	SWIFT		SWIFT	SWIFT	STANDARD	ACTIVA	ACTIVA	STANDARD
VPAP™ S (S9) / VPAP™ ST (S9)	X	Pillows	Pillows	Pillows		Pillows	Pillows	Nasal	Nasal	Nasal	Nasal
VPAP™ ST-A (S9)	X	N/V	Pillows	Pillows		Pillows	Pillows	Nasal	Nasal	Nasal	Nasal
VPAP™ ST (S9 iVAPS) / VPAP™ ST-A (S9 iVAPS)	X	N/V	Pillows	Pillows		Pillows	Pillows	Nasal	Nasal	Nasal	Nasal
S9 VPAP™ Tx	X	Pillows ²	N/V	N/V		Pillows	Pillows	Nasal	Nasal	Nasal	Nasal
VPAP™ Tx	X	N/V	N/V	N/V		PILLOWS	PILLOWS	NASAL	NASAL	NASAL	NASAL
Stellar™	X	N/V	PILLOWS	PILLOWS		PILLOWS ¹	PILLOWS ¹	NASAL	NASAL	NASAL	NASAL
VS Ultra™ / VS Integra™	X	N/V	N/V	N/V		N/V	N/V	N/V	N/V	N/V	N/V
VS III™	X	N/V	N/V	N/V		N/V	N/V	N/V	✓	N/V	✓
Astral™	X	N/V	Pillows	Pillows		N/V	N/V	Nasal	N/V	Nasal	Nasal

¹ SmartStop may not operate as described when using this mask with some CPAP or bi-level devices.

² Excluding PAC and iVAPS modes. Only compatible for pressure range 4 - 25 cm H2O.

N/V = The mask has not been validated for this device.

✓ = Compatible X = Not Compatible

Nasal Masks 2

Device	Pixi™	Mirage Kidsta™	Mirage Micro™	Mirage Vista™	Ultra Mirage II	Mirage	Mirage Activa™	Modular	Papillon™	Meridian
AirMini™	X	X	X	X	X	X	X	X	X	X
AirStart™ 10 / AirSense™ 10	N/V	Nasal	Nasal	Nasal	Nasal	N/V	Nasal	N/V	N/V	Nasal
AirCurve™ 10 VAuto / AirCurve™ 10 S / AirCurve™ 10 ST	N/V	Nasal	Nasal	Nasal	Nasal	N/V	Nasal	N/V	N/V	Nasal
AirCurve™ 10 CS Pacewave	N/V	Nasal	Nasal	Nasal	Nasal	N/V	Nasal	N/V	N/V	Nasal
AirCurve™ 10 CS-A Pacewave	N/V	Nasal	Nasal	Nasal	Nasal	N/V	Nasal	N/V	N/V	Nasal
Lumis™	N/V	Nasal	Nasal	Nasal	Nasal	N/V	Nasal	N/V	N/V	Nasal
S9 AutoSet™ / S9 Elite™	N/V	N/V	Nasal	Nasal	Nasal Ultra	N/V	Nasal	Nasal	Nasal	Nasal
S9 Escape™	N/V	N/V	Nasal	Nasal	Nasal Ultra	N/V	Nasal	Nasal	Nasal	Nasal
AutoSet CS™ 2	N/V	N/V	VISTA¹	VISTA	ULTRA	N/V	ACTIVA	N/V	PAPILLON	VISTA
AutoSet CS™ (S9)	N/V	N/V	Nasal	Nasal	Nasal Ultra	N/V	Nasal	Nasal	Nasal	Nasal
AutoSet CS™ (S9 PaceWave) / AutoSet CS-A™ (S9 PaceWave)	N/V	N/V	Nasal	Nasal	Nasal Ultra	N/V	Nasal	Nasal	Nasal	Nasal
C-Series Tango™	N/V	N/V	N/V		No setting required. Consult your clinician before changing to a different mask type.					N/V
S7™ Lightweight	N/V	N/V	N/V		No setting required. Consult your clinician before changing to a different mask type.					N/V
S7 Elite™ / S7 AutoSet Spirit	N/V	N/V	STANDARD	STANDARD	ULTRA	MIRAGE	STANDARD	STANDARD	STANDARD	N/V
S7 AutoSet Respond™	N/V	N/V	STANDARD	STANDARD	ULTRA	MIRAGE	STANDARD	STANDARD	STANDARD	N/V
S8™ Auto 25	N/V	N/V	STANDARD	STANDARD	ULTRA	MIRAGE	ACTIVA	N/V	N/V	STANDARD
S8 AutoSet Spirit™	N/V	N/V	STANDARD	STANDARD	ULTRA	MIRAGE	ACTIVA	STANDARD	STANDARD	STANDARD
S8 AutoSet Vantage™	N/V	N/V	STANDARD	STANDARD	ULTRA	MIRAGE	ACTIVA	STANDARD	STANDARD	STANDARD
S8 Escape™	N/V	N/V			Requires pressure calibration at 12 cm H2O. Consult your clinician and check the user guide.					
S8 Escape™ II	N/V	N/V			Requires pressure calibration at 12 cm H2O. Consult your clinician and check the user guide.					
S8 Elite™ / S8 Elite™ II	N/V	N/V	STANDARD	STANDARD	ULTRA	MIRAGE	ACTIVA	STANDARD	STANDARD	STANDARD
S8 AutoSet Spirit™ II	N/V	N/V	STANDARD	STANDARD	ULTRA	MIRAGE	ACTIVA	STANDARD	STANDARD	STANDARD
S8 Compact™	N/V	N/V			No setting required. Consult your clinician before changing to a different mask type.					
VPAP™ III / VPAP™ III ST	N/V	N/V	STANDARD	STANDARD	ULTRA	MIRAGE	STANDARD	STANDARD	N/V	N/V
VPAP™ III ST-A	N/V	STANDARD	STANDARD	STANDARD	ULTRA	MIRAGE	STANDARD	STANDARD	N/V	N/V
VPAP™ IV / VPAP™ IV ST	N/V	N/V	STANDARD	STANDARD	ULTRA	MIRAGE	ACTIVA	N/V	N/V	STANDARD
VPAP™ Adapt SV	N/V	N/V	VISTA¹	VISTA	ULTRA	N/V	ACTIVA	N/V	N/V	VISTA
VPAP™ Adapt (S9)	N/V	N/V	Nasal	Nasal	Nasal Ultra	N/V	Nasal	Nasal	Nasal	Nasal
VPAP™ Adapt (S9 PaceWave)	N/V	N/V	Nasal	Nasal	Nasal Ultra	N/V	Nasal	Nasal	Nasal	Nasal
VPAP™ Malibu	N/V	N/V	MERIDIAN	VISTA	ULTRA	N/V	N/V	N/V	N/V	MERIDIAN
VPAP™ Auto 25	N/V	N/V	STANDARD	STANDARD	ULTRA	MIRAGE	ACTIVA	N/V	N/V	STANDARD
Auto 25 (S9)	N/V	N/V	Nasal	Nasal	Nasal Ultra	N/V	Nasal	Nasal	Nasal	STANDARD
VPAP™ S / VPAP™ ST	N/V	N/V	STANDARD	STANDARD	ULTRA	MIRAGE	ACTIVA	N/V	N/V	STANDARD
VPAP™ S (S9) / VPAP™ ST (S9) / VPAP™ ST-A (S9)	Pediatric	N/V	Nasal	Nasal	Nasal Ultra	N/V	Nasal	Nasal	Nasal	Nasal
VPAP™ ST (S9 iVAPS) / VPAP™ ST-A (S9 iVAPS)	Pediatric	N/V	Nasal	Nasal	Nasal Ultra	N/V	Nasal	Nasal	Nasal	Nasal
S9 VPAP™ Tx	Pediatric	N/V	Nasal	Nasal	Nasal Ultra	N/V	Nasal	N/V	N/V	Nasal

Device	Pixi™	Mirage Kidsta™	Mirage Micro™	Mirage Vista™	Ultra Mirage II	Mirage	Mirage Activa™	Modular	Papillon™	Meridian
Stellar™	Pediatric	NASAL	NASAL	NASAL	ULTRA	N/V	NASAL	N/V	N/V	N/V
VS Ultra™ / VS Integra™	N/V	N/V	N/V	N/V	N/V	N/V	N/V	N/V	N/V	N/V
VS III™	N/V	✓	✓	✓	✓	N/V	✓	N/V	N/V	N/V
Astral™	Pediatric	Nasal	Nasal	Nasal	Nasal Ultra	N/V	N/V	N/V	N/V	N/V

N/V = The mask has not been validated for this device.

✓ = Compatible X = Not Compatible

Nasal Pillows

Device	AirFit™ P30i ¹	AirFit™ P10 For AirMini™	AirFit™ P10 / AirFit™ P10 for Her	Mirage Swift™	Mirage Swift™ II	Swift™ FX / Swift™ FX for Her / Swift™ FX Bella / Swift™ FX Bella Gray	Swift™ LT	Swift™ LT for Her / Swift™ LT-F
AirMini™	X	✓	X	X	X	X	X	X
AirStart™ 10 / AirSense™ 10	Pillows	X	Pillows ²	N/V	Pillows	Pillows	Pillows	Pillows
AirCurve™ 10 VAuto / AirCurve™ 10 S / AirCurve™ 10 ST	Pillows	X	Pillows ²	N/V	Pillows	Pillows	Pillows	Pillows
AirCurve™ 10 CS Pacewave	Pillows	X	Pillows ²	N/V	Pillows	Pillows	Pillows	Pillows
AirCurve™ 10 CS-A Pacewave	N/V	X	Pillows ²	N/V	Pillows	Pillows	Pillows	Pillows
Lumis™	N/V	X	Pillows ²	N/V	Pillows	Pillows	Pillows	Pillows
S9 AutoSet™ / S9 Elite™	Pillows	X	Pillows	Pillows	Pillows	Pillows	Pillows	Pillows
S9 Escape™	Pillows	X	Pillows	Pillows	Pillows	Pillows	Pillows	Pillows
AutoSet CS™ 2	N/V	X	N/V	N/V	N/V	N/V	N/V	N/V
AutoSet CS™ (S9)	Pillows	X	Pillows ²	Pillows	Pillows	Pillows	Pillows	Pillows
AutoSet CS™ (S9 PaceWave)	Pillows	X	Pillows ²	Pillows	Pillows	Pillows	Pillows	Pillows
AutoSet CS-A™ (S9 PaceWave)	N/V	X	Pillows ²	Pillows	Pillows	Pillows	Pillows	Pillows
C-Series Tango™	N/V	X		No setting required. Consult your clinician before changing to a different mask type.				
S7™ Lightweight	N/V	X		No setting required. Consult your clinician before changing to a different mask type.				
S7 Elite™ / S7 AutoSet Spirit	N/V	X	MIRAGE	MIRAGE	MIRAGE	MIRAGE	MIRAGE	MIRAGE
S7 AutoSet Respond™	N/V	X	MIRAGE	MIRAGE	MIRAGE	MIRAGE	MIRAGE	MIRAGE
S8™ Auto 25	N/V	X	SWIFT	SWIFT	SWIFT	SWIFT	SWIFT	SWIFT
S8 AutoSet Spirit™	N/V	X	SWIFT	SWIFT	SWIFT	SWIFT	SWIFT	SWIFT
S8 AutoSet Vantage™	N/V	X	SWIFT	SWIFT	SWIFT	SWIFT	SWIFT	SWIFT
S8 Escape™	N/V	X		Requires pressure calibration at 12 cm H2O. Consult your clinician and check the user guide.				
S8 Escape™ II	N/V	X		Requires pressure calibration at 12 cm H2O. Consult your clinician and check the user guide.				
S8 Elite™ / S8 Elite™ II	N/V	X	SWIFT	SWIFT	SWIFT	SWIFT	SWIFT	SWIFT
S8 AutoSet Spirit™ II	N/V	X	SWIFT	SWIFT	SWIFT	SWIFT	SWIFT	SWIFT
S8 Compact™	N/V	X		No setting required. Consult your clinician before changing to a different mask type.				
VPAP™ III / VPAP™ III ST	N/V	X	MIRAGE ²	MIRAGE	MIRAGE	MIRAGE	MIRAGE	MIRAGE
VPAP™ III ST-A	N/V	X	MIRAGE ²	MIRAGE	MIRAGE	MIRAGE	MIRAGE	MIRAGE
VPAP™ IV / VPAP™ IV ST	N/V	X	SWIFT ²	SWIFT	SWIFT	SWIFT	SWIFT	SWIFT
VPAP™ Adapt SV	N/V	X	N/V	N/V	N/V	N/V	N/V	N/V
VPAP™ Adapt (S9)	Pillows	X	Pillows ²	Pillows	Pillows	Pillows	Pillows	Pillows
VPAP™ Adapt (S9 PaceWave)	Pillows	X	Pillows ²	Pillows	Pillows	Pillows	Pillows	Pillows

Device	AirFit™ P30i ¹	AirFit™ P10 For AirMini™	AirFit™ P10 / AirFit™ P10 for Her	Mirage Swift™	Mirage Swift™ II	Swift™ FX / Swift™ FX for Her / Swift™ FX Bella / Swift™ FX Bella Gray	Swift™ LT	Swift™ LT for Her / Swift™ LT-F
VPAP™ Malibu	N/V	X	SWIFT ²	SWIFT	SWIFT	SWIFT	SWIFT	SWIFT
VPAP™ Auto 25	N/V	X	SWIFT ²	SWIFT	SWIFT	SWIFT	SWIFT	SWIFT
Auto 25 (S9)	Pillows	X	Pillows ²	Pillows	Pillows	Pillows	Pillows	Pillows
VPAP™ S / VPAP™ ST	N/V	X	SWIFT ²	SWIFT	SWIFT	SWIFT	SWIFT	SWIFT
VPAP™ S (S9) / VPAP™ ST (S9)	Pillows	X	Pillows ²	Pillows	Pillows	Pillows	Pillows	Pillows
VPAP™ ST-A (S9)	N/V	X	Pillows ²	Pillows	Pillows	Pillows	Pillows	Pillows
VPAP™ ST (S9 iVAPS) / VPAP™ ST-A (S9 iVAPS)	N/V	X	Pillows ²	Pillows	Pillows	Pillows	Pillows	Pillows
S9 VPAP™ Tx	Pillows ³	X	Pillows ²	N/V	Pillows	Pillows	Pillows	Pillows
VPAP™ Tx	N/V	X	PILLOWS ²	PILLOWS	PILLOWS	PILLOWS	PILLOWS	PILLOWS
Stellar™	N/V	X	PILLOWS ^{1,2}	PILLOWS ¹	PILLOWS ¹	PILLOWS ¹	PILLOWS ¹	PILLOWS ¹
VS Ultra™ / VS Integra™	N/V	X	N/V	N/V	N/V	N/V	N/V	N/V
VS III™	N/V	X	N/V	N/V	N/V	N/V	N/V	N/V
Astral™	N/V	X	N/V	Pillows	Pillows	Pillows	Pillows	Pillows

¹ SmartStop may not operate as described when using this mask with some CPAP or bi-level devices.

² Accuracy of tidal volume/minute ventilation may reduce using the AirFit P10 mask in bilevel modes, particularly at higher pressure supports.

³ Excluding PAC and iVAPS modes. Only compatible for pressure range 4 - 25 cm H2O.

N/V = The mask has not been validated for this device.

✓ = Compatible X = Not Compatible

Disinfection and sterilization guide

Clinical use only

This guide is intended for the multi-patient use of a ResMed full face mask, nasal mask or nasal pillows system ("mask") in a sleep lab, clinic or hospital. If you use the mask as a single user in the home, refer to the User Guide for cleaning instructions. This guide describes ResMed's recommended and validated procedures for cleaning, disinfection and sterilization of the mask in accordance with ISO17664.

Note: Only masks listed in the table below have been validated for reprocessing between patients.

Mask ¹	Replace these parts with new parts between patients	High-level thermal disinfection	High-level chemical disinfection		Sterilization		Validated number of cycles ²
		ISO 15883-1 158°F (70°C)-100 min; 167°F (75°C)-30 min 176°F (80°C)-10 min; 194°F (90°C)-1 min	CIDEX™ OPA Orthophthal-aldehyde 0.55%	CIDEX™ Plus Glutaraldehyde 3.4%	STERRAD™		
					100S	NX	
Full face masks							
Mirage Liberty	Swivel, inlet tube	✓	✓		✓	✓	20
Mirage Quattro ³	None	✓	✓		✓	✓	20
Quattro FX ⁷ / Quattro FX for Her ⁷	None	✓	✓		✓	✓	20
Quattro FX NV ⁷	None	✓	✓		✓	✓	20
Ultra Mirage	Valve membrane	✓	✓	✓	✓		15
Ultra Mirage Non-vented	None	✓	✓	✓	✓	✓	15
Quattro Air ⁷ / Quattro Air for Her ⁷	None	✓	✓			✓	20
AirFit F30i AirFit F30 AirTouch F20 / AirTouch F20 for Her AirFit F20 / AirFit F20 for Her	See ResMed.com/MaskSupport for the Disinfection and Sterilization guidelines.						
AirFit F10 ⁷ / AirFit F10 for Her ⁷	None	✓	✓			✓	20
Nasal masks							
Mirage Activa LT ³	None	✓	✓		✓	✓	20
Mirage FX ⁷ / Mirage FX for Her ⁷	None	✓	✓				20
Mirage Micro ³ / Mirage Micro for Kids ³	None	✓	✓		✓	✓	20

Disinfection and sterilization guide, continued

Clinical use only

Mask¹	Replace these parts with new parts between patients	High-level thermal disinfection	High-level chemical disinfection		Sterilization		Validated number of cycles²
		ISO 15883-1 158°F (70°C)-100 min; 167°F (75°C)-30 min 176°F (80°C)-10 min; 194°F (90°C)-1 min	CIDEX™ OPA Orthophthal-aldehyde 0.55%	CIDEX™ Plus Glutaraldehyde 3.4%	STERRAD™		
					100S	NX	
Ultra Mirage II	None	✓	✓	✓			15
Ultra Mirage Non-vented	None	✓	✓	✓	✓		15
Swift FX Nano/ Swift FX Nano for Her	None	✓	✓		✓⁵	✓⁵	20
AirFit N30	See ResMed.com/MaskSupport for the AirFit N30 Disinfection and Sterilization guidelines.						
AirFit N30i	See ResMed.com/MaskSupport for the AirFit N30i Disinfection and Sterilization guidelines.						
AirTouch N20/AirTouch N20 for Her	See ResMed.com/MaskSupport for the AirTouch N20/AirTouch N20 for Her Disinfection and Serilization guidelines.						
AirFit N20 / AirFit N20 for Her	See ResMed.com/MaskSupport for the AirFit N20/AirFit N20 for Her Disinfection and Sterilization guidelines.						
AirFit N10 / AirFit N10 for Her	None	✓	✓			✓	20
Nasal pillows systems							
Swift FX / Swift FX for Her / Swift FX Bella / Swift FX Bella Gray	None	✓	✓		✓⁵	✓⁵	20
Swift LT	None	✓	✓		✓⁵	✓⁵	20
Swift LT for Her/ Swift LT-F	None	✓	✓		✓⁵	✓⁵	20
Mirage Swift II	None	✓⁴	✓	✓	✓		15
AirFit P10 / AirFit P10 for Her	None	✓	✓			✓	20
AirFit P30i	See ResMed.com/MaskSupport for the AirFit P30i Disinfection and Sterilization guidelines.						

¹ Not all masks are available in all regions. For full details regarding the correct use of these masks, please refer to the specific User Guide. For a list of available replacement parts for each mask system, check the Components Card on ResMed.com.

² If a healthcare facility requires an additional disinfection or sterilization cycle after reassembly, the number of validated cycles must be halved.

³ Alternative disassembly available. See "Disassembling the masks before cleaning."

⁴ Some components cannot be treated. Alternative cleaning must be used. See "Cleaning other components."

⁵ Cushion cannot be treated via STERRAD sterilization.

⁶ Short tube assembly cannot be treated via STERRAD sterilization.

⁷ Short tube assembly validated to 158°F (70°C).

Validated disinfection and sterilization procedures

In the procedures below, only **one** of the following three disinfection or sterilization procedures needs to be performed. Refer to ResMed.com/MaskSupport for AirFit F30i, AirFit F30, AirFit F20, AirFit N30. AirFit N30i, AirFit N20, AirFit P30i and AirTouch F20 disinfection and sterilization procedures.

	High-level thermal disinfection	High-level chemical disinfection	Sterilization
Disassembly	Disassemble the mask according to the instructions in the User Guide.		
Cleaning and drying	1. Thoroughly clean the separated mask components (excluding headgear) with a soft-bristle brush for one minute while soaking them in the detergent Alconox™ (diluted at 1%) according to the manufacturer's instructions. Pay particular attention to all crevices and cavities.		
	2. Inspect each component and, if required, repeat washing until visually clean.		
	3. Rinse the components twice by agitating them vigorously in drinking quality water (1.3 gallons [5 liters] per mask).		
	4. Allow the mask components to air-dry out of direct sunlight.		
Pre-disinfection and drying	Not validated.		
Disinfection or sterilization and drying	1. Using a certified hot water disinfection system, soak the disinfectable mask components using a temperature-time combination: ¹ ISO 15883-1: • 158°F (70°C) for 100 minutes • 167°F (75°C) for 30 minutes • 176°F (80°C) for 10 minutes • 194°F (90°C) for 1 minute 2. On completion, remove the mask component from the hot water disinfection system. 3. Allow the mask components to air-dry out of direct sunlight.	1. Fully immerse and soak the disinfectable mask components in a commercially available solution according to the manufacturer's instructions: • Ortho-phthalaldehyde 0.55% (e.g., CIDEX OPA) for 20 minutes or • Glutaraldehyde 3.4% (e.g., CIDEX Plus) for 20 minutes 2. Rinse the mask components in drinking-quality water (1.3 gallons [5 liters] per mask). 3. Allow the mask components to air-dry out of direct sunlight.	1. Package the mask components prior to sterilization as described in the manufacturer's instructions for the STERRAD Sterilization System. ResMed masks have been validated using the tray method. Note: The use of pouches is not recommended. 2. Sterilize the mask following the manufacturer's instructions for the STERRAD Sterilization System. Note: Drying is achieved as part of the sterilization process.
Inspection	Perform a visual inspection of each mask component. If any visible deterioration of a mask component is apparent (i.e., cracking, crazing, tears), the mask component should be discarded and replaced. Slight discoloration of the silicone components may occur and is acceptable.		
Reassembly	Reassemble the mask according to the instructions in the User Guide.		
Packaging and storage	Store in a dry, dust-free environment away from direct sunlight. Storage temperature: -4F–140°F (-20C–60°C).		

¹ Calculated and predicted from known thermal inactivation kinetics of vegetative micro-organisms subjected to thermal disinfection (ISO 15883 1) and they are inclusive of the time-temperature combination recommended by the APIC (Associations for Professionals in Infection Control and Epidemiology) and the RKI (Robert Koch Institute).

Disassembling the masks before cleaning

These instructions apply only to:

- Mirage Quattro
- Mirage Micro
- Mirage Activa LT

As an alternate to the disassembly instructions shown in the relevant user guide for removing the forehead support and dial from the mask frame, you can disassemble these components as follows:

1. Unscrew the dial so that the forehead support protrudes about halfway out from the mask frame.
2. Firmly push the forehead support against the mask frame until the dial clicks and detaches from the mask frame.
3. Fully unscrew and remove the dial; then pull the forehead support out of the mask frame.

Cleaning the headgear

The headgear does not require disinfection and should not be washed in disinfection chemicals. Thoroughly washing the headgear between patients is acceptable.

1. Handwash the headgear in warm water (86°F/30°C) using mild soap.
2. Rinse well and allow it to air-dry out of direct sunlight.

Cleaning other components

These instructions apply only to:

- AirFit N10/AirFit N10 for Her – soft sleeves
- AirFit F10/AirFit F10 for Her – Frame, soft sleeves
- Quattro Air/Quattro Air for Her – Frame
- Quattro FX/Quattro FX for Her – Spring frame, soft sleeves
- Quattro FX NV – Spring frame, soft sleeves
- Mirage FX/Mirage FX for Her – Frame

The component does not require disinfection and should not be washed in disinfection chemicals. Thoroughly washing the component between patients is acceptable.

1. Handwash the component in warm water (86°F/30°C) using mild soap.
2. Rinse well and allow it to air-dry out of direct sunlight.



GENERAL WARNINGS AND CAUTIONS

- ResMed cannot give any assurance that deviations from the procedures listed in this guide (e.g., exceeding the number of reprocessing cycles), and their effect on the performance of the product, will be acceptable.
- Mask components should not be subjected to autoclave or ethylene-oxide gas sterilization.
- When using detergents, disinfectants or sterilization agents, always follow the manufacturer's instructions.
- Do not iron the headgear as the material is heat-sensitive and will be damaged.

Refer to [ResMed.com/MaskSupport](https://www.resmed.com/MaskSupport) for AirFit F30i, AirFit F30, F20, AirFit N30. AirFit N30i, AirFit N20, AirFit P30i and AirTouch F20 cleaning and disassembly instructions.



9001 Spectrum Center Blvd.
San Diego, CA 92123
(800)-424-0737

Patient identity and corresponding data used are fictional.

ResMed Corp San Diego, CA, USA (858) 836-5000 or (800) 424-0737 (toll free). ResMed Ltd Bella Vista, NSW, Australia +61 (2) 8884 1000 or 1 800 658 189 (toll free). See [ResMed.com](https://www.resmed.com) for other ResMed locations worldwide. Activa, ActiveCell, AcuCare, Adapt SV, Air10, AirCurve, AirFit, AirMini, AirSense, AirStart, AirTouch, AirView, Astral, AutoRamp, AutoSet, ClimateLine, ClimateLineAir, EPR, H5i, HumidAir, HumidX, HumidX Plus, InfinitySeal, iVAPS, Kidsta, Mirage, Mirage Activa LT, Mirage FX, Mirage Kidsta, Mirage Liberty, Lumis, Mirage Micro, Mirage Quattro, Mirage Swift, Mirage Vista, Mobi, myAir, Pixi, Quattro, QuickFit, QuietAir, S8 Escape, S9, SlimLine, SmartStart, SoftEdge, Spring Air, Spring Air Technology, Stellar, Swift, Swift FX, Swift FX Bella, Swift FX Nano, Swift LT, TiControl, U-Sleep, Ultra Mirage, UltraSoft, Vista, VPAP, VS and VS Ultra are trademarks and/or registered trademarks of the ResMed family of companies. Android and Google Play are trademarks of Google, Inc., registered in the U.S. and other countries. Apple, App Store and iPhone are trademarks of Apple Inc., registered in the U.S. and other countries. Amara, DreamWear and Phillips are trademarks of Koninklijke Philips N.V, registered in the US and other countries. Specifications may change without notice. For patent and other intellectual property information, see [ResMed.com/ip](https://www.resmed.com/ip). © 2020 ResMed. 1013499/25 2020-08

ResMed.com