



A lighter choice  
for a feeling of freedom



Now there's an extra choice for everyone in the FX range with Swift™ FX Bella Gray—an unobtrusive design that's easy to fit, lighter, softer and more streamlined than ever.

Benefits for you ...

**Intuitive, quick to fit**

Once the right pillow size is selected, Swift FX Bella Gray requires minimal or often no adjustment. So it saves you time, money and additional patient visits.

**Interchangeable with Swift FX**

Swift FX Bella Gray is fully interchangeable with Swift FX (both use the same pillows and tubing) for an increased level of patient compliance.

**Easier therapy acceptance**

Designed for instant appeal, Swift FX Bella Gray is not daunting, helping users accept and continue therapy. And with very few parts, the system is easy to clean and assemble.

**Stable seal**

The Swift FX Bella Gray system stays in place and maintains seal effectively.

... and your patients

**Instant appeal**

With an unobtrusive design and minimal contact on the head, the entire Swift FX Bella Gray system is designed to enhance comfort and ease of use.

**Interchangeable with Swift FX**

It is fully interchangeable with Swift FX (both have the same pillows and tubing) which provides users with another headgear option.

**Choice of sleeping positions**

Swift FX Bella Gray is stable during movement as the user sleeps—so they can sleep on their side or their back, positioning the tube as desired.

**Quiet performance**

Innovative vent design diffuses airflow for quiet mask performance—only 25 dBA.\*



**Swift™ FX Bella Gray**  
NASAL PILLOWS SYSTEM

\* Declared dual numbered noise emission values (with uncertainty 3dBA) in accordance with ISO 4871.e 1




## Easy, comfortable style

With its unobtrusive design and sophisticated cushion technology, Swift FX Bella Gray provides a comfortable fit with a feeling of freedom.

Pillows (S, M and L) come packaged with the mask system to fit a wide range of users.

### Pillows

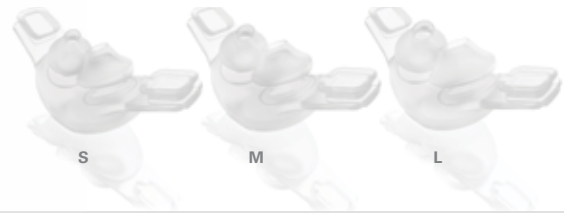

**Flexible pillows with trampoline action**  
accommodate movement in multiple directions while maintaining seal

**Dual-wall nasal pillows**  
seal on contact for a gentle, comfortable seal, offer a wide fit range and disperse air gently to the nose

**Soft, narrow cushion with integrated flexible chamber**  
positions the tubing for ease of use

**Quiet vents**  
diffuse and direct airflow quietly away from the user and bed partner

### Product Codes

Swift FX Bella Gray mask system*	<b>61569</b>
Swift FX Bella Gray headgear**	<b>61582</b>
Swift FX Bella Gray headgear and Swift FX headgear Combo**	<b>61583</b>

\* Available in Europe 1 only. Includes S, M and L pillows, Swift FX Bella Gray headgear assembly, Swift FX headgear assembly, gray backstrap and gray soft wraps.

\*\* Available in Europe 1, 2, 3.


**Europe 1:** English (UK), French, German and Italian (Switzerland)

**Europe 2:** English, Danish, Finnish, Norwegian and Swedish

**Europe 3:** English, Czech, Dutch (Belgium), French (Belgium), German (Belgium), Greek, Hungarian, Italian, Polish, Portuguese, Russian, Slovak, Spanish and Turkish



**Bella Gray headgear**  
fits comfortably around the ears with unobtrusive loops, for a streamlined look and feel

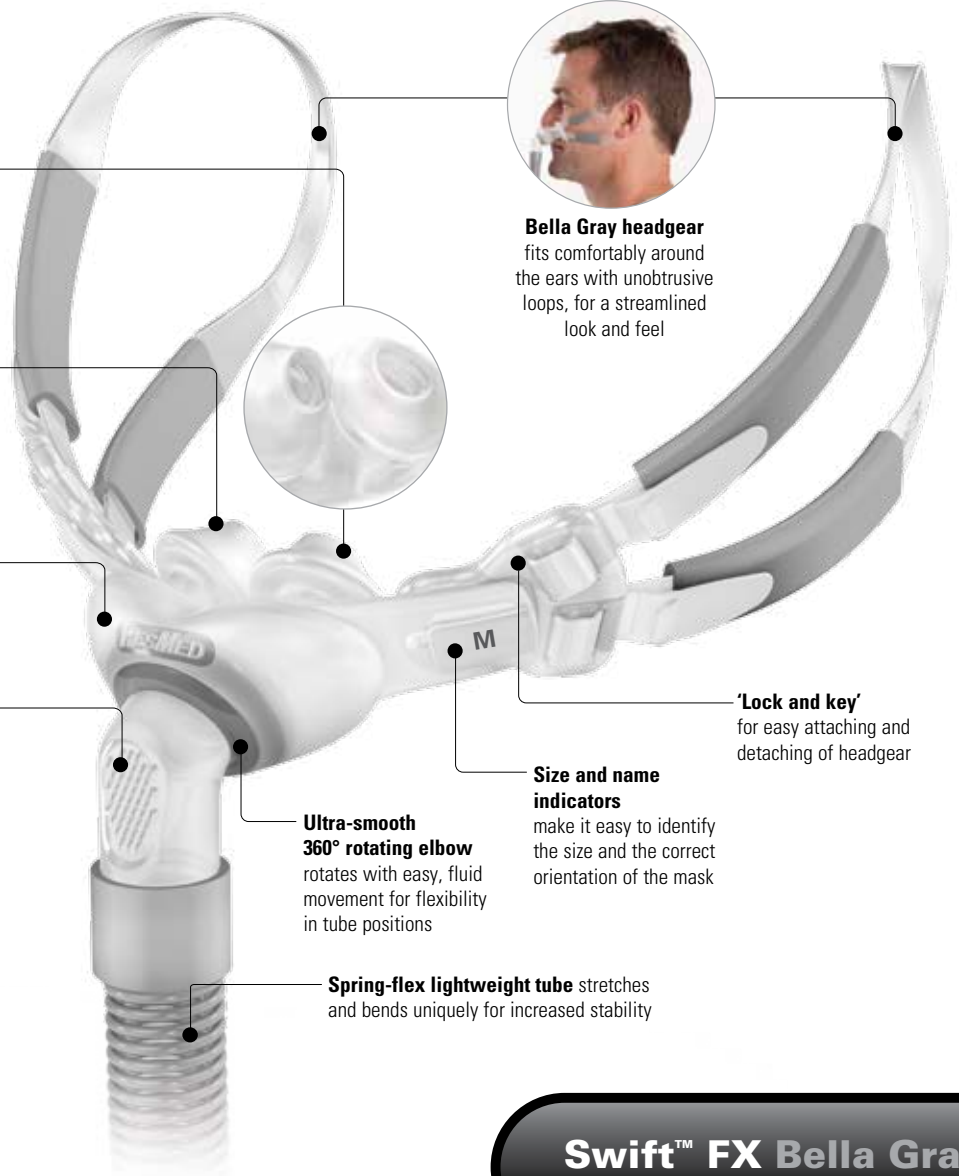

**'Lock and key'**  
for easy attaching and detaching of headgear

**Size and name indicators**  
make it easy to identify the size and the correct orientation of the mask

**Ultra-smooth 360° rotating elbow**  
rotates with easy, fluid movement for flexibility in tube positions

**Spring-flex lightweight tube** stretches and bends uniquely for increased stability

**Swift™ FX Bella Gray**  
NASAL PILLOWS SYSTEM